

HUMAN HEALTH

ENVIRONMENTAL HEALTH

CHROMERA TUTORIAL FOR NEW USERS

Chromera® Chromatography Data System



Chromera™
Tutorial for New Users
Basic Principles

Release History

Part Number	Release	Publication Date	Software Version
09936823	A	November 2009	2.0 – 2.1

Any comments about the documentation for this product should be addressed to:

User Assistance
PerkinElmer, Inc
710 Bridgeport Avenue
Shelton, Connecticut 06484-4794 U.S.A.

Or emailed to: info@perkinelmer.com

Notices

The information contained in this document is subject to change without notice.

Except as specifically set forth in its terms and conditions of sale, PerkinElmer makes no warranty of any kind with regard to this document, including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose.

PerkinElmer shall not be liable for errors contained herein for incidental consequential damages in connection with furnishing, performance or use of this material.

Copyright Information

This document contains proprietary information that is protected by copyright.

All rights are reserved. No part of this publication may be reproduced in any form whatsoever or translated into any language without the prior, written permission of PerkinElmer, Inc.

Copyright © 2009 PerkinElmer, Inc.

Trademarks

Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are protected by law.

PerkinElmer is a registered trademark of PerkinElmer, Inc.

Table of Contents

Tutorial Introduction

About this Tutorial	4
Before You Begin	5
Checklist.....	5
Other Helpful Hints.....	5

Tutorial Part 1: Chromera Manager

About the Chromera Manager	9
Starting Chromera Manager	9
Create an Instrument.....	10
About Database Management.....	12
Create a Database.....	13
Launching Chromera	14

Tutorial Part 2: Working with Example Data

About Working with Example Data.....	19
Importing Chromera Data and Methods	20
Using the Data Selector.....	24

Tutorial Part 3: Using Chromera

About the Chromera Method	29
Editing the Method Graphically.....	30
Selecting Data for Graphic Editing	30
Setting Integration Parameters	32
Setting the Bunching Factor (BF)	32
Setting the Noise Threshold (NT) and Area Threshold (AT)	32
Entering Component Information	34
Entering Baseline Timed Events	36
Specifying Component Calibration Information	38
Using the Report Format Wizard.....	41

Tutorial Part 4: Viewing Results and Reprocessing Data

Viewing Results in Post Run	51
Modes of Operation	53
Peak Identification Review	53
Go to Graphic Method Editor	53
General Views of Plot Displays	54
Single Plot	54
Stacked Plots	56
Matrix	57
Overlay	58
3-D Display.....	59
Reprocessing the Data.....	60

Tutorial Part 5: Running a Simple Analysis

Overview.....	67
Prepare Your System	67
Creating an Instrument and Database in Chromera Manager	68
Create an Instrument.....	68
Create a Database	69
Launch Chromera.....	70
Using Run Time	72
Manual Control Mode.....	72
Purge the Pump	72
Equilibrate the System.....	73
Flush the Autosampler	73
Verify that the Detector is Working	74

Tutorial Introduction

Create a New Method	74
Running the Method in the Single Run Mode.....	77
Building a Sequence	79
Using the Sequence Wizard	79
Editing Your Sequence	82
Running a Sequence	83
What If I Have Problems?	
Try This First	87
Internet Support Help	89
When All Else Fails and You Need Help	89
Register for a User's Class	89
PerkinElmer Product Specialists	89
Placing the Call	90
Contacting PerkinElmer.....	90
What Happens to My Problem or Suggestion?.....	91
What If My Problem Can't Be Resolved Immediately?.....	91
Appendix 1: Importing TotalChrom Data and Methods (Advanced Task)	
How to Import TotalChrom Data and Methods.....	95
Import TotalChrom Index Data Files.....	96
Import TotalChrom raw data Files.....	99
Import TotalChrom component name:.....	101
Copying Component Information in Chromera Method Editor	103

Tutorial Introduction

PerkinElmer would like help you to become quickly productive with this product. Before running the software, please take some time to review this tutorial. It is intended to provide you with an overview of how to get the most out of your investment and our software support organization. For information on specific Chromera functions, refer to the on-line help.

There is no magic in becoming proficient with pipettes, pH meters, and precision balances. It took practice and training for these tools to become valuable aids in your daily life. Chromera is your data-handling tool, and like any other tool in your lab, the time spent learning and practicing will make your experience much more productive and satisfying.

This tutorial does not cover all of the Chromera functions. For more training contact PerkinElmer at 1-800-762-4000 Option 4 to hear the available training options or visit PerkinElmer on the web at www.perkinelmer.com.

About this Tutorial

The Chromera Tutorial is designed for you to interact with Chromera to understand how to get started by performing basic tasks. To familiarize you with this process, the tutorial leads you through a series of procedures that are organized to illustrate your typical workflow. The workflow listed below is one way to help you get started with Chromera for the first time:

1. Using Chromera Manager, **Configure** an Instrument.
2. **Launch** the main Chromera application.
3. Investigate the **Run Time** Environment.
4. Build a **Method** for the current instrument configuration.
5. **Import** data and method information.
6. **Extract** a method from imported data.
7. Complete and optimize a method for reprocessing the imported data.
8. Batch **Reprocess** the imported data.
9. Interact with data batches in **Post Run**.
10. Create a sequence using the **Sequence Wizard**.
11. Create a report using the **Report Wizard**.

As you read through this tutorial you will receive information, including screen captures, to illustrate the workflow outlined above. Again, this documentation is not comprehensive, nor is it fully detailed. Its purpose is only to provide a suitable overview to guide you through the software and perform some of its basic functions. It does contain hints and cautions as appropriate.

After completing this tutorial you will be able to:

- Create a method - A *method* is a collection of parameters that determine how Chromera acquires and analyzes data.
- Optimize the method - *Optimizing* means adjusting method parameters to provide the best analysis of the raw data.
- Reprocess stored data with that sequence - *Reprocessing* is the re-analysis of stored data.
- Create a report format for printing the results of the analysis - A *report format* contains the specifications for the analysis report and some of the report's formatting characteristics.
- Build a sequence that utilizes the optimized method - A *sequence* is a collection of one or more rows of information that defines analysis cycles. The sequence that you create controls how data is acquired, analyzed, and/or reprocessed.

Before You Begin

The tutorial assumes that Chromera has been installed as described in the manual *Installation Instructions for Chromera Software Kits*. It also assumes that you know how to operate your instruments and that you know the basics of using Microsoft Windows.

Checklist

Complete the following preparations before you start the tutorial:

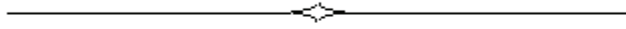
- ☐ Make sure that at the Edgeport is properly connected to the instruments.
- ☐ Turn on your printer and check the paper supply.
- ☐ Ensure that Windows is running.
- ☐ In the Windows Explorer, you may want to create a training directory on your local drive.

Other Helpful Hints

- Follow the directory and file name instructions carefully.
- When you collect data, you must provide values that are appropriate for the analysis that you are performing and for the equipment you are using.
- The defaults on your system may be different from what the tutorial describes or illustrates. The tutorial may direct you to enter or select something that is already present or selected, or where your instruments do not provide a choice.
- If there are no directions for a parameter, just leave it as is.
- When there are multiple ways to do something, the choices are listed and separated by **OR**.
- Operation of Chromera is **not compatible** with Fast User Switching and therefore **Fast User Switching must be turned off** under User Accounts in any workstation on which Chromera will be run.

Fast User Switching, is a feature of Windows XP Home Edition and Windows XP Professional (when it is not joined to a domain) that makes it possible to switch between users without actually logging off from the computer. Multiple users can share a computer and use it simultaneously, switching back and forth without closing the programs they are running.

We are ready to begin the interactive part of the tutorial.



Tutorial Part 1: **Chromera Manager**



About the Chromera Manager

Chromera Manager is the starting point for your interaction with instruments (devices) in the Chromera environment. More specifically in Chromera Manager...

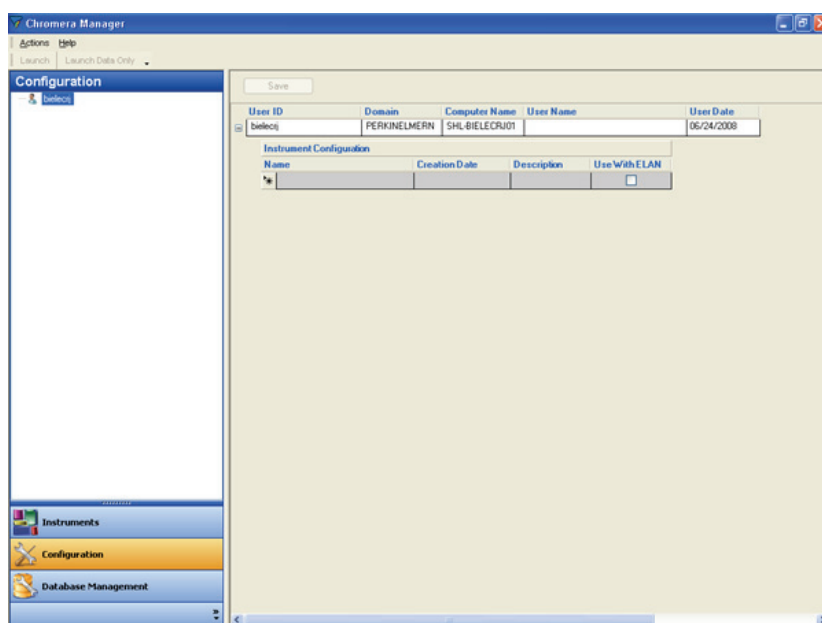
- Instrument Configurations are created here.
- Overall instrument Status can be monitored here.
- Database Management utilities are accessed here.
- Instances of Chromera, either for instrument control or for **Data Only**, are launched from the **Chromera Manager**. An instance of Chromera for an HPLC cannot be started directly from Windows.
- When using Chromera with the **ELAN** ICP-MS, the Chromera software can only be started from the **ELAN** software.

Starting Chromera Manager

- To open Chromera Manager; Click the **Windows Start** button, then click **Chromera Manager**.



The **Chromera Manager** dialog appears:




Create an Instrument

An instrument is actually a collection of devices (such as a pump module, a detector, or an autosampler) which function together as a unit (instrument). Databases are also defined as part of the instrument.

NOTE: Before beginning, make sure the Edgeport is installed and your devices are connected as described in the Chromera Installation Instructions.

To configure an instrument, follow this procedure:

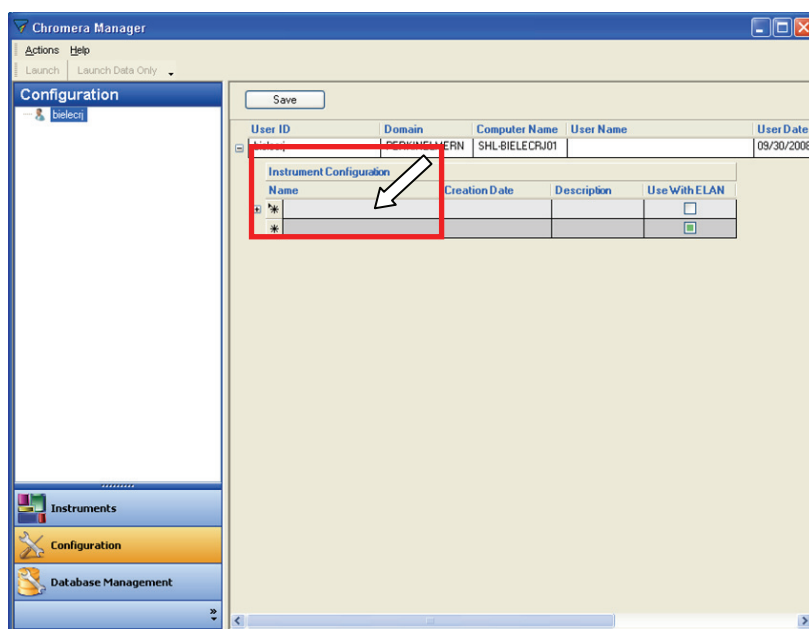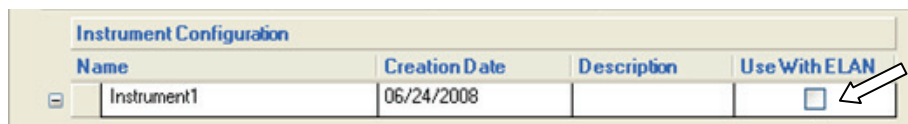
1. Click on the **Configuration** button  to display the initial Configuration screen.




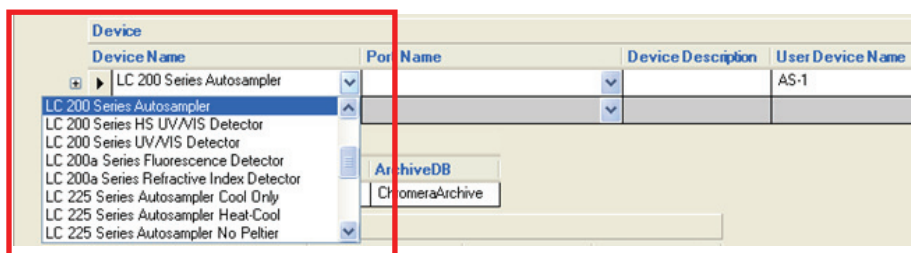
2. Click in the **Name** field and type a name to define the name of the Instrument.

This name can be up to 50 characters and must be unique for this user. The **Creation Date** displayed is the creation date of this Instrument Configuration. You can also type a full **Description** of this Instrument Configuration.

NOTE: Click the **Use With ELAN** checkbox if this Instrument Configuration is used when Chromera is run from the ELAN software.

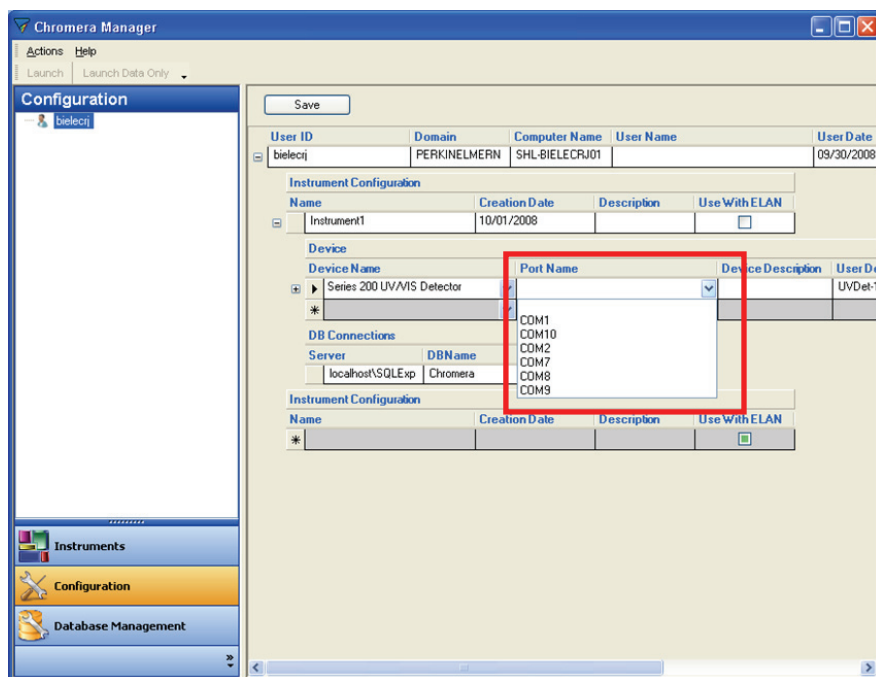


3. Click the plus button  to the left of **Name** to expand the entry fields so that you can enter the **Device Names**.
4. Click in the **Device Name** field and select a Device Name from the drop-down list.




5. Click in the **Port Name** field and select a Port Name from the drop-down list.


All instruments and the Fluorescence detector will have a COM port. For example if the Autosampler is connected to the Computer, it could be COM1. Depending on the Edgeport, COM ports could be labeled from COM1, COM2, etc.



6. Click in the **Device Description** field and type a description of this device.
Initially the **User Device Name** field displays a short name for the device. This is an abbreviation for the device type, with a number added if more than one device of that type is present in the configuration. You can modify this by typing a name (using up to 50 characters) by which the device will be known within the software. Each user device name must be unique within the instrument. This name is used (for example) in the Method Editor.
7. If necessary, click in the **Serial Number** field and type the serial number of the device (using up to 50 characters). This step is an optional.
8. If necessary, click in the **Firmware Version** field and type the version number of the firmware installed in the device (using up to 50 characters). This step is an optional.
9. To add another device, click in the next **Device Name** field and select a Device Name from the drop-down list. Continue to enter the information for the device.
10. Leave the default **DB Connections** settings for Server, DBName, and ArchiveDB. These settings will appear as the Server, DBName, and Archive Database for this instrument.
The default DBName is Chromera. Each user can have their own database. Each instrument can have its own database. Instruments and/or users can share a common database.

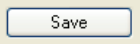12 . Chromera Tutorial for New Users

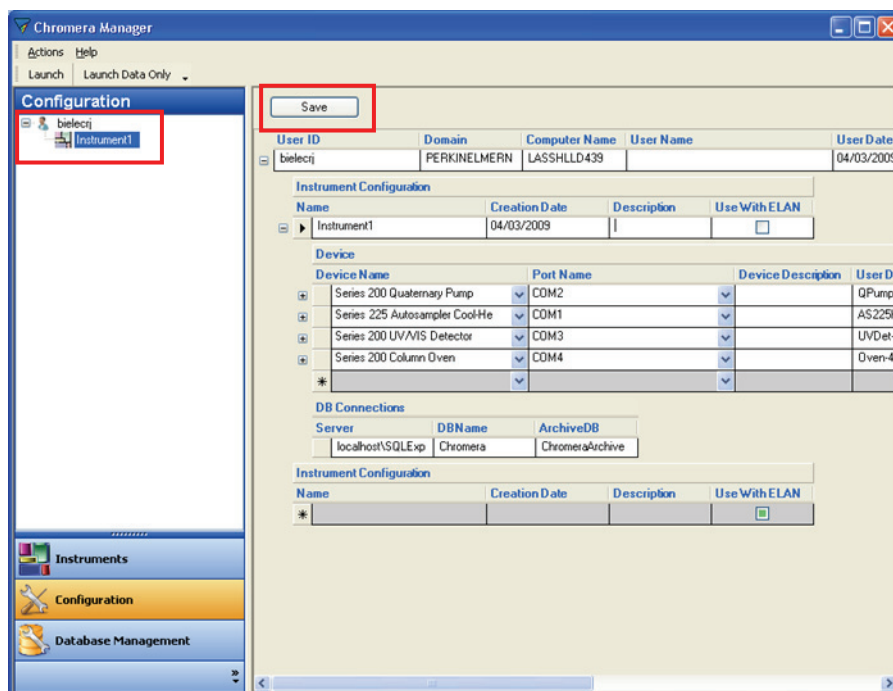11. When done, click the **Save** button  to save this Instrument Configuration. The Instrument is displayed under the user name in the Navigation Pane.



About Database Management

The following general principals apply when using a database.

- Each instrument requires a database.
- You can create the new database for the configuration.
- Each person using Chromera can have their own database. Each instrument can have its own database. Also, instruments and/or users can share a common database.
- Searching in a database is a database function. A database has no files; however, there are database records.
- If the database already exists, a message will be displayed alerting you. Switch to the Delete tab to delete the existing database, if desired.

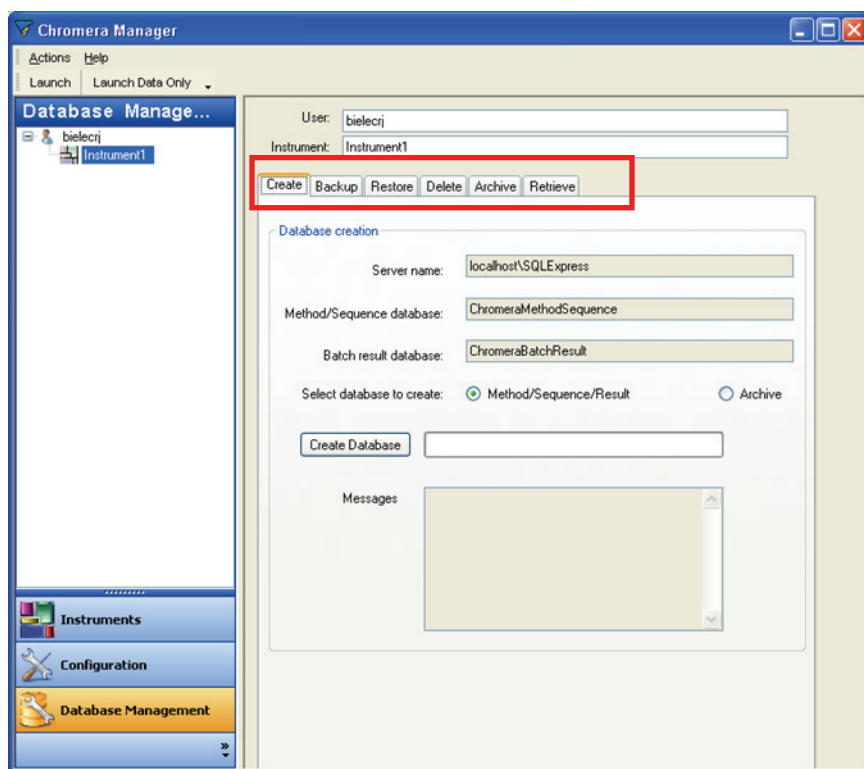
IMPORTANT: Once you delete a database it **CANNOT** be recovered.

Create a Database

Before running Chromera you must create a database. To create a database, follow this procedure:

1. Click on the **Database Management** button  to display the database screens.

The database functions (**Create, Backup, Restore, Delete, Archive, and Retrieve**) are displayed as tabs. Click on the tab to display more detail of your desired function. The first screen that appears is **Create**.



Notice that the **Server name** is filled in. This name corresponds to the default **DB Connection** settings that were made when you configured an Instrument in **Configuration**.

2. Select one of two types of databases you want to create. Select either **Method/Sequence/Result** or **Archive**.

In this example we are creating the **Method/Sequence/Result** database.

3. Click the **Create Database** button.

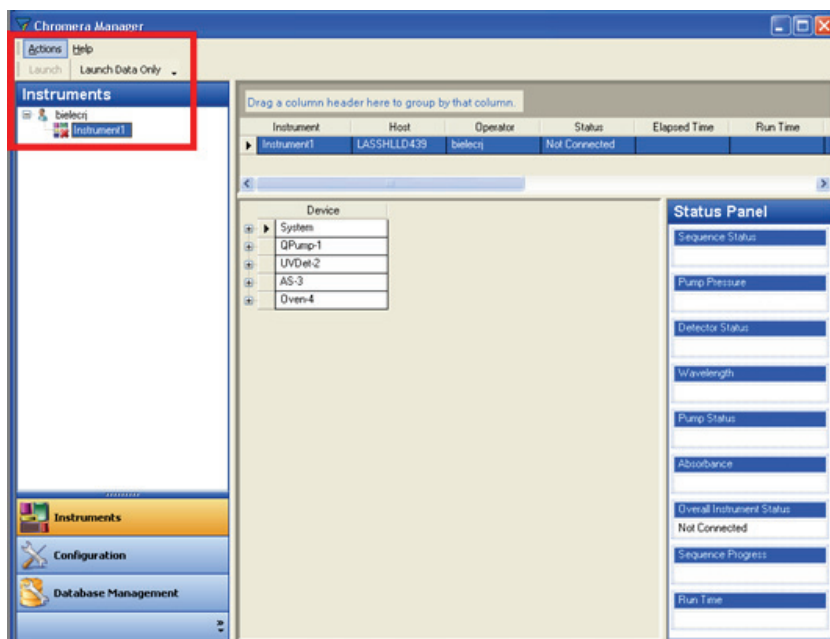
Observe the progress bar to the right of the Create Database button as the database is created. When the database creation is complete a message appears in the **Messages** box indicating this.

Launching Chromera


The Chromera application is different from Chromera Manager. To launch Chromera, follow this procedure:

1. Select the **Instrument** you want to initialize and click the **Launch** button ... or simply double-click the instrument name.

Each instrument needs to be launched in its own session.



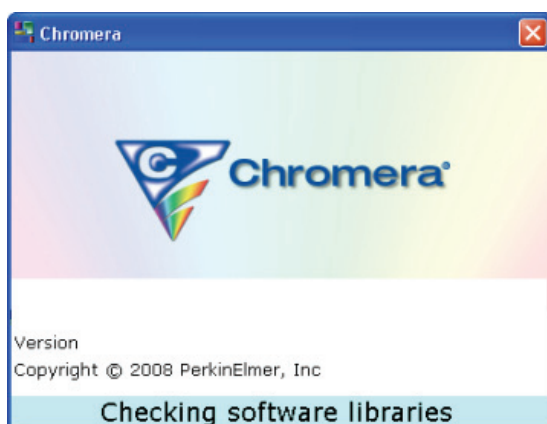
The system does extensive checking of the Chromera installation to ensure that the database is intact, proper program files have been installed, etc.




Any normal errors or warnings encountered during startup are displayed.

2. Once the system is running, instrument errors are displayed in the fade-in notification area, and will be written to the **Error Log**, accessible from the Chromera **Tools** menu.

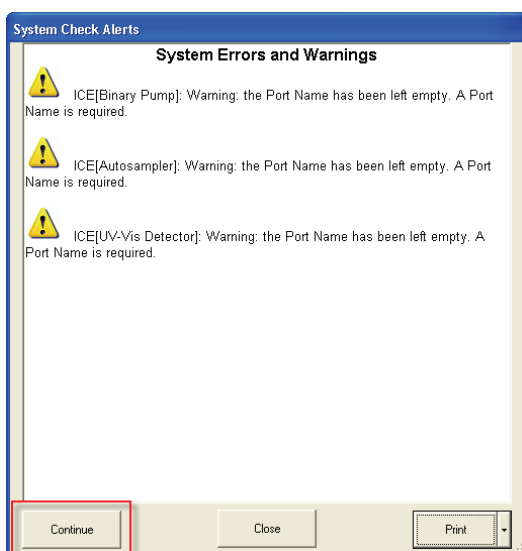
Unexpected errors are trapped at the OS level and recorded in the Windows **Event Log** under **Application Errors**.



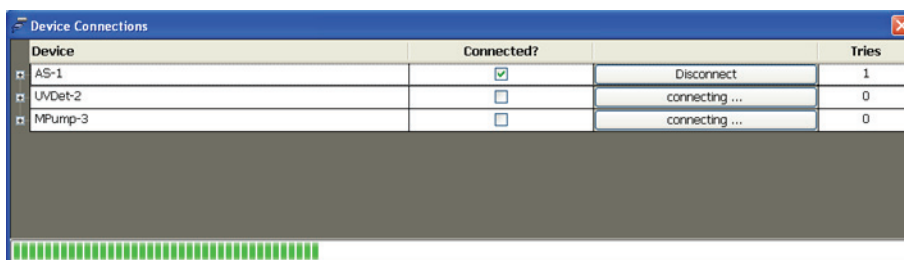
- Once Chromera has been launched, the **Device Connections** dialog displays showing the status of the connections. As the Device tries to connect, **Connecting ...** displays in the button.



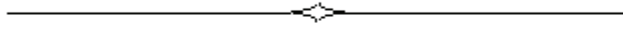If a device does not connect, click the **Connect** button to retry making a connection to the instrument module. You can display the **Device Connections** dialog at any time from the Chromera **Tools** menu.



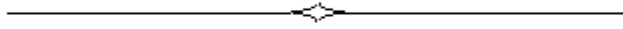Tutorial Part 2: Working with Example Data




About Working with Example Data


This session describes how you can **Import** example data into Chromera and work with the Data Selector. We are using a specific example data file so you can familiarize yourself with the specific features within Chromera.

NOTE: *If you are somewhat familiar with Chromera and you have acquired your own data you may want to use it in place of the Example Data.*

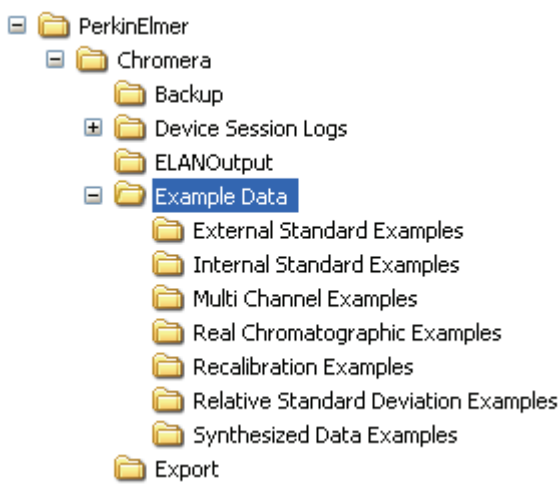1. Upon installation of Chromera, the **Example Data** are installed into in the **Chromera** directory as shown below.



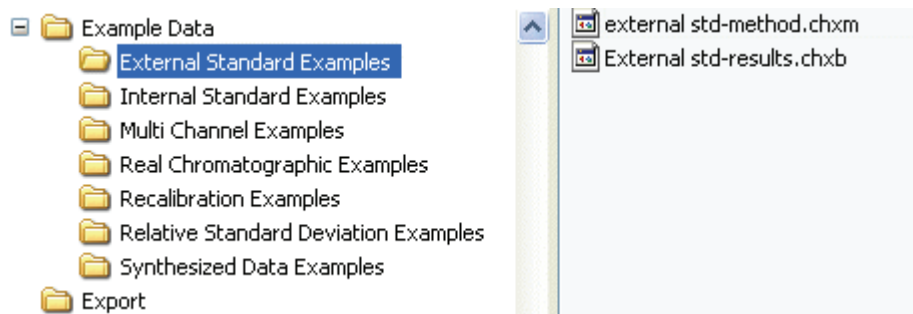
2. Expand the **External Standard Examples** directory to display the two examples.




Importing Chromera Data and Methods

To import the Chromera Results file **External std-method.chxb**

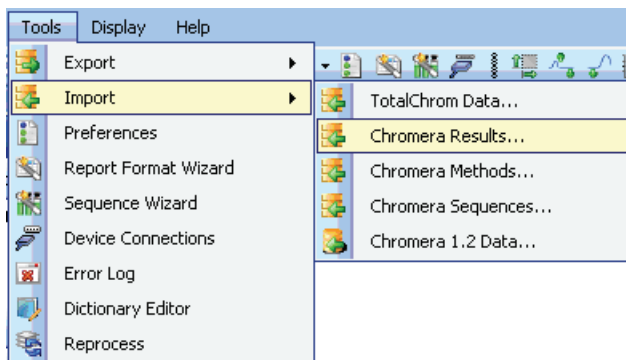1. Select **Import** from the **Tools** menu then select **Chromera Results...**




The **Import Results** dialog appears:

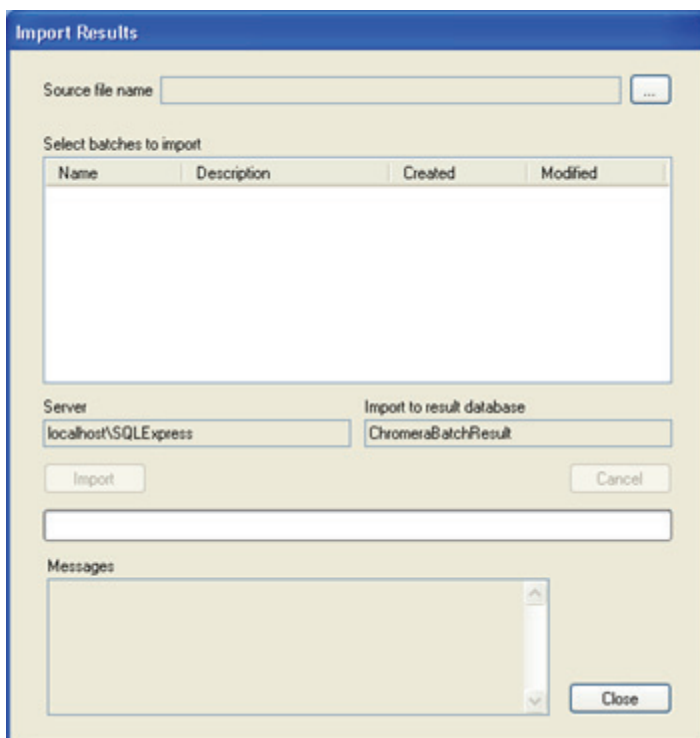


2. Click the browse button  to the right of **Source file name**.

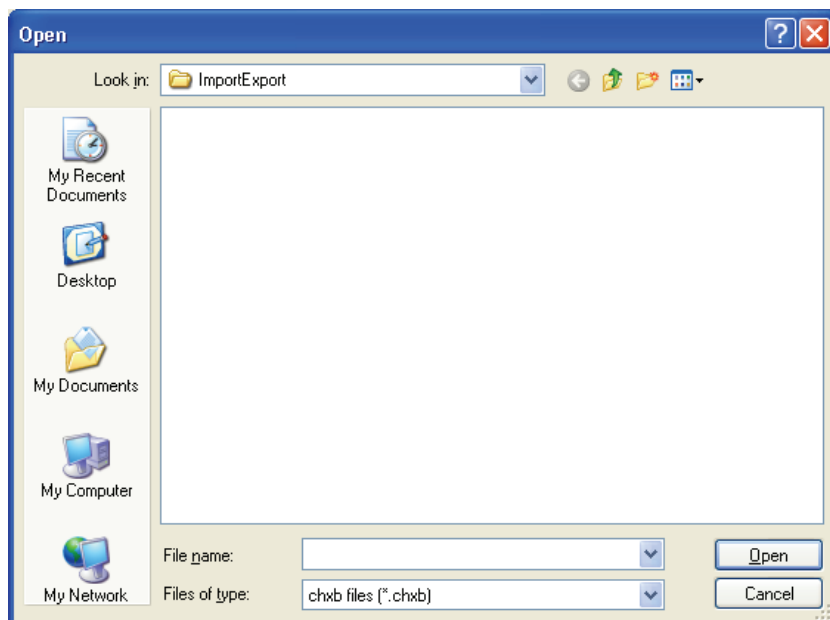
The **Open** dialog appears:



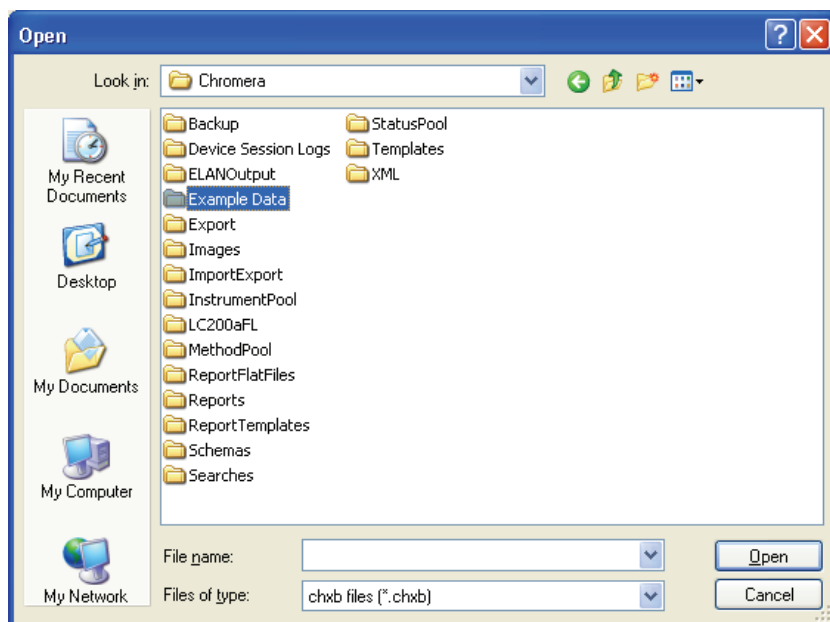
3. Navigate to the **Example Data** directory.



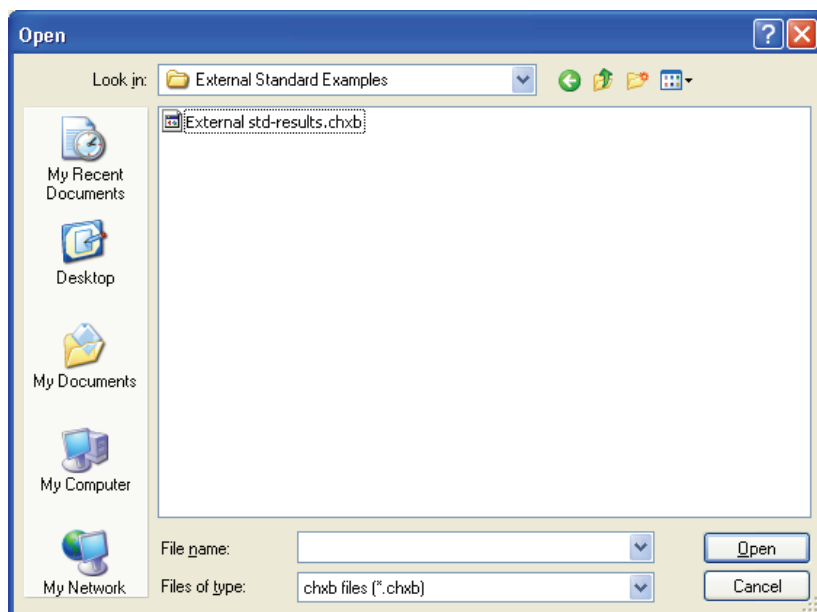4. Double-click on the **Example Data** directory then double-click on the **External Standard Examples** directory.



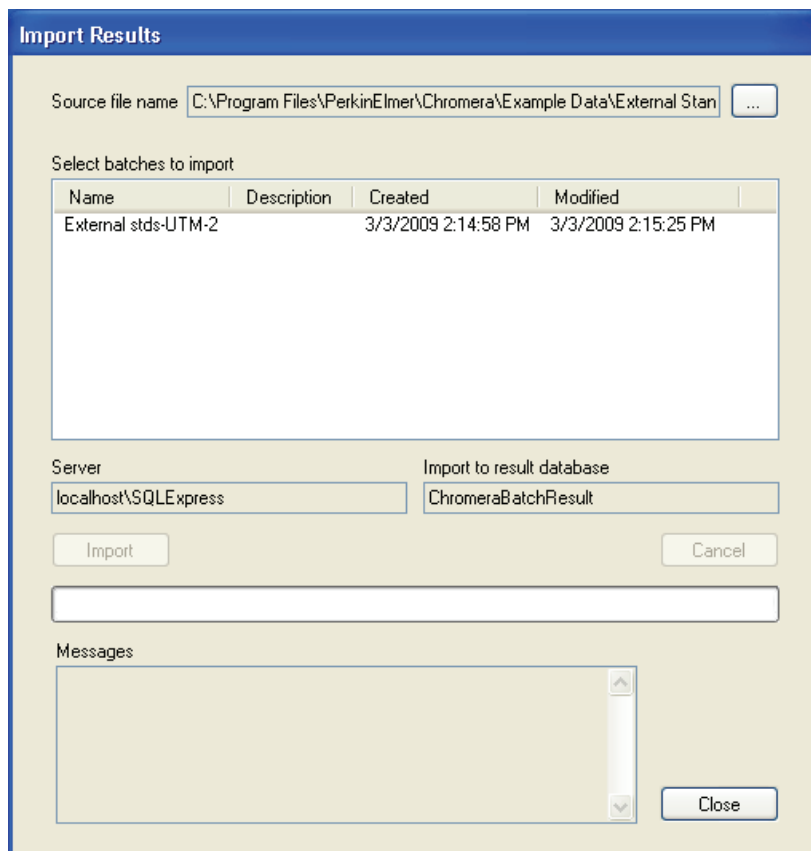5. Select **External std-results.chxb**, then click **Open**.
The file name appears in the **Select batches to import** list.



6. Select **External stds-UTM-2** then click the **Import** button.

The progress bar shows the import progress. Upon completion, the message **Import of results successful** appears in the **Messages** box.

Import Results

Source file name: C:\Program Files\PerkinElmer\Chromera\Example Data\External Stan ...

Select batches to import

Name	Description	Created	Modified
External stds-UTM-2		3/3/2009 2:14:58 PM	3/3/2009 2:15:25 PM

Server: localhost\SQLEXPRESS Import to result database: ChromeraBatchResult

Import Cancel

Progress bar: 100% complete

Messages: Import of results successful.

Close

7. Click the **Close** button.

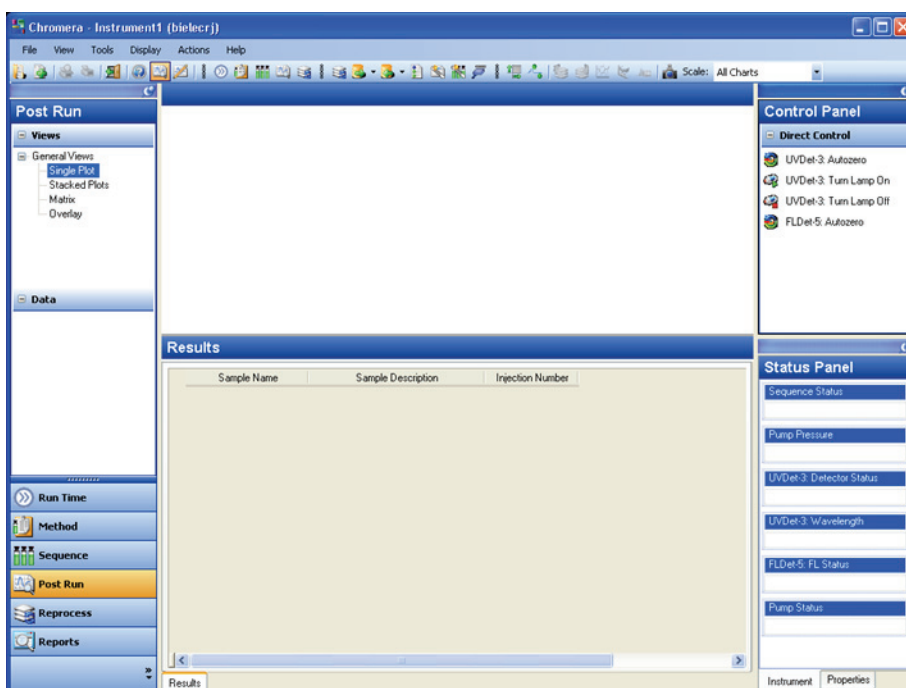
Using the Data Selector

Now that data has been imported into Chromera you can open it using the Data Selector. To open example data file using the Data Selector follow this procedure:


1. Click on the **Post Run** button




The **Post Run** environment opens:



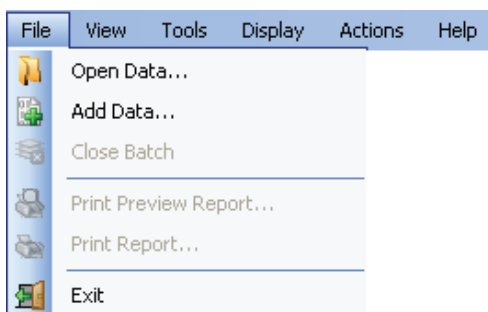2. Select **Open Data...** from the **File** menu.




The **Data Selector** opens.

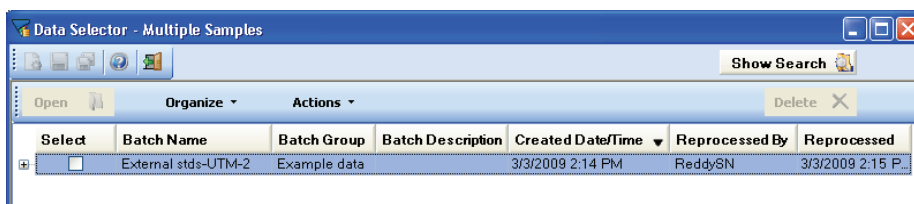

3. Click on the plus sign to expand the batch, then choose one or more samples from a batch.

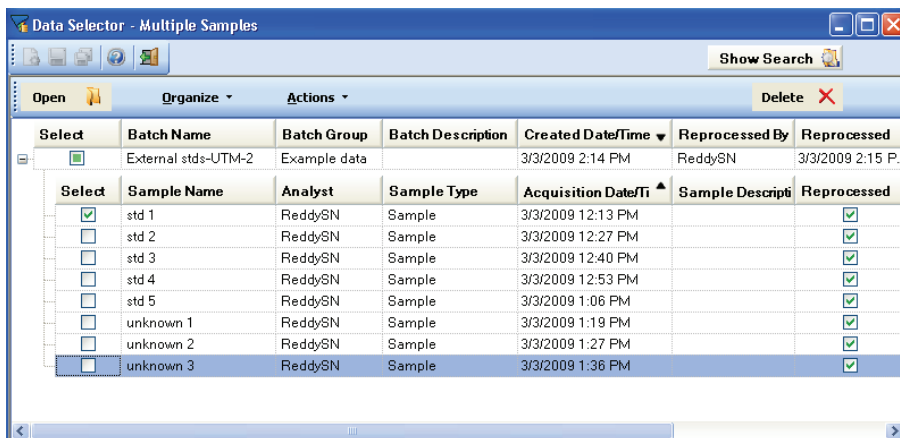The data are stored in Batches, which are, in turn, stored in groups. Use the **Data Selector** to choose data for guiding the method editing. The operation of the **Data Selector** is the same throughout the software.



A **Batch** contains one or more **Samples**. Each Sample can have one or more **Replicate** injections. Each Replicate contains one or more Data **Channels**. For example:

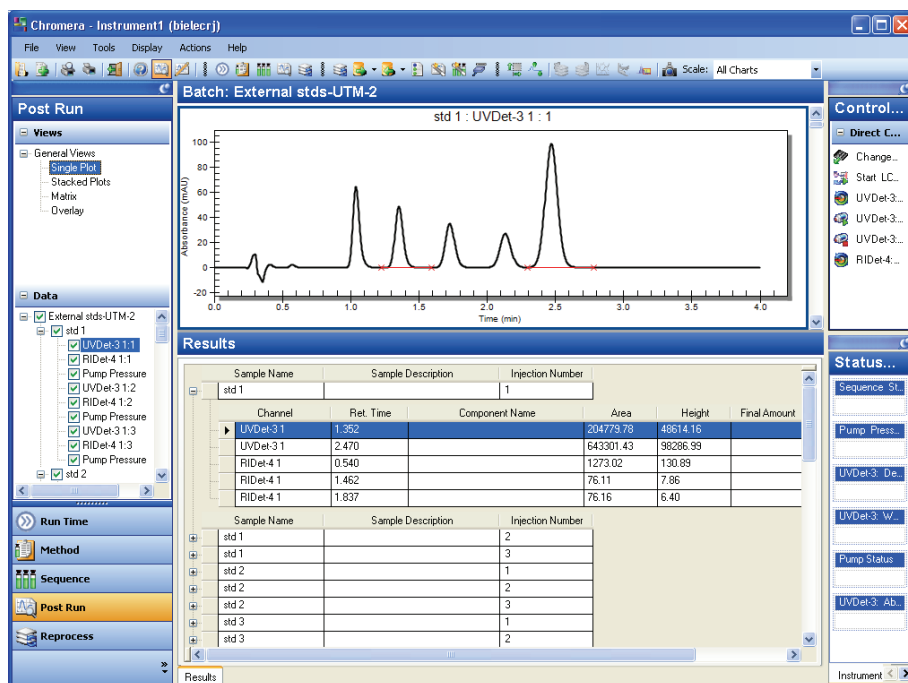
- **Batch**
 - ◆ **Sample**
 - **Replicate**
 - Channel
 - Channel
 - **Replicate**
 - Channel
 - Channel
 - ◆ **Sample**
 - **Replicate**
 - Channel
 - Channel
 -
 - **Replicate**
 - Channel
 - Channel

4. Click in the check box to **Select** the **External stds-UTM-2**, then click the **Open** button.




This selects all of the samples under the batch. The **Post Run** screen displays the **Example Data**.


26 . Chromera Tutorial for New Users





Tutorial Part 3: **Using Chromera**

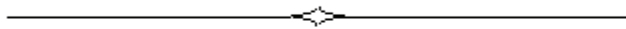



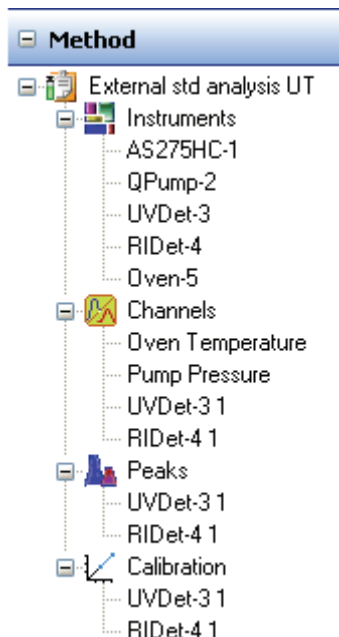

About the Chromera Method

A method contains fundamental parameters to provide unattended control of an instrument during acquisition, calibration, quantitation, and post run. The Method combines Instruments, data Channels, Peak detection, and Calibration parameters into a single entity. These appear as distinct sections in the Method Editor. Clicking the Method button  in the Navigation Pane displays the Method Identification page in the client area (center) of the screen.



Since you are in the **Method** environment, the upper section of the Navigation Pane displays the Method name (along with the general method icon ), whose contents are based on the configuration of the displayed method (this is the configuration of the instrument defined in Chromera Manager). The client area (the large area in the middle of the screen) displays the part of the method corresponding to the tree selection. For example, selecting QPump-2 displays the method parameters for the Quaternary Pump.

The second level: Instruments , Channels , Peaks , and Calibration  represent the main sections of the method and always appear in the Navigation pane. You can reorder these (by right-clicking on a selection using the Move Up/Move Down commands in the pop-up menu) but they will remain at the same level.

The names displayed under each of the second level sections depends on the configuration of the instrument as defined in Chromera Manager. A name appears under Instruments for each device configured. The text displayed is the name you give to the device, as defined in the Chromera Manager's Instrument Configuration.

Each device that generates a data stream that is subject to peak detection will appear under Peak Detection and Calibration. The text displayed is the name of the channel, as defined in the Channels section. Channels are used to control the type of plot on your printed report. Only channels that can be integrated display under Peaks and Calibration. The part of the method displayed in the client area of the screen depends on the node selected in the method tree.

Method editing consists of modifying existing method parameters. This includes a modification to the method name, since in Chromera the name is a method parameter like any other (the method is uniquely identified by a database index that is not visible to a user). After making changes and saving the method, the previous version of the method is replaced by the modified version. Only a single version of a method is saved in the database. Also, a method cannot be saved with the same name as another method in the Group. To create a copy of an existing method but with a different name, leaving the original method unchanged, use the Save As command.

Editing the Method Graphically

In this part of the tutorial you will use some of the basic features of the Graphic Method Editor to:

- Change the plot display
- Adjust peak detection parameters
- Create a component list

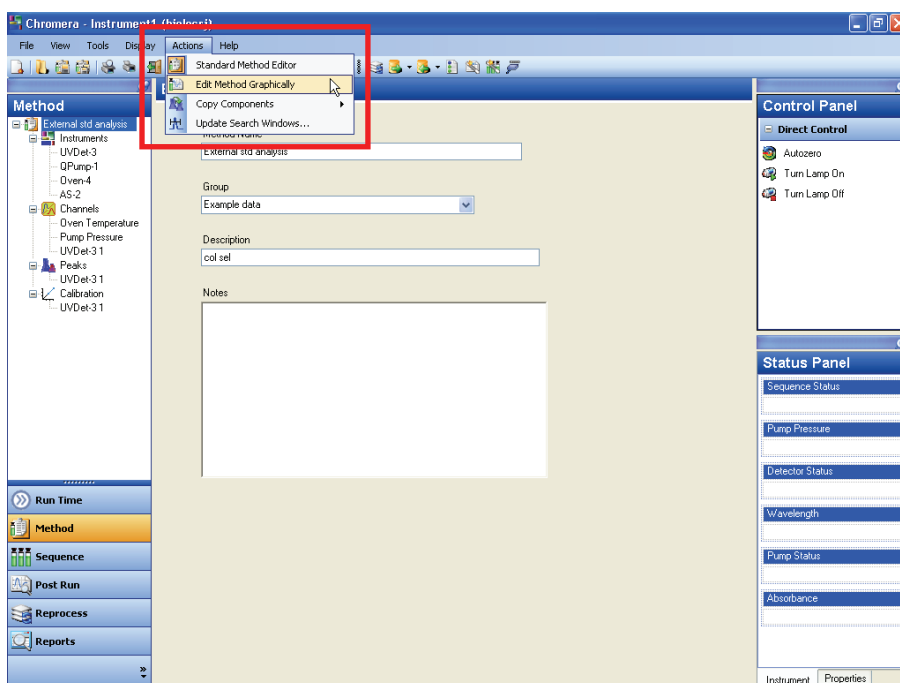
A method contains the parameters that determine how Chromera analyzes raw data. With the Graphic Method Editor, you can adjust those method parameters and see the effects of the change immediately.

Viewing the data in graphic form helps you determine if peaks are being detected correctly, if baselines are being placed properly, and if expected retention times and search windows for components are identifying peaks accurately.


To edit a method in the Graphic Method Editor, select raw data that represents a chromatogram of the analysis that you want to optimize. Chromera processes the raw data with the method currently open and then displays the chromatogram, showing the resulting integration and peak identification.


Selecting Data for Graphic Editing


1. Select **Edit Method Graphically** from the **Actions** menu.




The **Data Selector** appears:

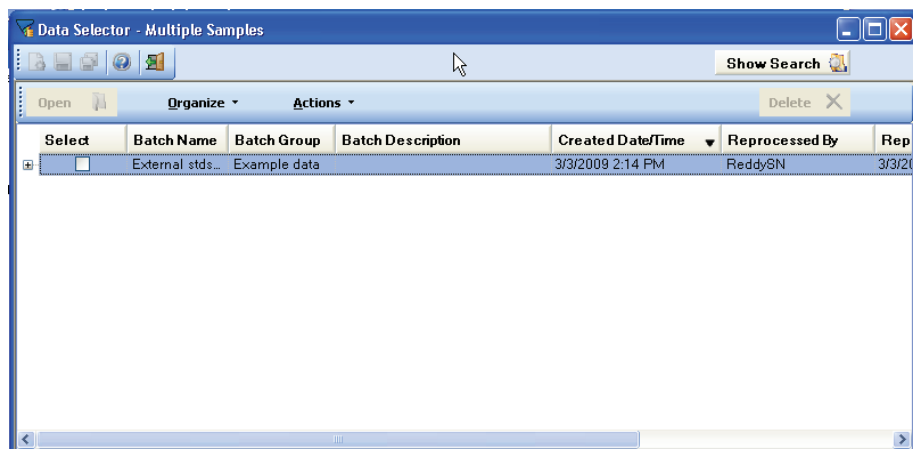

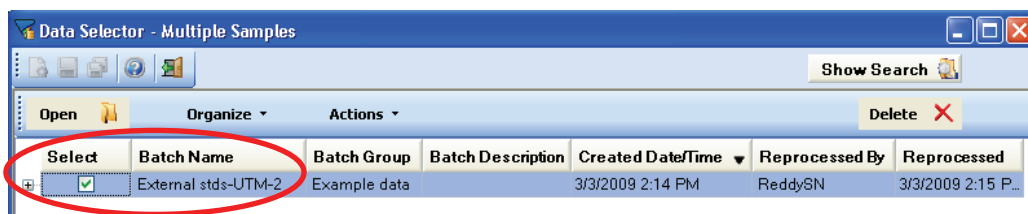- Click in the check box to **Select** the **External stds-UTM-2**, then click the **Open** button. You can also click on the plus sign  to expand the batch, then choose one or more individual samples from a batch, but in this tutorial we are selecting the entire batch **External stds-UTM-2**.



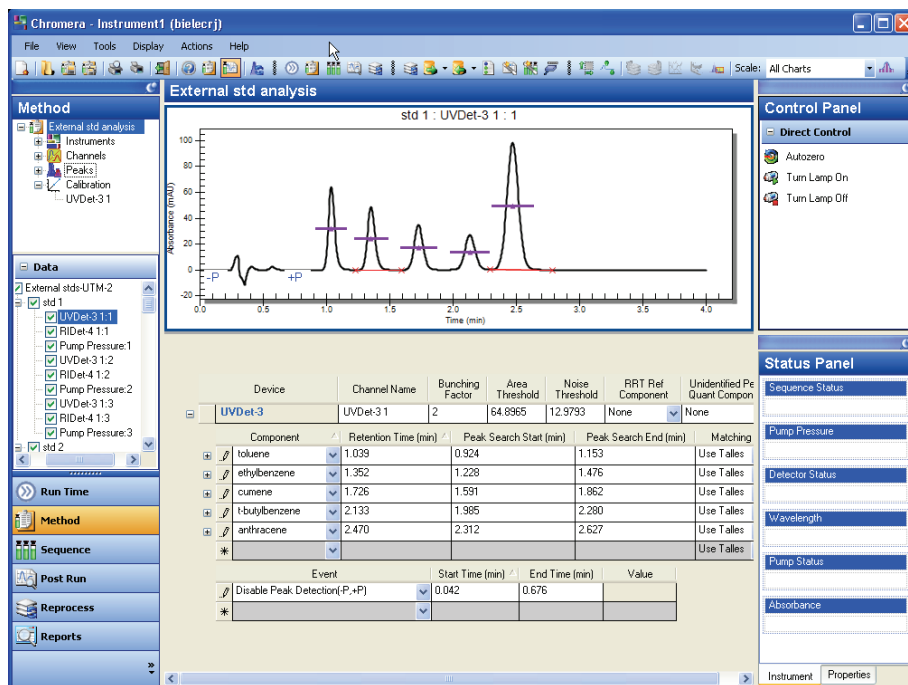The method appears for you to edit graphically.



- Use the currently selected data to optimize the method settings for processing similar samples.

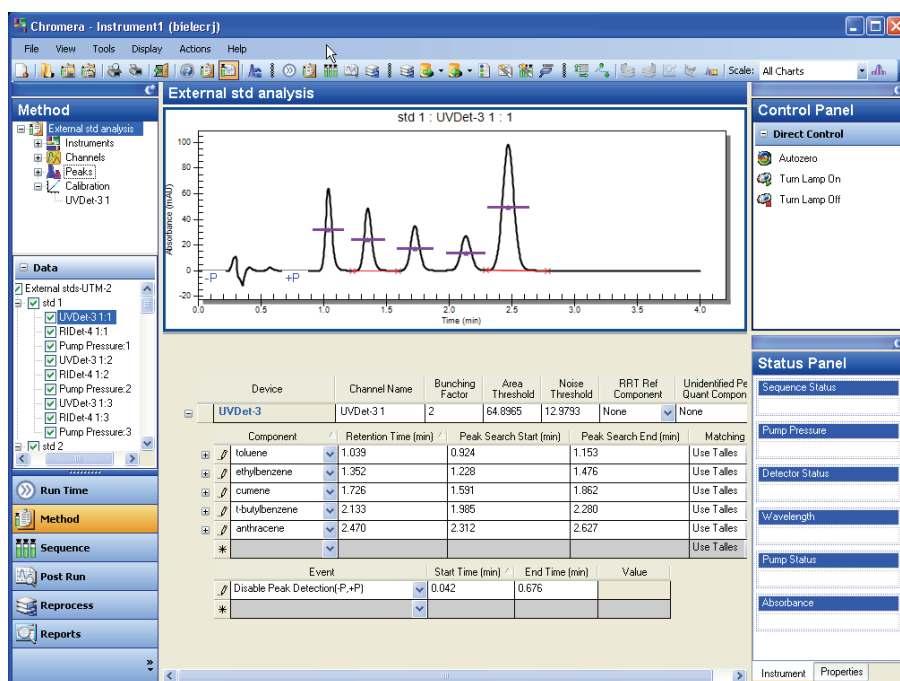
Setting Integration Parameters

Chromera uses the following parameters for optimizing:

- Bunching Factor
- Noise Threshold (Uses decimal digits, not just integers.)
- Area Threshold (Uses decimal digits, not just integers.)

Setting the Bunching Factor (BF)


1. Click in the **Bunching Factor** field.
2. Hold down the **Ctrl** key, click and hold the left mouse button, and drag the pointer across the width of the narrowest peak (this is called Ctrl, click, and drag) to set the **Bunching Factor** for the method.




3. The appropriate **Bunching Factor** value is automatically filled in.

Setting the Noise Threshold (NT) and Area Threshold (AT)

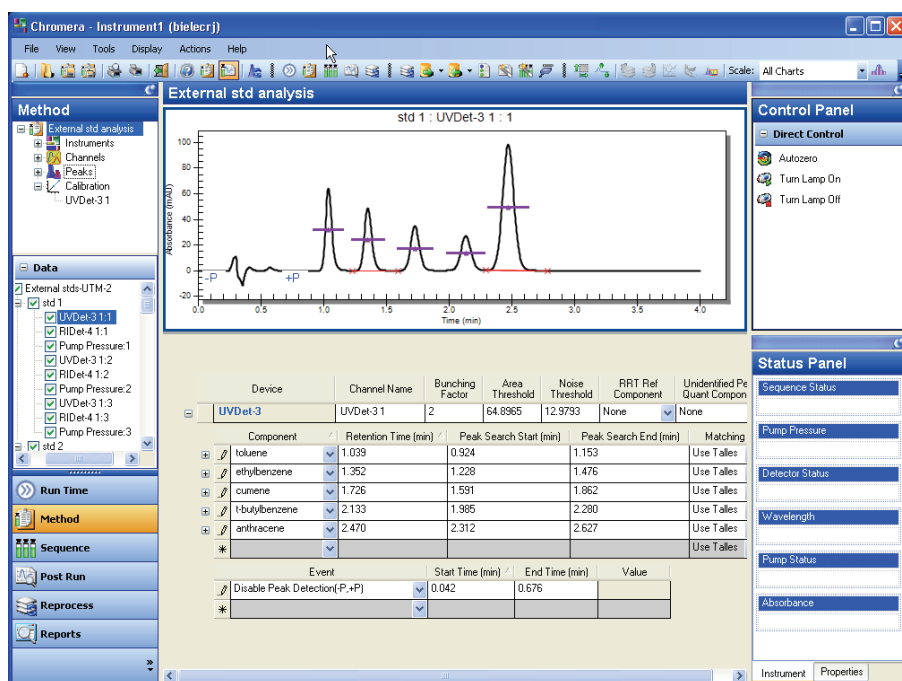
1. Clicking in the **Noise** or **Area Threshold** field.
2. Hold down the **Ctrl** key click and hold the left mouse button and drag the pointer across the region of baseline noise (Ctrl, click, and drag) to set the **Noise** or **Area Threshold** for the method.
3. Select a representative section of baseline to use for calculating the **Noise** and **Area Threshold** values for the chromatogram.
4. The appropriate **NT** and **AT Threshold** values will be filled in automatically.



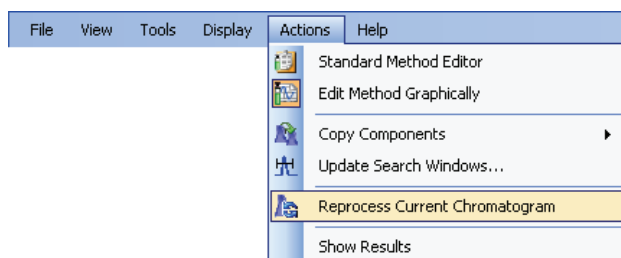5. Select **Reprocess Current Chromatogram** from the **Actions** menu.
Or select the **Reprocess Current Chromatogram** icon on the toolbar.



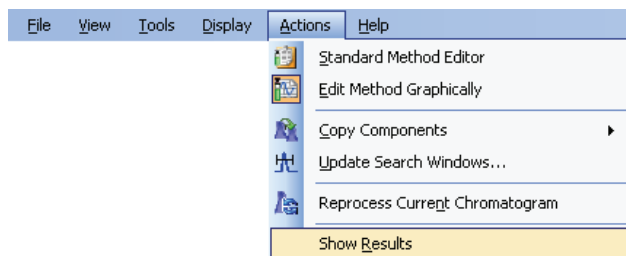6. Then select **Show Results** from the **Actions** menu to view a peak table of the newly reprocessed chromatogram.



GME Results

Sample Name	Sample Description	Injection Number
std 1		1

Channel	Ret. Time	Component Name	Area
UVDet-3 1	1.039	toluene	228525.41
UVDet-3 1	1.352	ethylbenzene	204449.67
UVDet-3 1	1.726	cumene	173214.01
UVDet-3 1	2.133	t-butylbenzene	156255.37
UVDet-3 1	2.470	anthracene	643301.43
RIDet-4 1	0.540		1273.02
RIDet-4 1	1.462		76.11
RIDet-4 1	1.837		76.16

Sample Name	Sample Description	Injection Number
std 1		2

Channel	Ret. Time	Component Name	Area
UVDet-3 1	1.344		435244.20
UVDet-3 1	2.463		645975.07
RIDet-4 1	0.527		1290.40

The current peak information displays in a pop-up window that is always on top and is automatically updated with every new reprocessing.

Entering Component Information

To enter component information, follow this procedure:


1. Select the **Peaks** in the **Method** tree.

Device	Channel Name	Bunching Factor	Area Threshold	Noise Threshold	RRT Ref Component	Unidentified Peak Quant Component
UVDet-3	UVDet-3 1	2	66.7828	13.3566	None	None


Component	Retention Time (min)	Peak Search Start (min)	Peak Search End (min)	Matching
toluene	1.030	0.877	1.198	Use T-Allies
ethylbenzene	1.352	1.228	1.476	Use T-Allies
cumene	1.726	1.591	1.862	Use T-Allies
t-butylbenzene	2.133	1.965	2.280	Use T-Allies
anthracene	2.470	2.312	2.627	Use T-Allies

Event	Start Time (min)	End Time (min)	Value
Disable Peak Detection(P,AP)	0.042	0.676	

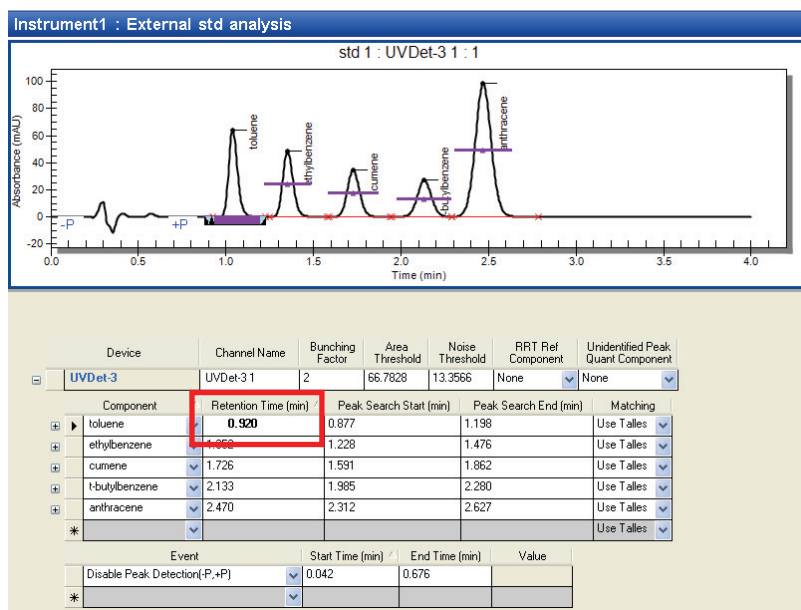
2. In the **Component** column, enter the component names for each peak in the sample.
Auto-Complete can be used to minimize typing when entering names for components already in the dictionary.
3. Click in the **Retention Time** field for a component.



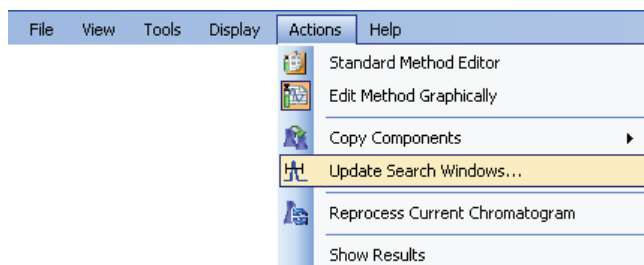
- Click within the desired peak on the chromatogram to enter the expected **Retention Time** for that peak. This value is automatically filled in.
- Select **Update Search Windows** from the **Actions** menu.



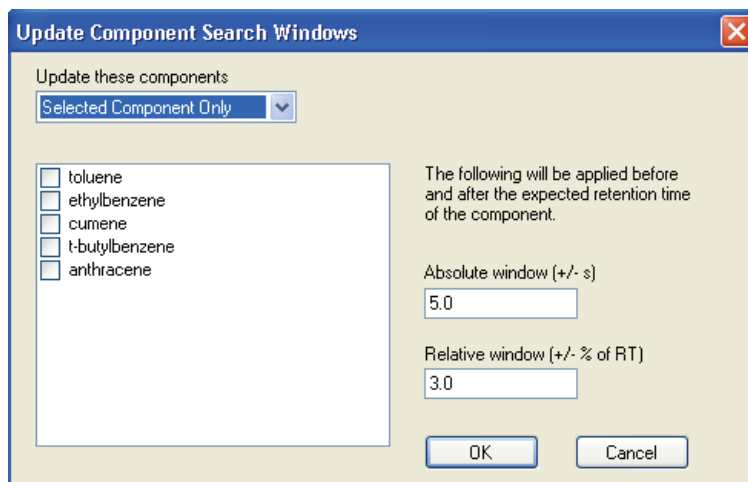A dialog displays for you to change the default settings for one or more components.



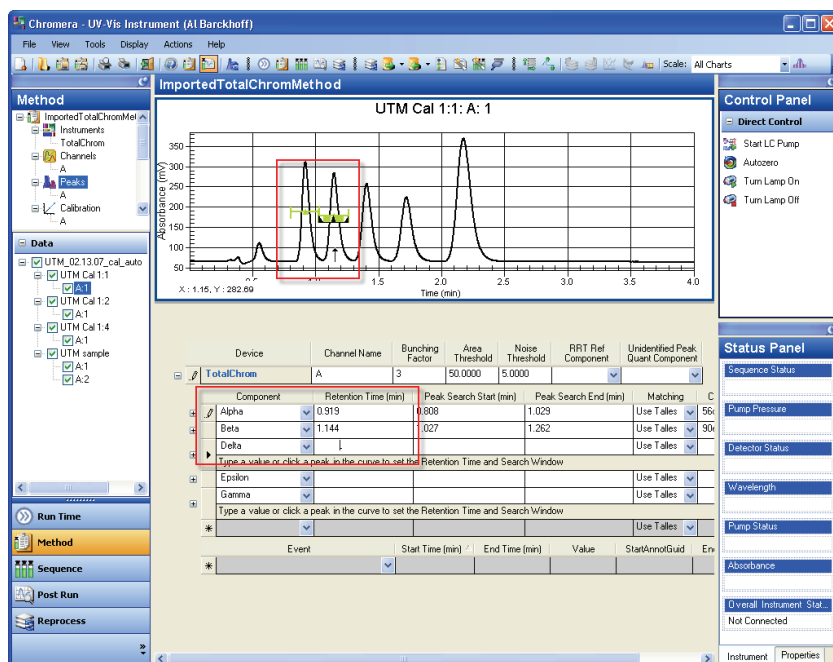
- When done click **OK**.
- Repeat for each component in the method.



Components also can be entered into the method through the **Calibration** section, but retention times and component windows still must be entered in the **Peaks** section.

Entering Baseline Timed Events

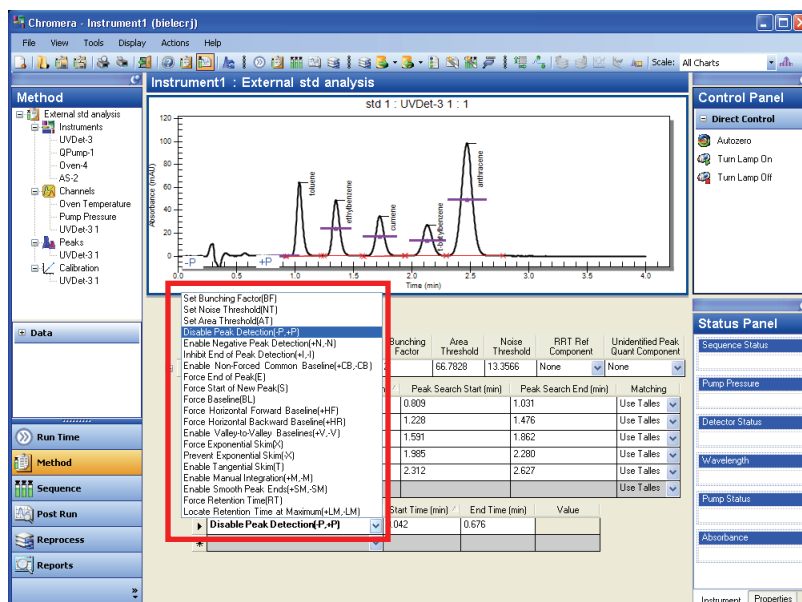
1. Select **Peaks** from the **Method** tree.
2. Select any **Baseline** or **Integration Timed Events** from the list box in the **Events Table**, or by choosing them from the **Toolbar**.



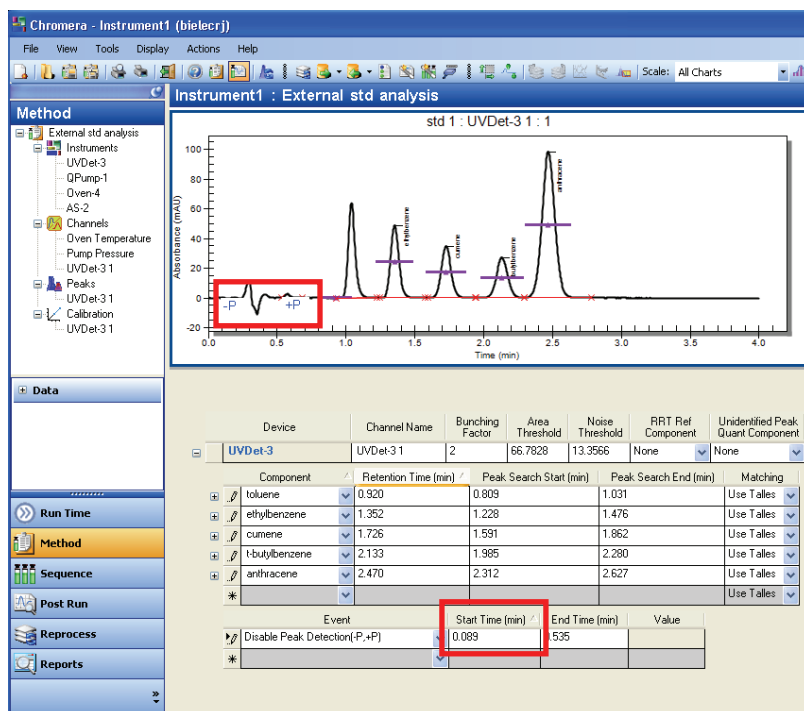The **Integration Timed Events** icons appear on the Toolbar.

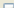


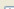




























3. Select the event to be added.

For example, select **Disable Peak Detection**.



4. Click in the **Start Time** field.
5. Ctrl, Click, and drag on the plot to indicate where the selected events should be applied.
 - If it is a single event (no \pm) then click on the chromatogram where the event is to occur.
 - If it is a \pm event, click on the chromatogram where you want the event to start and drag a box to where you want the event to end.

The height of the box does not matter, only the beginning and end matter.

Device		Channel Name	Bunching Factor	Area Threshold	Noise Threshold	RRT Ref Component	Unidentified Peak Quant Component
	UVDet-3	UVDet-3.1	2	66.7828	13.3566	None 	None 
Component		Retention Time (min)	Peak Search Start (min)		Peak Search End (min)		Matching
	 toluene 	0.920	0.809		1.031		Use Tallies 
	 ethylbenzene 	1.352	1.228		1.476		Use Tallies 
	 cumene 	1.726	1.591		1.862		Use Tallies 
	 t-butylbenzene 	2.133	1.985		2.280		Use Tallies 
	 anthracene 	2.470	2.312		2.627		Use Tallies 
	 						Use Tallies 
Event		Start Time (min)	End Time (min)		Value		
	Disable Peak Detection(-P,+P) 	0.089	0.535				
	 						

The **Start Time** and **End Time** values for the event are automatically filled in. The actual position of the event is at the left edge of the event annotation marker on the plot. Timed events can be repositioned on screen by clicking-and-dragging the event markers on the plot.

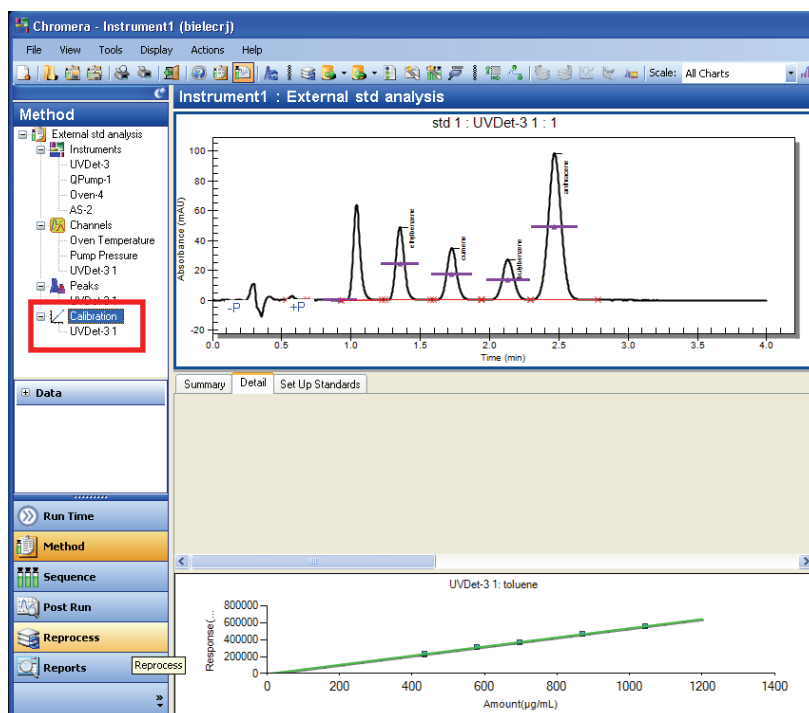
Specifying Component Calibration Information

Calibration information and calibration levels can be entered either from the **Details** tab or from the **Set Up Standards** tab. Calibration information in the method can be viewed as a component-by-component **Summary** of the current calibration curves.

The **Details** view shows the current calibration settings and replicate data, if any, on a component-by-component basis. Using **Set Up Standards** provides a convenient way to enter calibration standard information on a standard-by-standard basis and to visualize each calibration standard level.

1. Select **Calibration** in the **Method** tree.


Calibration opens with the **Detail** tab selected.



The component names should be automatically transferred into the **Calibration** table from the **Peaks** list.

2. Click on the device under **Calibration**. In this example the device is, **UVDet-3 1**

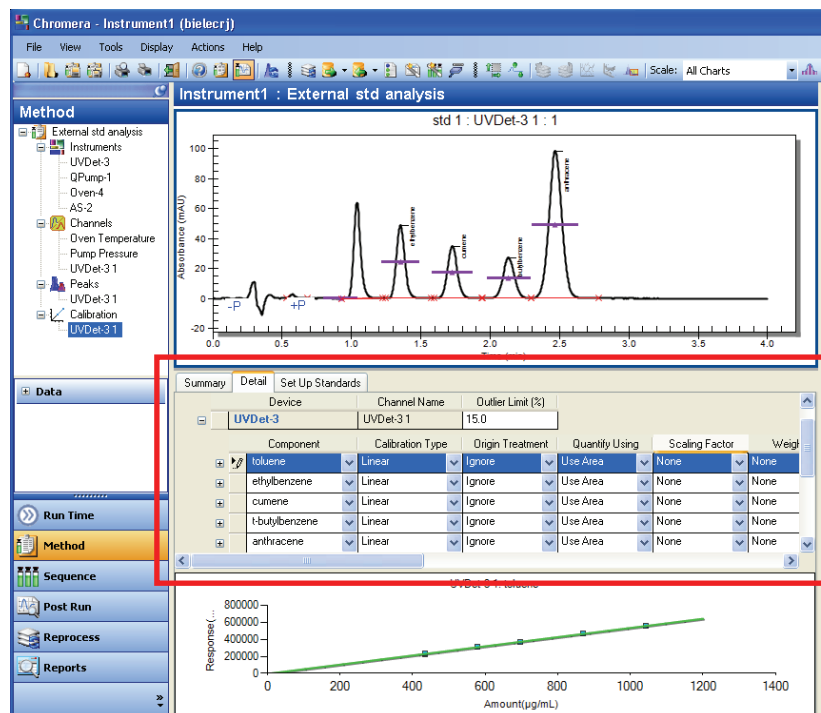The Details appear.




- Fill in the component and calibration information as desired.
- Expand the component record in the **Details** view (in this example, ethylbenzene is expanded) to display the current calibration average values, as well as any individual replicate values that constitute those current averages.

The **Details** view is a convenient way to review the current state of calibration for the components in the method.




5. Select the **Set Up Standards** tab to identify the calibration standards with their constituent components and concentrations.

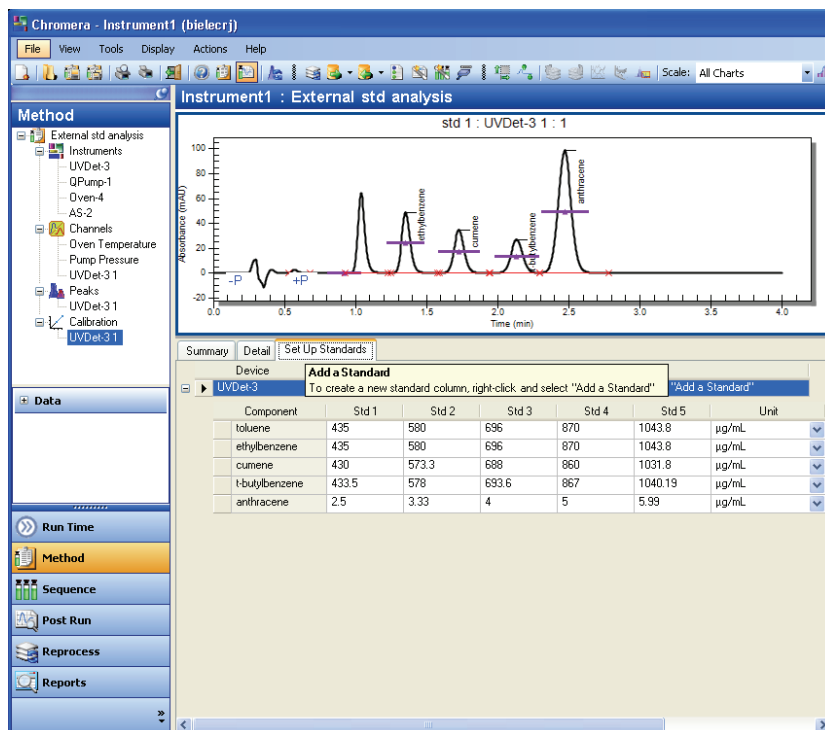NOTE: **Set Up Standards** is the more convenient way for entering Calibration Standards.



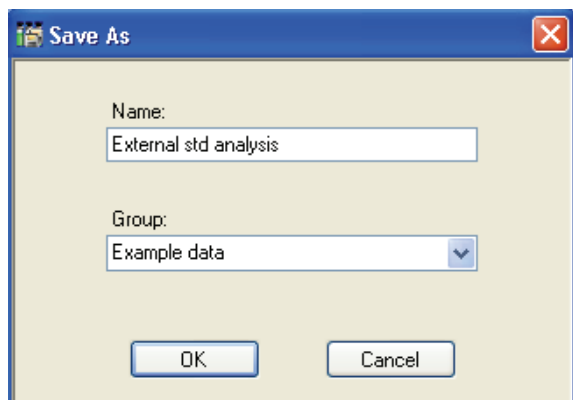6. Enter any number of calibration levels by using a right mouse click through the **Add a Standard** command on the context menu.
Context menu selections provide convenient editing functions to fill in the Standards table.
7. Change other calibration settings for the components in the method as required.
8. Save the changes to the current method or select **Save Method As** to create a new method with the current settings.




9. Give a new name to the method and assign it to a **Group**.

Using the Report Format Wizard


Use the **Report Format Wizard** tool to create single sample or summary reports, method printouts, calibration printouts, or sequence printouts. The wizard is the only mechanism to create report templates in Chromera. There is no Report Editor as there is a Sequence Editor with Sequence Wizard.

To create a new report format using the Report Format Wizard:

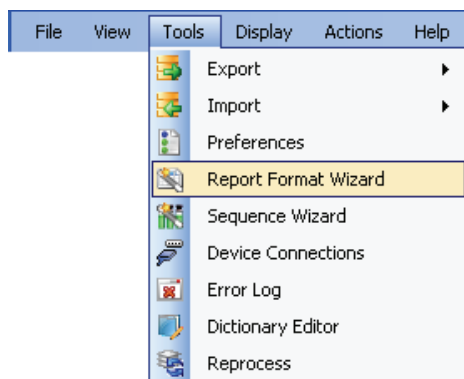1. Select the **Report Format Wizard** from the **Tools** menu.




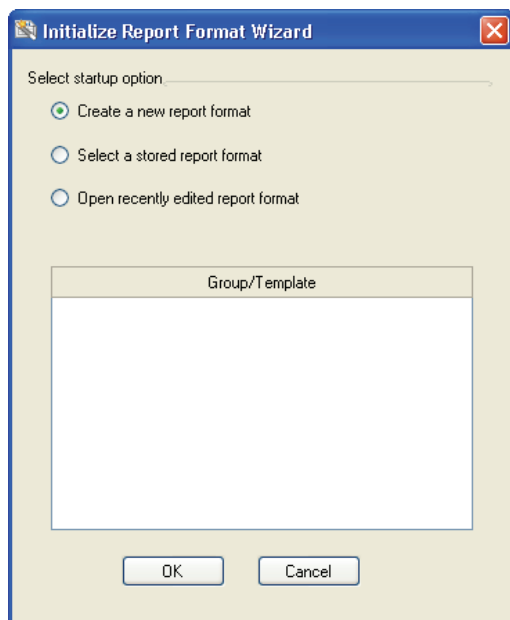The **Initialize Report Format Wizard** screen appears.




2. Under **Select startup option** select **Create a new report format** then click **OK**.

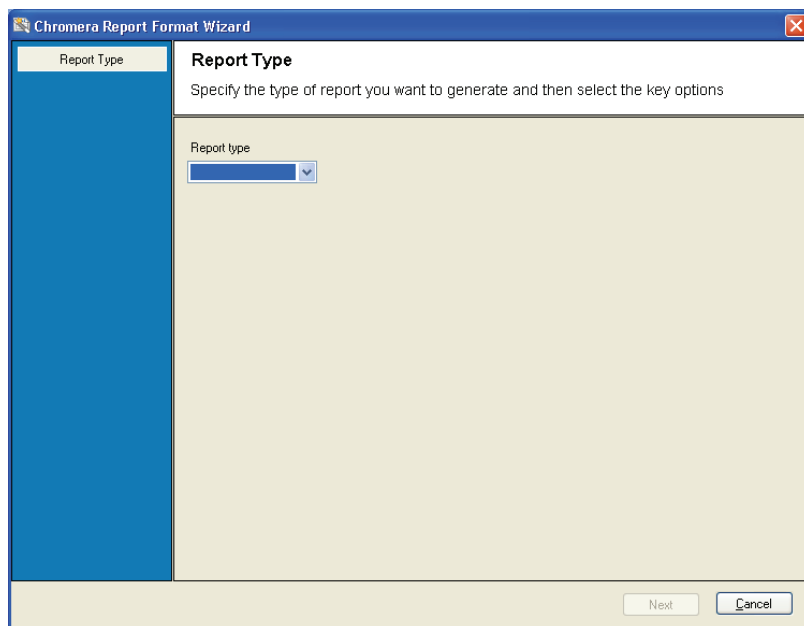The **Report Type** screen appears:




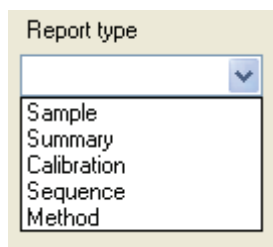3. Select a **Report Type** from the drop-down list.




You can create **Sample**, **Summary**, **Calibration**, **Sequence**, or **Method** reports.

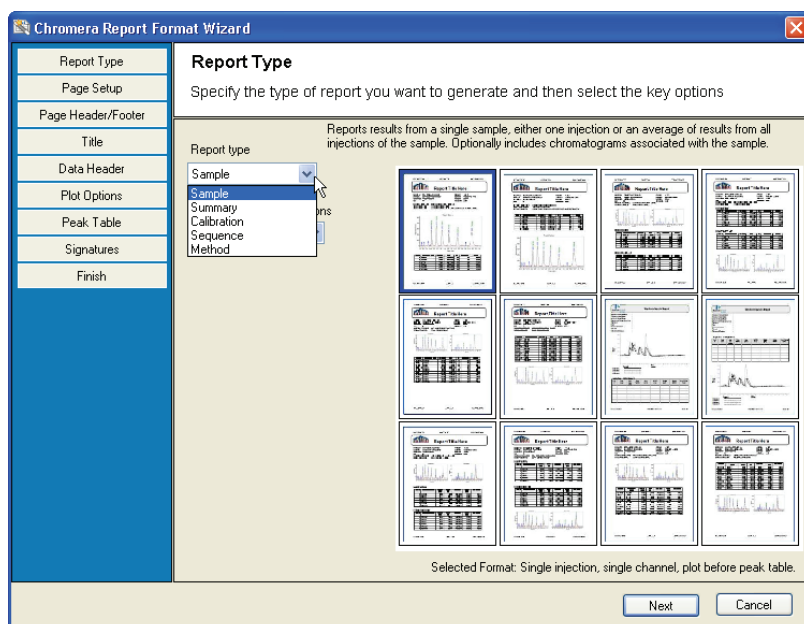4. Click **Sample** to display a gallery of report styles showing the possible layouts for the report type selected.




The gallery styles available are narrowed-down based upon the selections for the type of data to be reported. A **Sample Report** can include the data and plot from one or more replicate injections and one or more data channels.

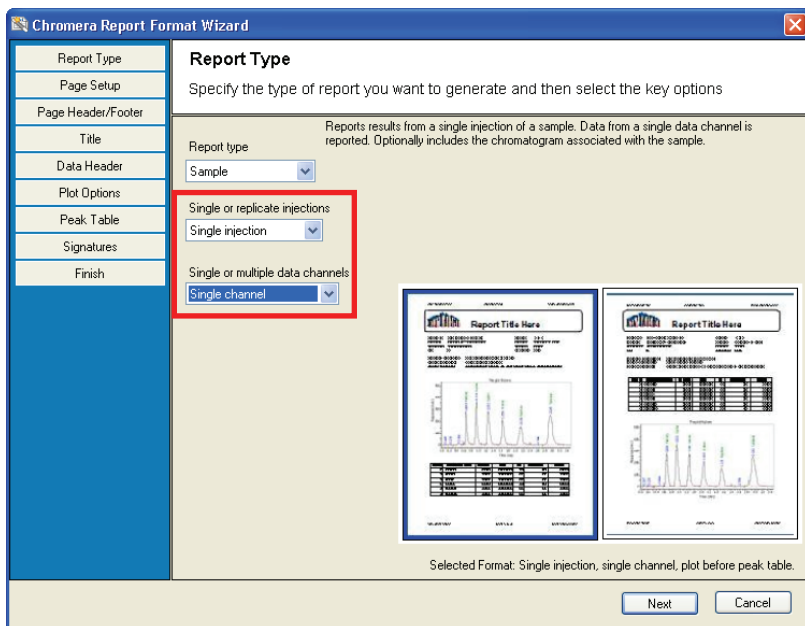- From the **Single or replicate injections** drop-down list select **Single Injection**.
- Then from the **Single or multiple data channels** drop-down list select **Single Channel**.
- Click on a **Report Format** from the two displayed.



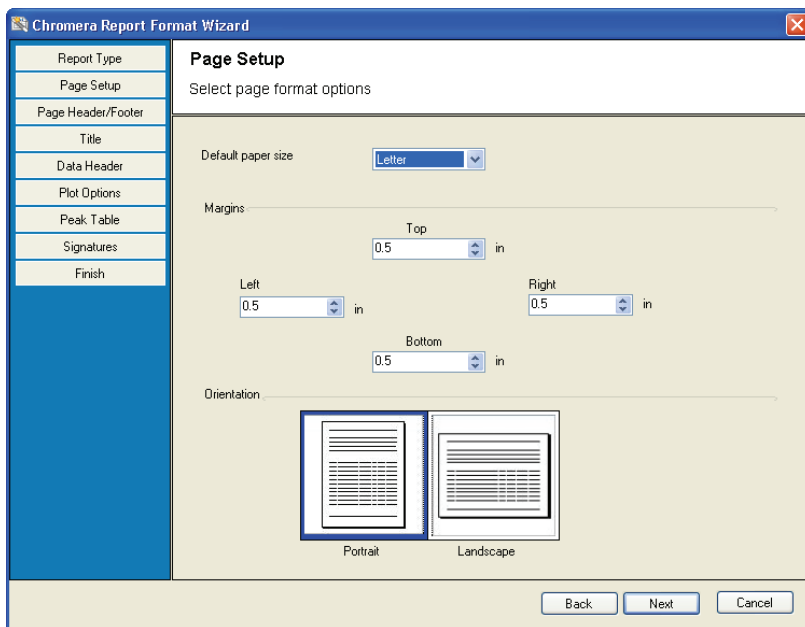- Click **Next** to display the **Page Setup** screen.



- Select the **paper size**, **margin** settings, and page **Orientation**. The units will change based upon the paper size selected.

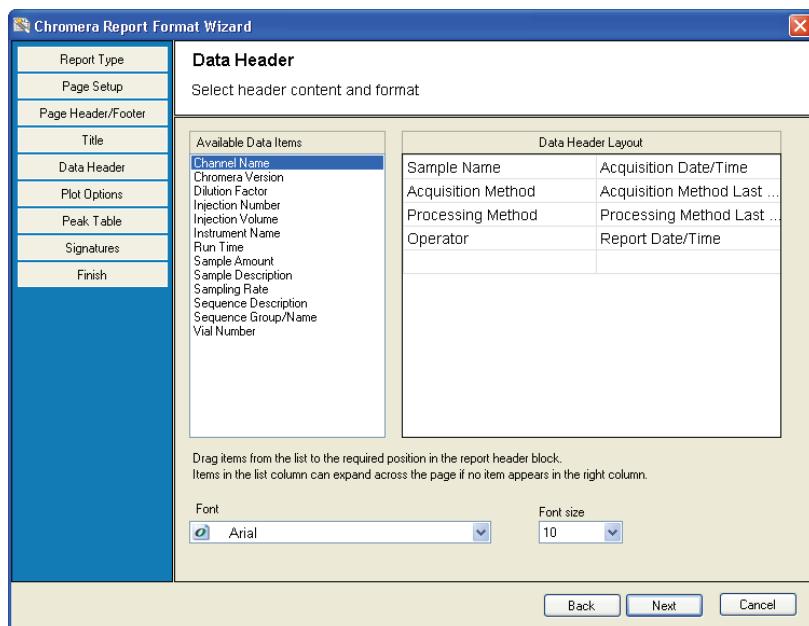
10. Click **Next** to display the **Page Header/Footer** screen.

Each page can have a 3-section **Header** at the top and a 3-section **Footer** at the bottom.

11. From each drop-down list select the information to be displayed in the left, right, and center sections of the **Header** and **Footer**, if desired.
12. Select or deselect the check box if the **Header** and/or **Footer** should be displayed on the first page.
13. Click **Next** to display the report **Title** screen.


14. Select if you want a report **Title** printed on every page or just the first page.

15. You can add a **Logo** image to the report. Any valid graphic file type can be selected. The image can either be embedded within the report template, or loaded from disk and rendered each time the report is printed. You can also add an optional **Border** drawn around the **Title**.
16. Click **Next** to display the **Data Header** screen.



The screenshot shows the 'Data Header' screen of the 'Chromera Report Format Wizard'. The left sidebar contains a list of steps: Report Type, Page Setup, Page Header/Footer, Title, Data Header (selected), Plot Options, Peak Table, Signatures, and Finish. The main area is titled 'Data Header' with the instruction 'Select header content and format'. It features two panes: 'Available Data Items' on the left and 'Data Header Layout' on the right. The 'Available Data Items' list includes Channel Name, Chromera Version, Dilution Factor, Injection Number, Injection Volume, Instrument Name, Run Time, Sample Amount, Sample Description, Sampling Rate, Sequence Description, Sequence Group/Name, and Vial Number. The 'Data Header Layout' pane shows a table with columns for Sample Name, Acquisition Date/Time, Acquisition Method, Acquisition Method Last ..., Processing Method, Processing Method Last ..., Operator, and Report Date/Time. Below these panes, there are 'Font' and 'Font size' dropdown menus, both set to 'Arial' and '10' respectively. At the bottom are 'Back', 'Next', and 'Cancel' buttons.

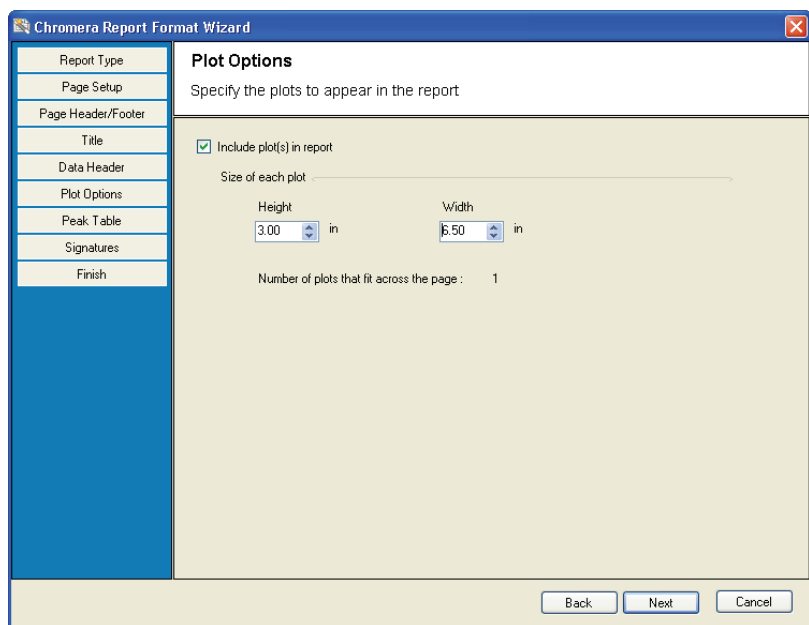17. Select items from the **Available Data Items** list. Position the items in the **Data Header Layout** box by click-and-drag. The Data Header can be constructed from available sample information.
18. Click **Next** to display the **Plot Options** screen.




The screenshot shows the 'Plot Options' screen of the 'Chromera Report Format Wizard'. The left sidebar is the same as in the previous screen, with 'Plot Options' now selected. The main area is titled 'Plot Options' with the instruction 'Specify the plots to appear in the report'. It contains a checkbox labeled 'Include plot(s) in report' which is checked. Below this is a 'Size of each plot' section with 'Height' and 'Width' dropdown menus, both set to '3.00' and '6.50' inches respectively. Below these is a label 'Number of plots that fit across the page : 1'. At the bottom are 'Back', 'Next', and 'Cancel' buttons.

Since a report can consist of one or more sample replicates and/or one or more data channels, then one or more **Plots** can be displayed as well. The **Plot Options** entry will determine how many can be drawn across a page.

19. Select if you want to include the plots in the report and enter the size of each plot.

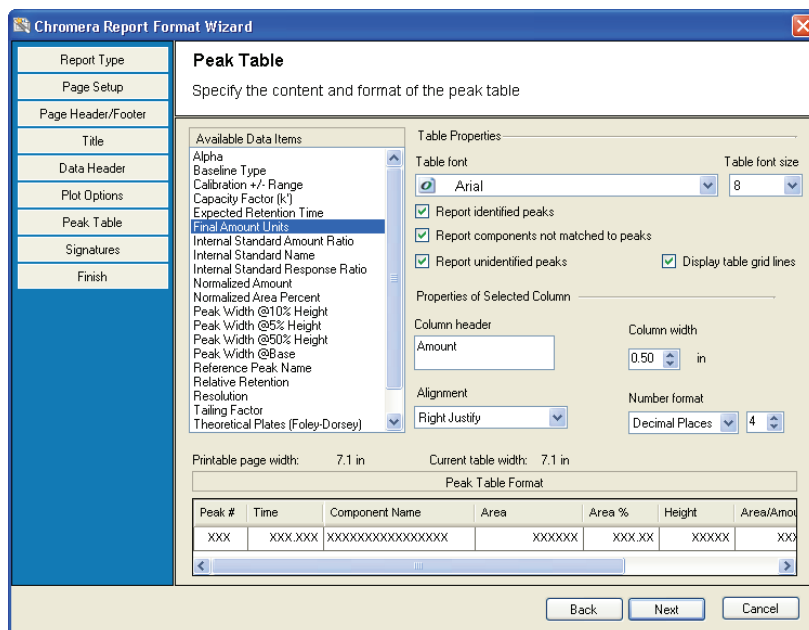20. Click **Next** to display the **Peak Table** screen.



The **Peak Table** screen in the Chromera Report Format Wizard allows users to specify the content and format of the peak table. It includes a list of available data items, table properties, and a preview of the table format.

Available Data Items:

- Alpha
- Baseline Type
- Calibration +/- Range
- Capacity Factor (k')
- Expected Retention Time
- Final Amount Units
- Internal Standard Amount Ratio
- Internal Standard Name
- Internal Standard Response Ratio
- Normalized Amount
- Normalized Area Percent
- Peak Width @10% Height
- Peak Width @50% Height
- Peak Width @90% Height
- Peak Width @Base
- Reference Peak Name
- Relative Retention
- Resolution
- Tailing Factor
- Theoretical Plates (Foley-Dorsey)

Table Properties:

- Table font: Arial, Table font size: 8
- ☒ Report identified peaks
- ☒ Report components not matched to peaks
- ☒ Report unidentified peaks
- ☒ Display table grid lines

Properties of Selected Column:

- Column header: Amount
- Column width: 0.50 in
- Alignment: Right Justify
- Number format: Decimal Places, 4

Printable page width: 7.1 in, Current table width: 7.1 in

Peak Table Format

Peak #	Time	Component Name	Area	Area %	Height	Area/Amount
XXX	XXX.XXX	XXXXXXXXXXXXXXXXXX	XXXXXX	XXX.XX	XXXXX	XX

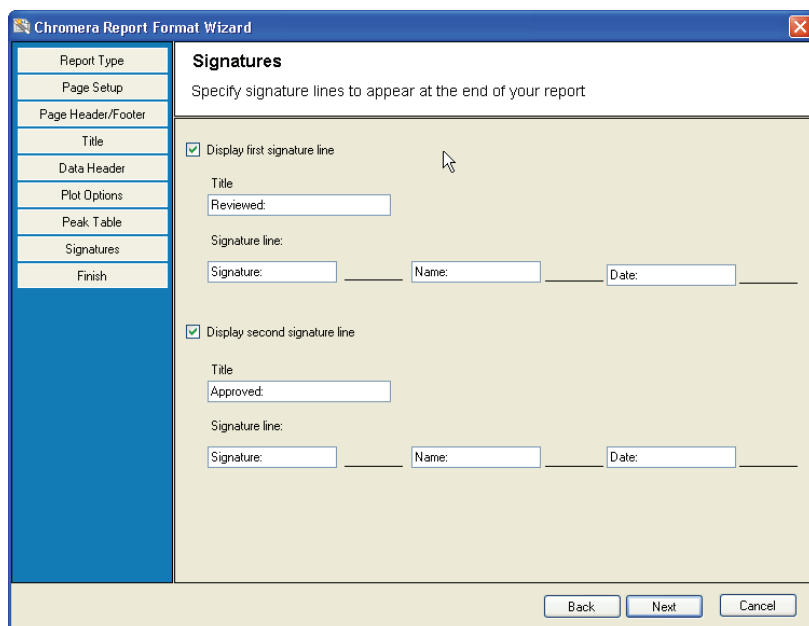
Buttons: Back, Next, Cancel

The **Peak Table** is where the results for each injection will be displayed.

21. Select the desired **Data Items** for reporting.

The selection of items for reporting will change based upon the type of report to be printed. The **Available Data Items** include system suitable parameters that can be calculated without required user input.

22. Select to format the Numeric values with fixed decimal digits or as significant figures.
23. Click **Next** to display the **Signatures** screen.



The **Signatures** screen in the Chromera Report Format Wizard allows users to specify signature lines to appear at the end of their report. It includes checkboxes for displaying the first and second signature lines, and fields for Title, Signature, Name, and Date.

Signatures

Specify signature lines to appear at the end of your report

☒ Display first signature line

Title: Reviewed:

Signature line:

Signature: _____ Name: _____ Date: _____

☒ Display second signature line

Title: Approved:

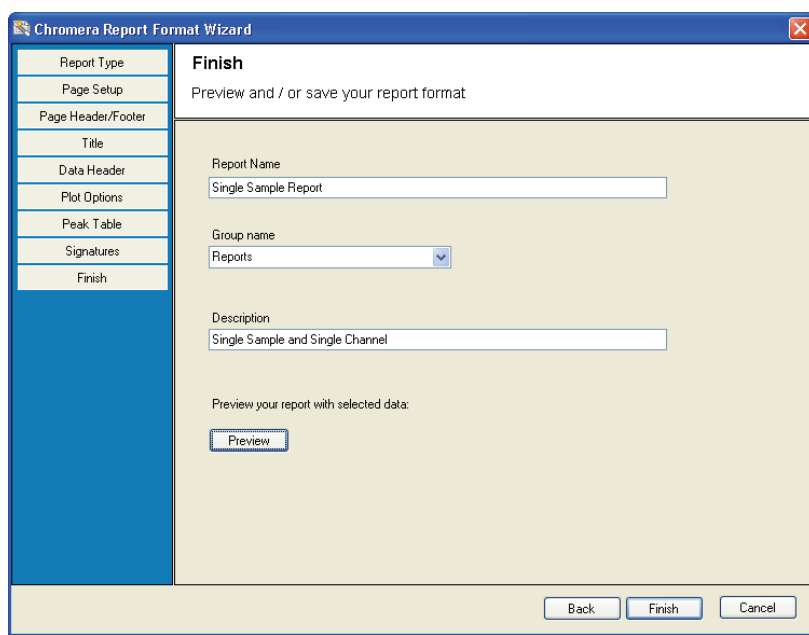
Signature line:

Signature: _____ Name: _____ Date: _____

Buttons: Back, Next, Cancel

Two different **Signature** lines can be displayed for hand-written signatures and timestamps.

24. Click **Next** to display the **Finish** screen.



Chromera Report Format Wizard

Finish
Preview and / or save your report format

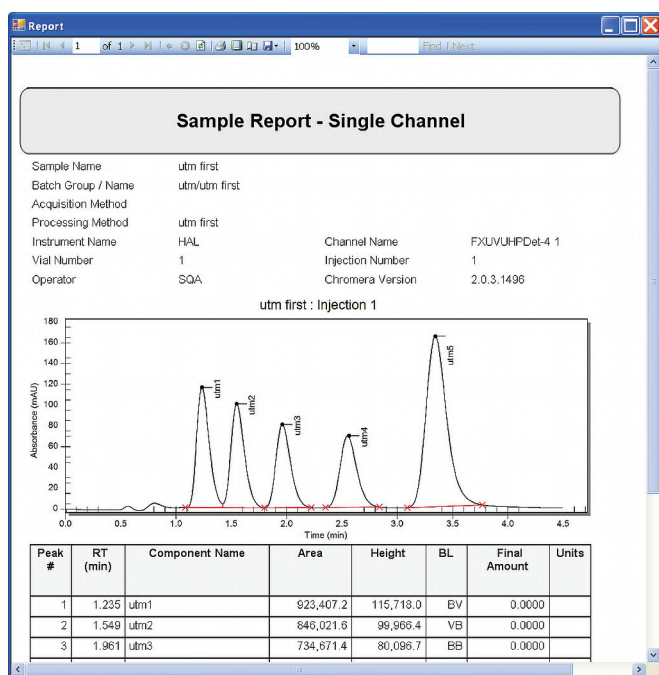
Report Name
Single Sample Report

Group name
Reports

Description
Single Sample and Single Channel

Preview your report with selected data:


25. Enter a **Report Name**, **Group name**, and **Description** for this report.
26. Click **Preview** to display a visual rendering of the report using mock data and graphics.




This is an example of the print **Preview** for the completed report. The preview contains only mock data and plots, but the layout of the report should be as specified through the wizard.

You have completed Part 3 of the Chromera tutorial.

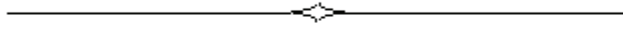NOTE: Every step in the rest of the tutorial is built on, and is a continuation of, the previous step. Therefore, you need to complete it as a whole.



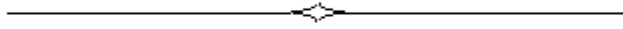Tutorial Part 4: Viewing Results **and Reprocessing Data**




Viewing Results in Post Run

Whether reprocessing existing data or acquiring new data, the completed samples will be displayed in the **Post Run** environment, and can be inspected by navigating through the Sample Tree and interacting graphically with the chromatographic display.


- Data can be treated as view only from the standard **Post Run** display.
- Individual results can be optimized graphically.
- The current version of the method can be graphically modified (**GME**, Graphic Method Editing) using the selected sample data.
- Data can be viewed in **Single Plot** mode, **Stacked Plot** mode, **Matrix** mode for multiple channels and replicate injections, or in **Overlay** and **3D** mode.

To view results in Post Run:

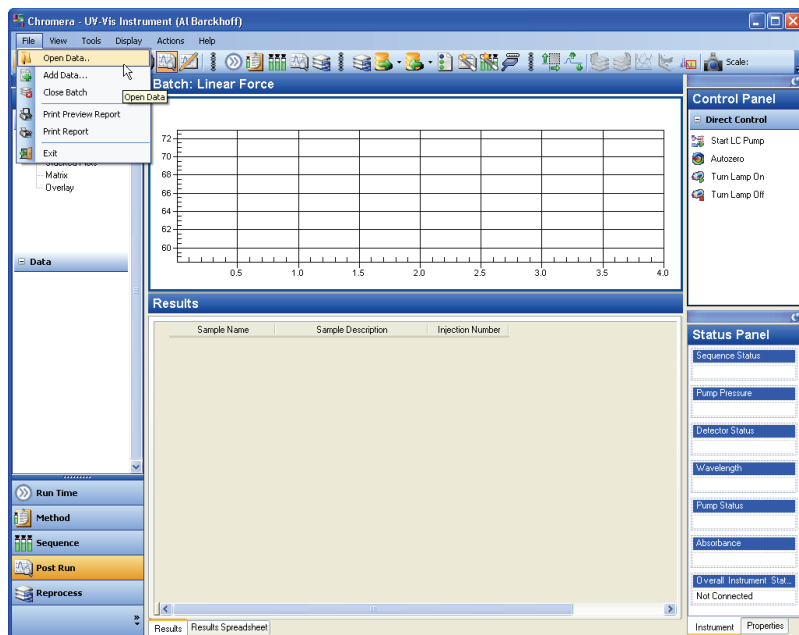1. Click the **Post Run** button in the navigation pane.
2. Select **Open Data** from the **File** menu.



The **Data Selector** opens.

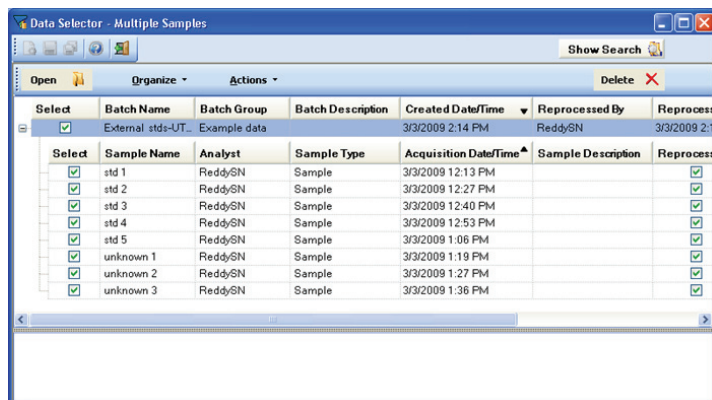

3. Select the **Batch** listed in the Data Selector you want to open.

A Batch in Chromera consists of the entire collection of data and metadata used to produce a set of results. A batch may represent a single injection of a sample on a single detector instrument, or it may represent many samples, with many replicate injections, on an instrument with multiple detectors and multiple data channels. The creation of a batch is directed by a sequence.

A batch contains:

- All the sample identification specified in the sequence.
- All the original raw data points generated from that sequence.
- All the original method settings used to control the instrument and acquire the original raw data. (When possible.)
- The current and complete set of integration parameters used to generate the current set of peaks.
- The current and complete set of calibration information used to generate the current set of results.
- All the results generated from the processing of the raw data using these "embedded" method parameters.

Data acquired from the **Sequence** mode in **Run Time** will be written to a new batch containing all the data from all the injections made from that sequence.

Data acquired from **Single Injection** mode in **Run Time** will cause a new record to be written to a new batch created by Chromera, containing all the data from each injection made, until the current **Run Time** mode is changed.

NOTE: *When TotalChrom data are imported into Chromera, a new batch is automatically created, containing the raw data and sample information from the TotalChrom RAW files, and as much of the TotalChrom Instrument Method as possible. Since results are not imported from TotalChrom, the imported batch will not contain any peak information.*

If data are processed more than once, the new results may be either appended to the existing batch, or written to a new sequence containing all of the information described above. The choice is driven by the nature of the selections and changes made in the Reprocessing environment.

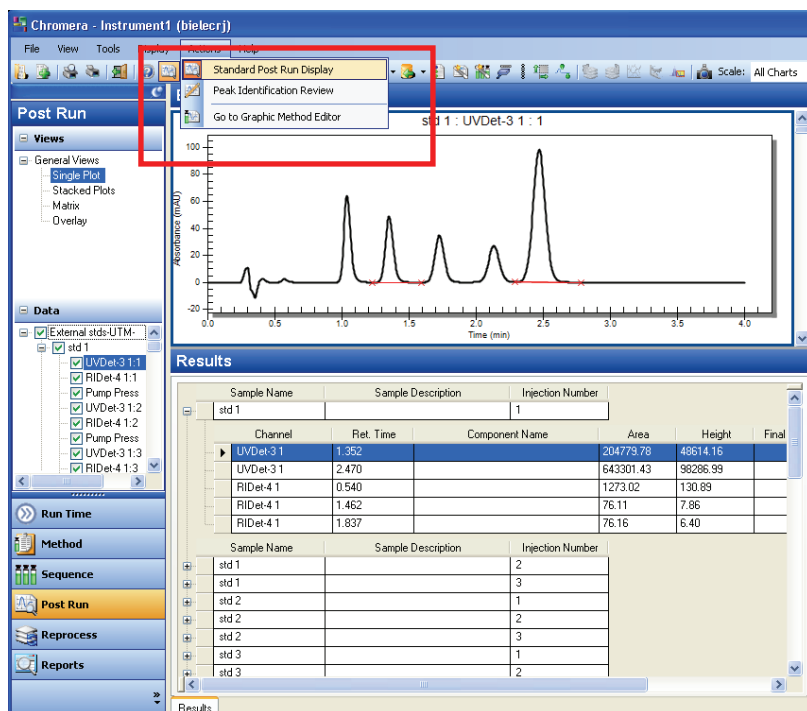
Results are never overwritten or modified in a batch. New results are just appended to the data in the existing batch. Open the desired Batch of data to view through the **Data Selector**.

- ♦ A database search query can be constructed to refine the selection of data.
- ♦ Searches can be saved and recalled for convenient retrieval of data fitting routine criteria.

4. Click **Open**.

Modes of Operation

In the **Post Run** environment depending upon the mode selected, you take one of three actions from the **Actions** menu. You can either view data without changing anything (**Standard Post Run Display**), optimize the integration of individual injections (**Peak Identification review**), or optimize the current version of the method in the database using one of the injections as an example (**Go to Graphic Method Editor**).



Peak Identification Review

The Peak Identification Review allows optimization of the currently selected chromatogram, by permitting changes to the integration parameters for that channel data, and that channel data only, thereby generating potentially new results for that data. These results are appended to the same batch.

Go to Graphic Method Editor

This functions the same as Graphic Method Editor, but only for the data on that Channel. It launches the method editor for current version of the method from the database, and displays the currently selected data in **Graphic Method Edit (GME)** mode. This action leaves the Post Run mode and automatically loads into the Method Editor the current version of the method for the selected channel data.

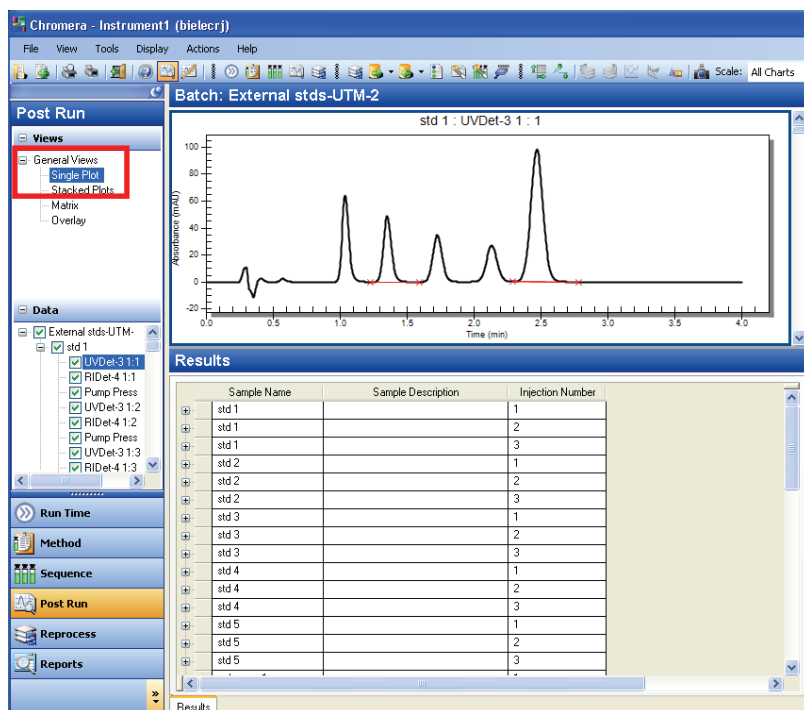
General Views of Plot Displays

Navigate through the **Sample Tree** in the **Navigation Bar** to select the **Channel Data** and **Sample Replicate** to view. The plot mode will display one or more chromatograms at a time. The results for the currently selected **Data Channel** will be displayed in the lower frame. Plot Options control the display of annotations, peak information, peak fills, etc.


Single Plot

To display a single plot:


1. Select Single Plot from the General Views.
2. View the chromatogram for each Data Channel one-at-a-time.
In this example, **UVDet-3 1:1**, **RIDet-4 1:1**, **Pump Pressure:1**, etc.
3. Navigate through the **Sample Tree** or scroll through the plots.



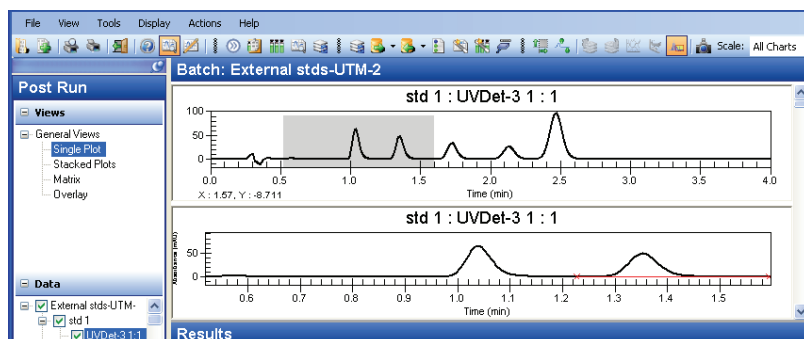
4. You can optionally display a **Radar Window** from the **Display** menu. The Radar Plot indicates the localized region of the main plot that is currently in view.



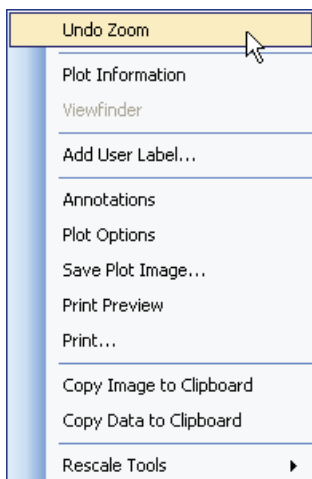- Clicking-and-dragging in the **Radar Window** will display that selected region of the chromatogram in the main plot.
- Zoom and un-zoom using click-and-drag and selections from the context menus.



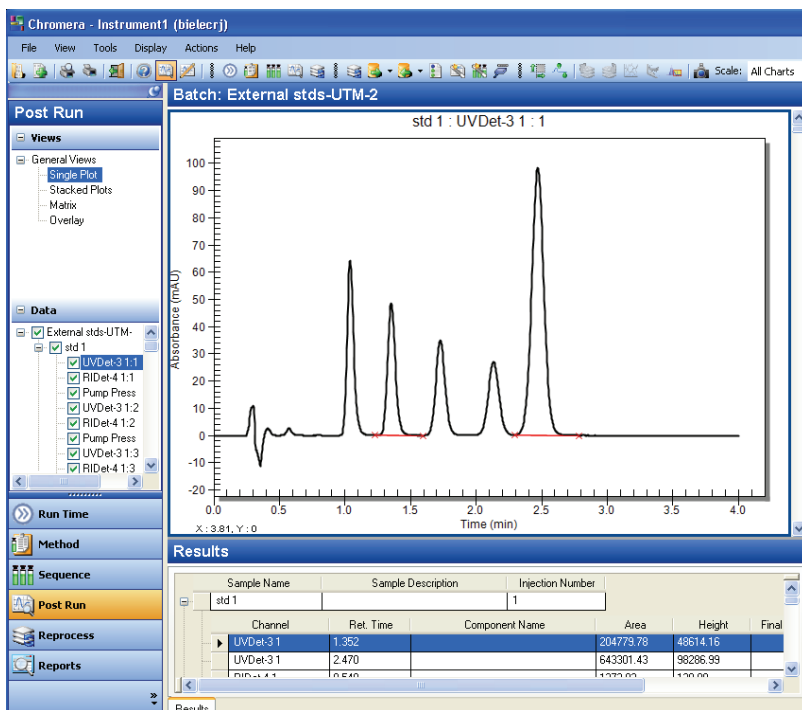
- Expand the size of the chromatogram frame either by dragging the splitter between it and the **Results** frame, or by double-clicking on the Title Bar for the chromatogram.




Stacked Plots

Stacked plots enable you to view more than one plot at a time in a single column. To display stacked plots:

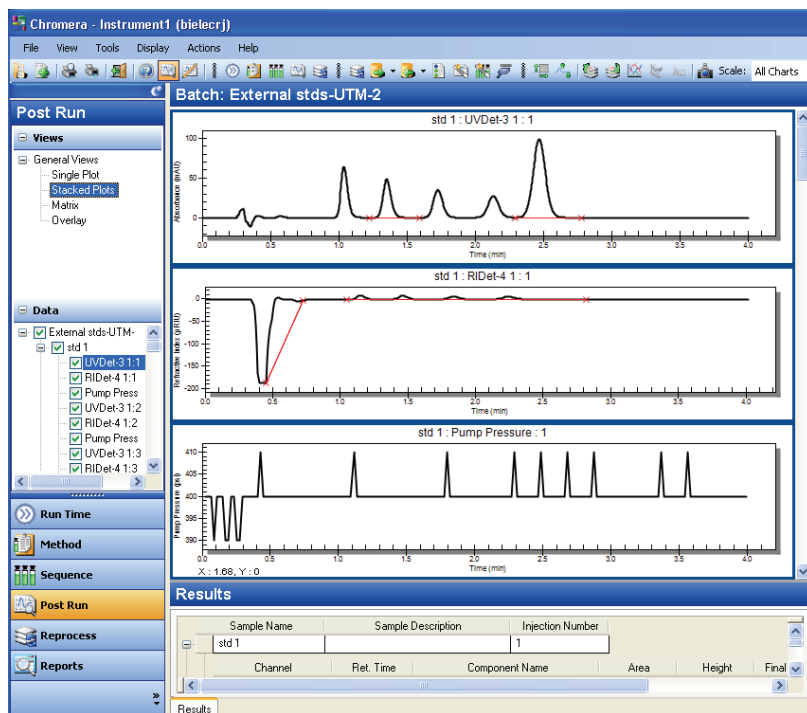1. Select **Stacked Plots** from the **General Views**.
2. Navigate through the Sample tree or scroll through the stacked plots.



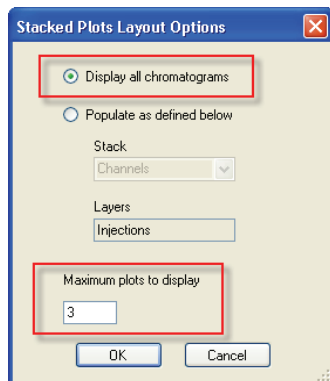
3. The number of chromatograms to display and the order of scrolling through Replicates and Channels can be controlled by selecting **Plot Layout** from the **Display** menu.




Matrix

The Matrix mode allows multiple Data **Channels** or **Replicates** to be shown simultaneously for each **Sample**. To display a plot matrix:

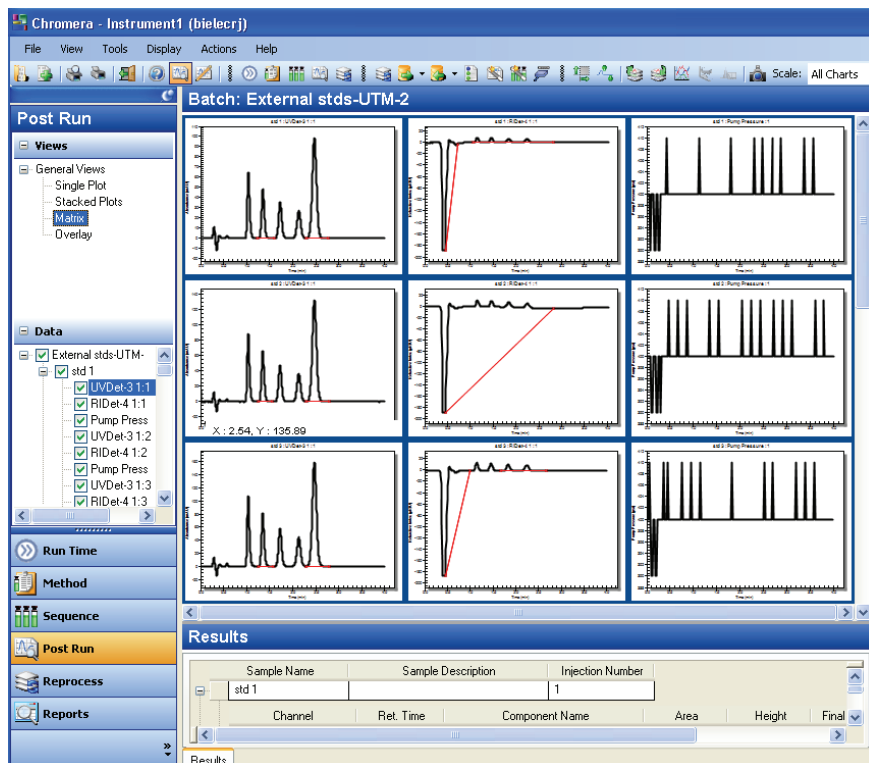1. Select **Matrix** from the **General Views**.



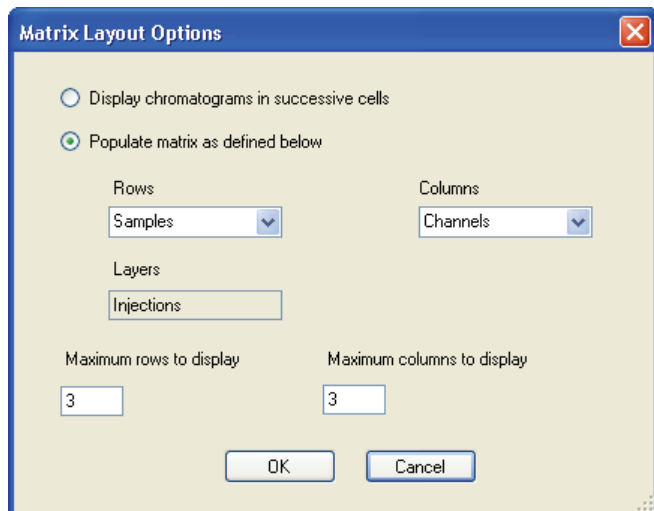2. Selecting the **Matrix Layout Options** control enables you to display and navigate through the **Layers** of **Samples** and Replicate **Injections** and Data **Channels**.




Overlay

The **Overlay** mode allows up to 30 individual chromatograms to be displayed in the same plot window. Use **Plot Options** to control the colors to be assigned to each of the chromatograms. You can select specific plots to be displayed or not displayed can be selected from the **Sample Tree**.

To overlay plots:

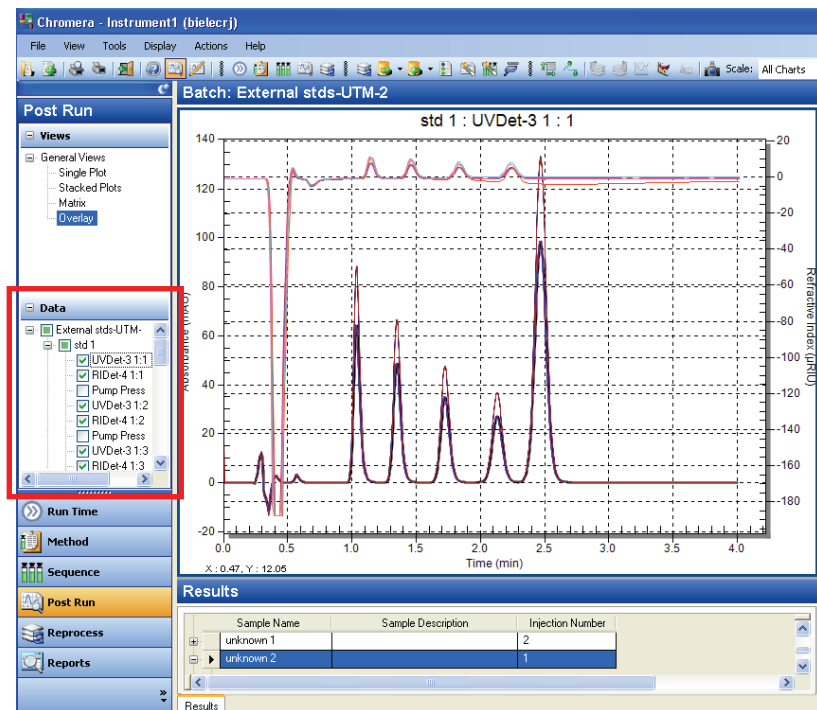1. Check and uncheck the data items to overlay.




3-D Display

When in **Overlay** mode, you can enable the **3-D** mode. Use the horizontal and vertical scroll bars to rotate the plots around the X, Y, and Z axes. You must be in the **Overlay** mode to undo the **3-D** mode.

To display plots in a 3-D display:

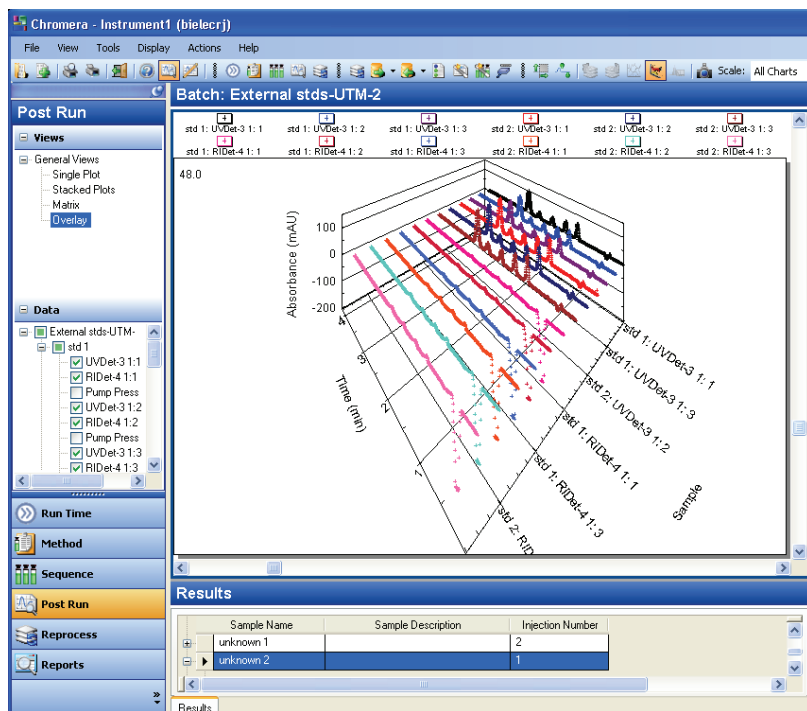1. Select **3-D** from the **Display** menu.



2. Double-click on the 3-D plot to start or stop auto-rotation.

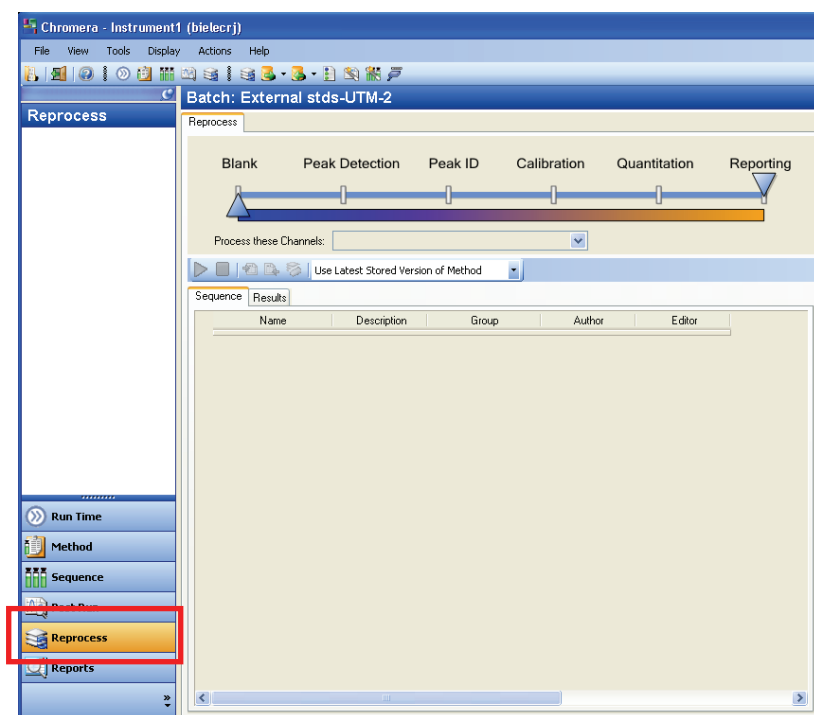
Reprocessing the Data

Once you have acquired data and stored it in a raw data file, you can process that data repeatedly. The Batch Reprocessing function enables you to reprocess a large number of files in succession.

One of the purposes of reprocessing data is to re-analyze the raw data files listed in a sequence. For our last exercise, we will use the sequence file we just created to reprocess data files with our optimized method.


To reprocess files in a sequence:

1. Click the **Reprocess** button in the Navigation pane.




NOTE: If you selected the Batch Reprocessing function when you had a sequence open, that sequence is supplied as the default. The **Reprocess** environment consolidates all the controls needed to selectively reprocess a batch of data.

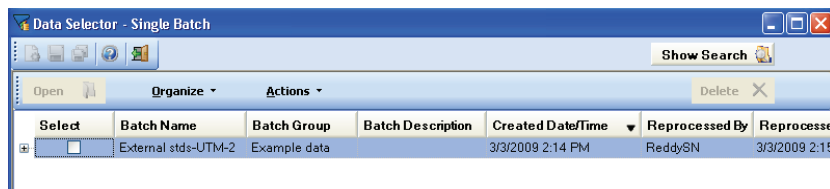
2. Select **Open Data...** from the **File** menu to use the **Data Selector** to open a **Batch** of data to reprocess. New **Searches** or saved **Searches** can be used to refine the criteria for retrieving data from the database.



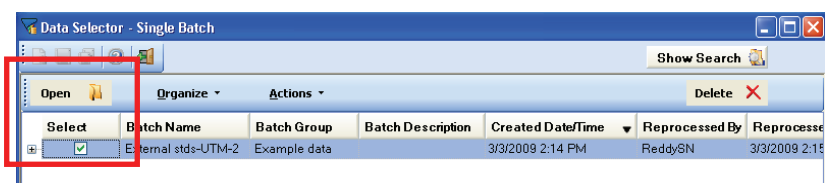
3. Click in the **Select** checkbox to select the batch, then click **Open**.



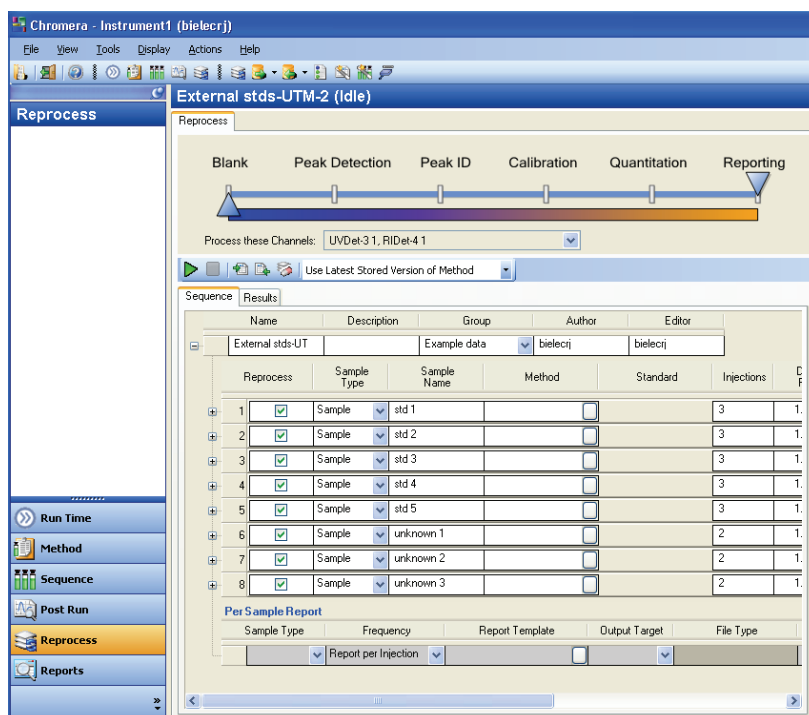The screen is populated with the data to be reprocessed.



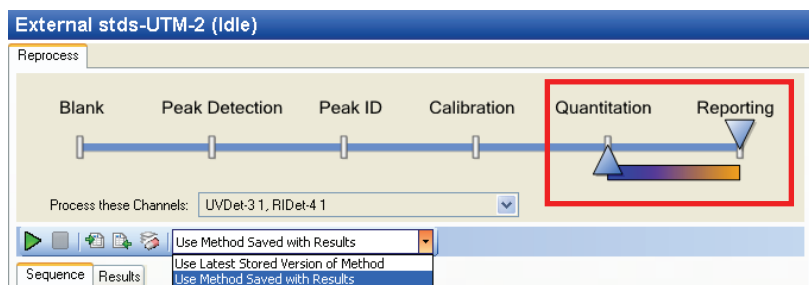The amount of processing to be performed is controlled by using the processing **Slider Bar**.



4. Select the method version from the drop-down list to use for processing. This is disabled when the batch is being processed. Options available are:

Use Latest Stored Version of Method: The current version of the method specified in the sample row is read from the method database and used for reprocessing.

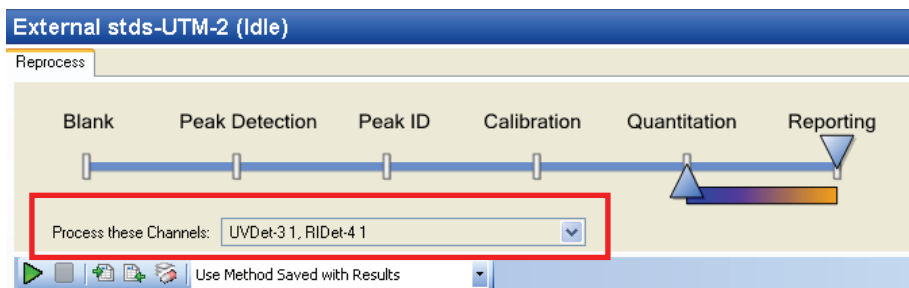
The starting point for reprocessing is automatically set to **Quantitation**, because **Peak Detection**, **Peak ID**, and **Calibration** information will be used from the existing results.

Use Method Saved with Results: The version of the method saved with the results are used for processing (processing must start with Quantify or later – see below).


NOTE: Rows in which a new method name is selected will always be processed with the current version of that method.

When selecting a different method for processing, a new batch will be created using either the **Name** of your choice, or the same name as the existing batch but with the timestamp of the reprocessing.


5. Select from the drop-down, one or more of the **Data Channels** to be reprocessed.



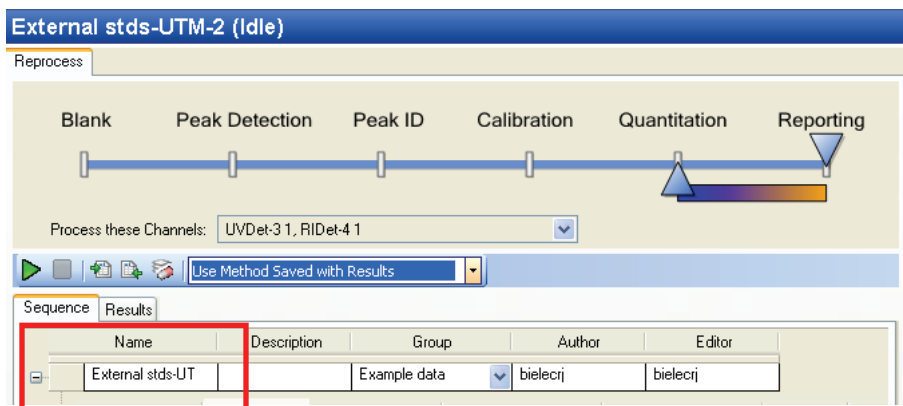
Changing any information in the sequence will cause a new batch to be created using either the **Name** of your choice, or the same name as the existing batch but with the timestamp of the reprocessing.



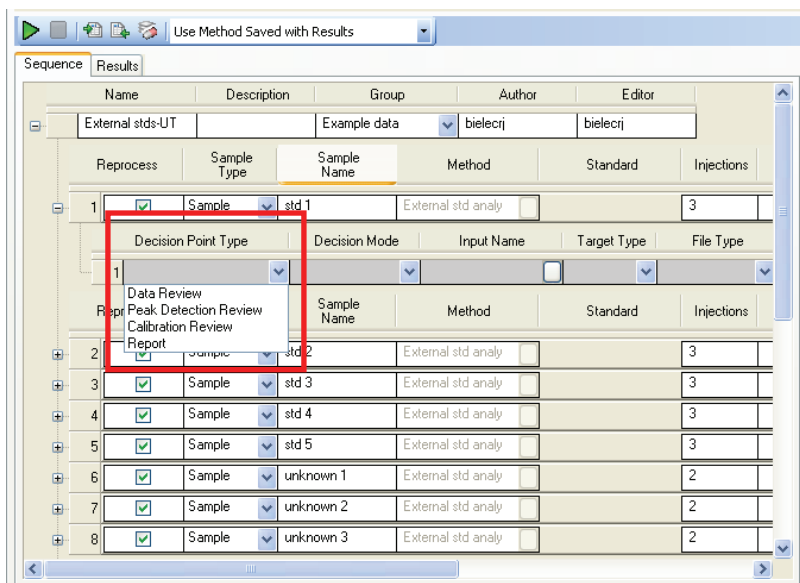
6. Set desired **Decision Points** between samples and injections to automatically pause the **Batch** to allow manual review of the data, calibration, or peak detection.



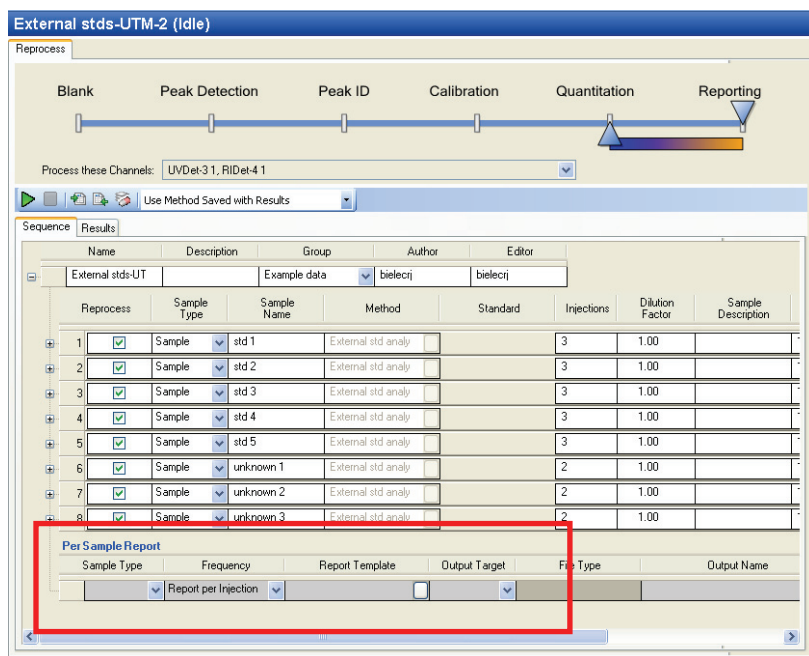
7. If report templates have been created for different Sample Types, **Per Sample Reports** can be generated on a **Sample Type** basis.



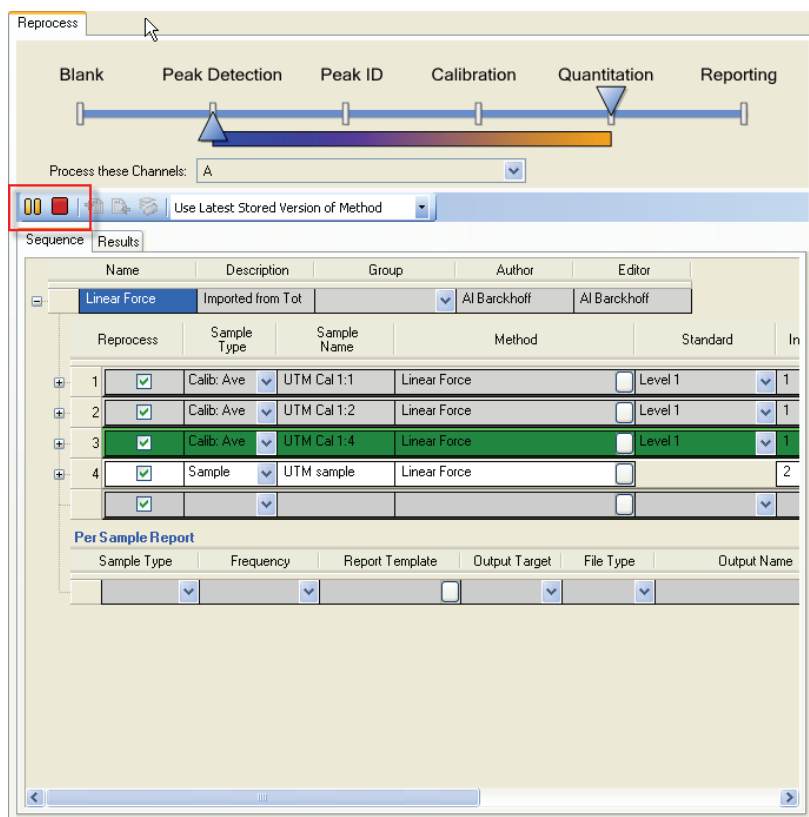8. Use the controls on the toolbar to **Start** and **Stop** the batch processing as desired.




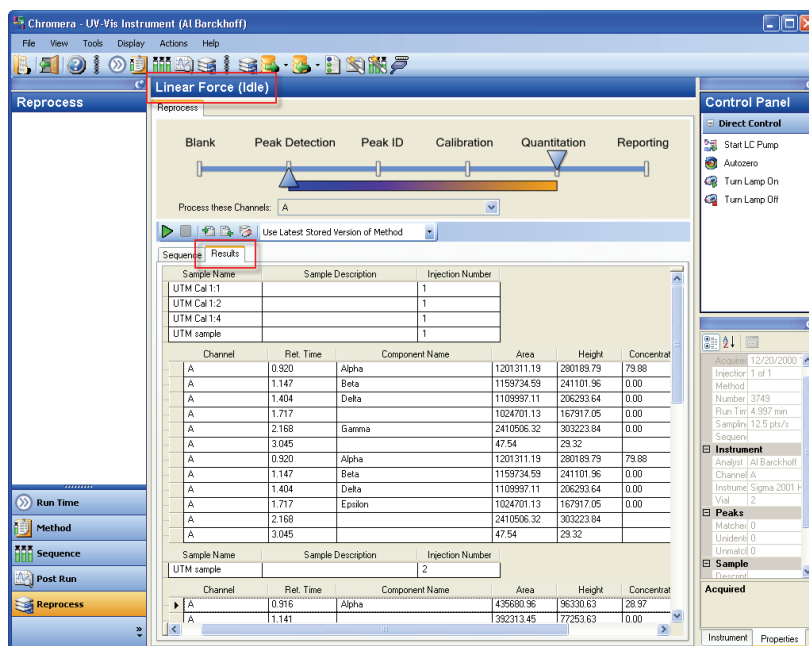9. Select the **Results** tab to view the data as they being generated.
No quantitated results are shown for the analysis of calibration standards.




10. The data are automatically loaded into the **Post Run** environment for either review or optimization. Now we are ready to use this sequence with the Batch Reprocessing function.

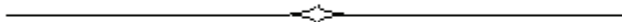The next section describes how to set up and run a simple sample analysis



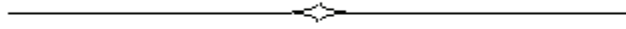Tutorial Part 5: Running a Simple Analysis



Overview

In this part of the tutorial you will:

- Configure an Instrument (pump, autosampler, detector, oven) in Chromera Manager
- Create and save a Method
- Run the Method in Single Run Mode
- Acquire data
- Create and save a Sequence Mode
- Run the Sequence

Prepare Your System

Prepare your system for this tutorial with the solvents (mobile phase), column, and sample for you're your analysis.



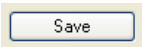
- LC grade methanol and water. Reservoir A = Water; Reservoir B = Methanol
- Universal Test Mix (P/N 00890893, 5 vial pack)
- 3x3 CR C18 column (P/N 02580195) and column holder (P/N 07150028)


Creating an Instrument and Database in Chromera Manager

This example briefly describes how to configure an instrument to run a simple analysis. In this tutorial example the instrument consists of a Series 200 Quaternary Pump, Series 200 Autosampler, Flexar UHPLC detector, and Series 200 Column Oven. For more details refer to **Create an Instrument** on page 10.

Create an Instrument

To create an instrument, follow this procedure:

1. Click on the **Configuration** button  to display the initial Configuration screen.
2. Click in the **Name** field and type a name to define the name of the Instrument (for example LC Tutorial).
3. Click the plus button  to the left of **Name** to expand the entry fields so that you can enter each of the **Device Names** to make the instrument. Click in the **Device Name** field and select a Device Name from the drop-down list.
4. Click in the **Port Name** field and select a Port Name from the drop-down list.
5. In this example we are not going to enter a **Device Description**, **Serial Number** or **Firmware Version**.
6. To add another device, click in the next **Device Name** field and select a Device Name from the drop-down list. Continue to enter the information for the device.
7. Set the **DB Connections** settings. Leave the Server setting but change the DBName to Tutorial, and ArchiveDB to TutorialArchive. Each instrument can have its own database.
8. When done, click the **Save** button  to save this Instrument Configuration. The Instrument is displayed under the user name in the Navigation Pane.

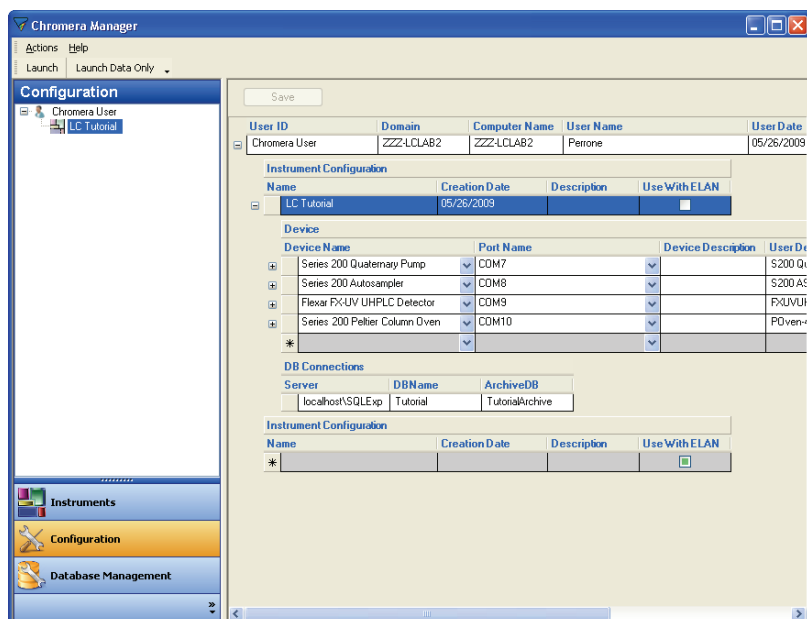


Create a Database

Before running Chromera you must create a database. To create a database, follow this procedure:

1. Click on the **Database Management** button  to display the database screens.

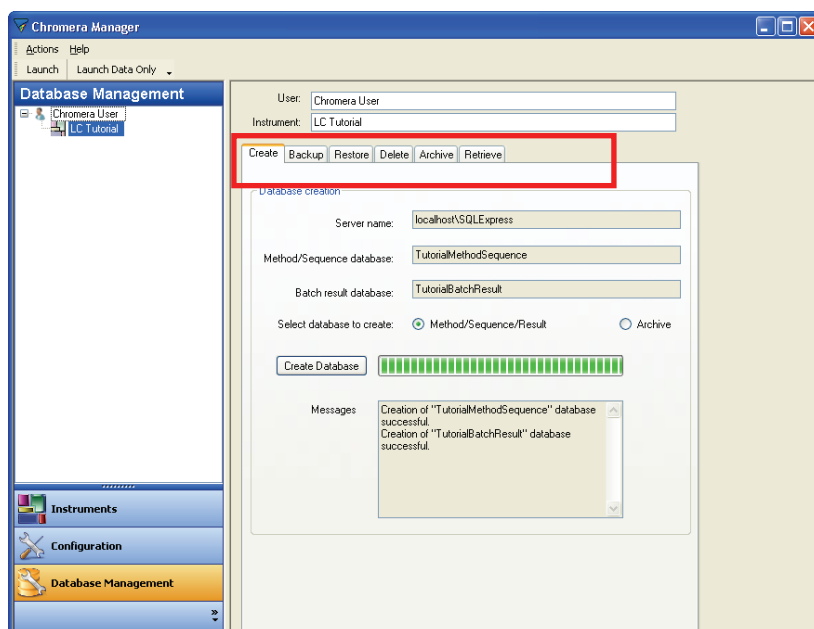The database functions (**Create, Backup, Restore, Delete, Archive, and Retrieve**) are displayed as tabs. Click on the tab to display more detail of your desired function. In this tutorial click the first tab to **Create** a database.



Notice that the **Server name** is filled in. This name corresponds to the default **DB Connection** settings that were made when you configured an Instrument.

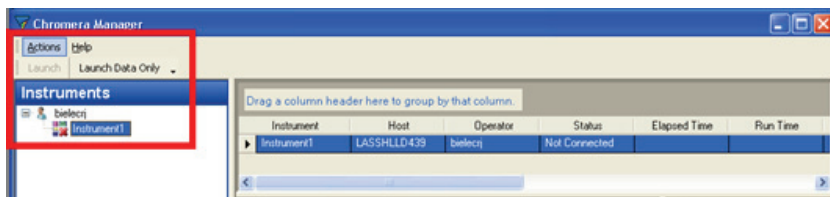
2. In this example we are creating the **Method/Sequence/Result** database.
3. Click the **Create Database** button.

Observe the progress bar to the right of the Create Database button as the database is created. When the database creation is complete a message appears in the **Messages** box indicating this.


Launch Chromera

The Chromera application is a separate application from Chromera Manager and must be launched (started). To launch Chromera, follow this procedure:

1. Select the **Instrument** you want to initialize and click the **Launch** button ... or simply double-click the instrument name.




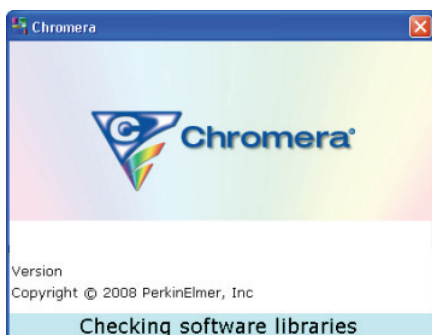The system does extensive checking of the Chromera installation to ensure that the database is intact, proper program files have been installed, etc.




Any normal errors or warnings encountered during startup are displayed.

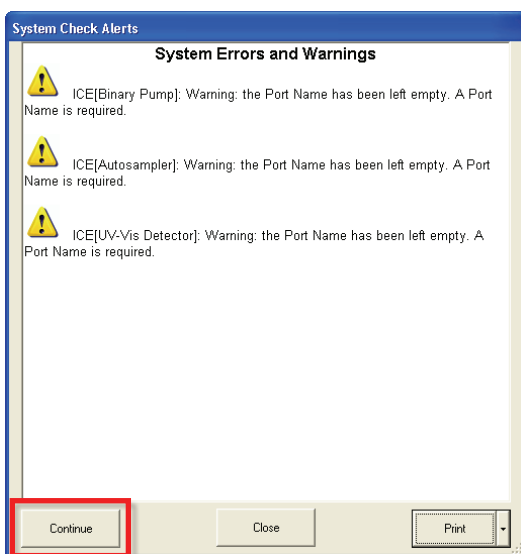2. Once the system is running, instrument errors are displayed in the fade-in notification area, and will be written to the **Error Log**, accessible from the Chromera **Tools** menu. Unexpected errors are trapped at the OS level and recorded in the Windows **Event Log** under **Application Errors**.



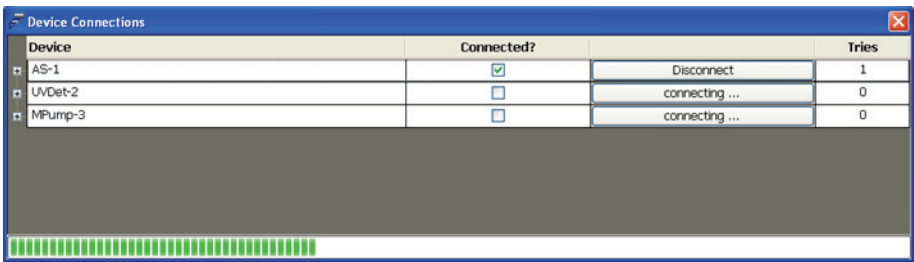3. Once Chromera has been launched, the **Device Connections** dialog displays showing the status of the connections. As the Device tries to connect, **Connecting ...** displays in the button.



If a device does not connect, click the **Connect** button to retry making a connection to the instrument module. You can display the **Device Connections** dialog at any time from the Chromera **Tools** menu.

After all Device Connections are made Chromera displays in the **Run Time** mode.

Using Run Time

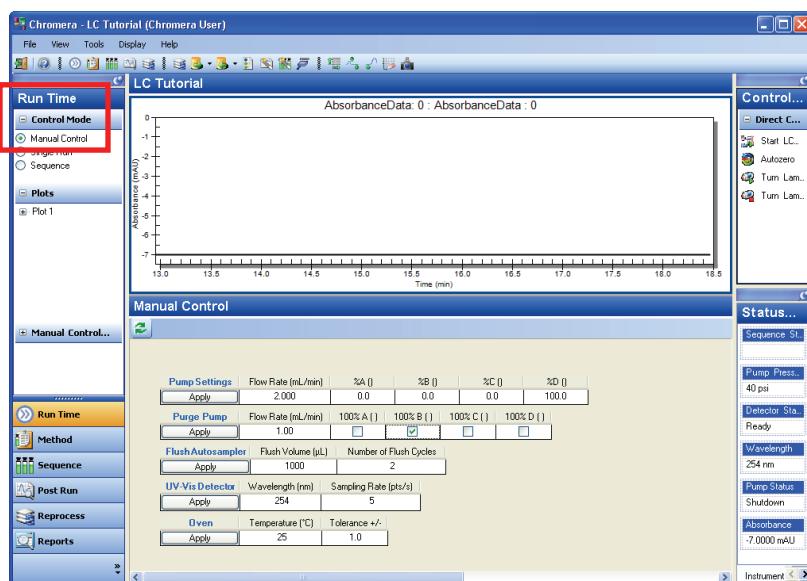
The Run Time mode allows you to collect data without developing a comprehensive method. You enter only the parameters that are essential for performing a run. This is helpful when you are first starting to use Chromera; for example, if you are performing a survey run on a type of sample that you have not analyzed before or when you are developing instrument settings for an analysis.

In this tutorial we will create a simple **Method** to use in **Single Run** mode. We will later use this method when we build a sequence. Before creating a method let's set initial conditions using the Manual Control Mode.

Manual Control Mode

Use Manual Control to directly set initial conditions for system startup, equilibration, and detector monitoring. Here you will purge the pump, set the initial conditions and equilibrate the system, flush the autosampler, and verify that the detector is properly working. This is essentially the Direct Control mode of the software.

- Select the **Manual Control** mode.



Purge the Pump

You need to purge the pump to remove any trapped air from the solvent lines and fill the lines from the solvent reservoirs to the pump with the selected mobile phase.

1. Make sure the pump solvent inlet lines are in **Reservoir A: Water** and **Reservoir B: Methanol**.
2. Locate the Prime/Purge valve on the front of your pump.
You will have to either open the pump door on a Series 200 pump or on a Flexar pump remove the front cover.
3. Attach the 30ml priming syringe to the Prime/Purge valve on the pump.
4. On the **Purge Pump** settings, type **10** for the **Flow Rate** and click in the **100% A** checkbox.

Purge Pump	Flow Rate (mL/min)	100% A ()	100% B ()	100% C ()	100% D ()
Apply	10	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

- Open the Prime/Purge valve by turning it counterclockwise.
- Click **Apply** to start the pump.
- Collect about 25 to 30 mL in the syringe then click **Stop LC Pump** in the **Control Panel**.

Control Panel	
Direct Control	
	Stop LC Pump
	Autozero
	Turn Lamp On
	Turn Lamp Off

- Close the Prime/Purge valve by turning it clockwise, then remove and empty the syringe.
- Attach the 30ml priming syringe to the Prime/Purge valve on the pump.
- On the **Purge Pump** settings, click in the **100% B checkbox**.
- Open the Prime/Purge valve by turning it counterclockwise.
- Click the **Apply** button to start the pump.
- Collect about 25 to 30 mL in the syringe then click **Stop LC Pump** in the **Control Panel**.
- Close the Prime/Purge valve by turning it clockwise, then remove and empty the syringe.

Equilibrate the System

Equilibrate the system with your mobile phase concentration to provide a stable analysis.

- On the **Pump Settings**, type **2** for the Flow Rate, **30 for %A** (Water) and **70 for %B** (Methanol).

Pump Settings	Flow Rate (mL/min)	%A ()	%B ()	%C ()	%D ()
Apply	2.000	0.0	0.0	0.0	100.0

- Click the **Apply** button to start the pump with these settings.
Let the system equilibrate for 10 minutes.

Flush the Autosampler

While the system is equilibrating, flush the autosampler.

- Click the **Apply** button to flush the autosampler.

Flush Autosampler	Flush Volume (μL)	Number of Flush Cycles
Apply	1000	2


The autosampler will flush using the displayed **Flush Volume** and the **Number of Flush Cycles**. In this example, it will flush **1000 μL two** times.

Verify that the Detector is Working

1. While the system is equilibrating, change the wavelength to **260** nm then click the **Apply** button.

UV-Vis Detector	Wavelength (nm)	Sampling Rate (pts/s)
Apply	254	5

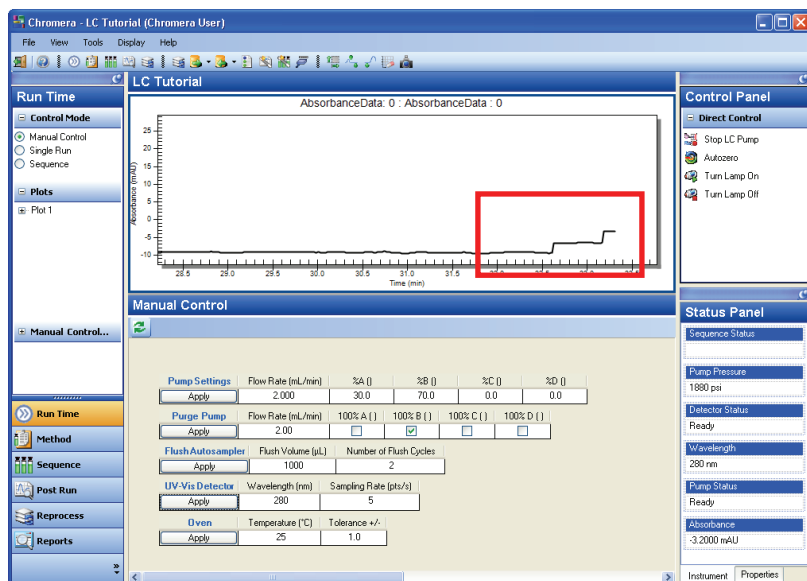The baseline should change on the display.




Create a New Method

While the system is equilibrating, create a method to be used with your analysis.

To create a new method for the current instrument configuration:

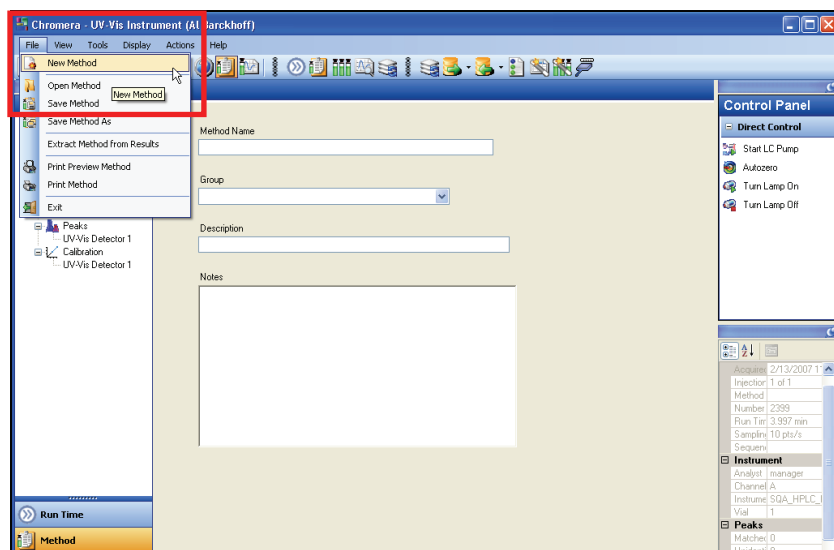
1. Click the **Method** button to select the method environment. Then select **New Method** from the **File** menu.



2. Type a **Method Name**, **Group** name, **Description**, and **Notes**.


The **Group** name is like a folder, and is a **required field**; whereas **Description** and **Notes** are optional.

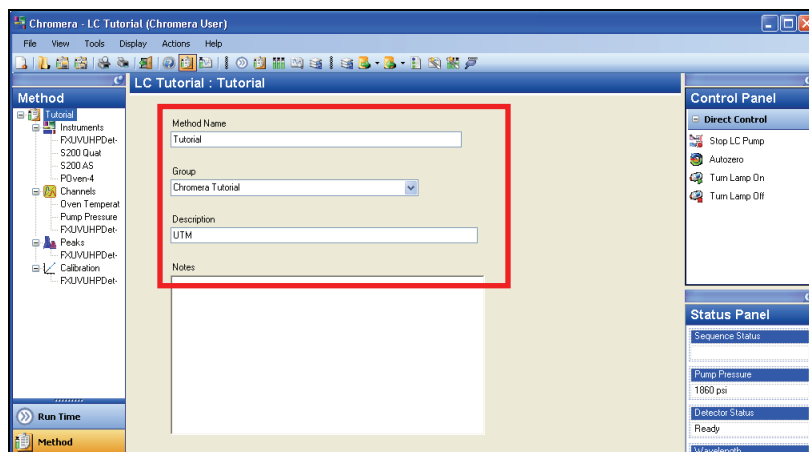In this example the Method Name is **Tutorial**, the Group is **Chromera Tutorial**, and the Description is **UTM**.




The Method Tree (in the **Method** pane on the left) provides access to all the instrument settings or to each module separately.

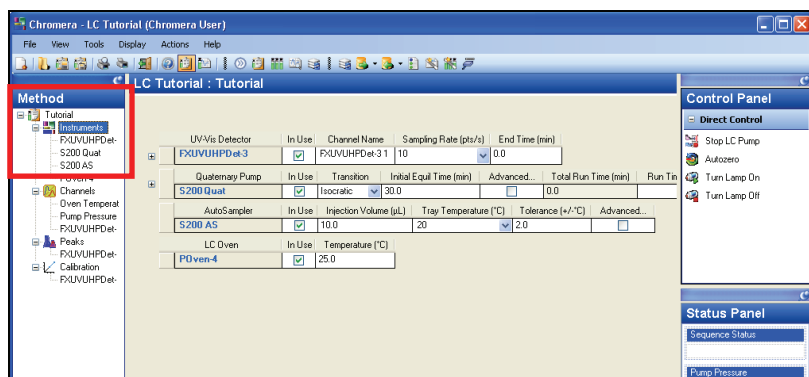3. Click **Instruments** in the Method pane to display the devices and their parameters. Enter the required HPLC parameters.



This example has a quaternary pump configured for this Instrument in Chromera Manager.

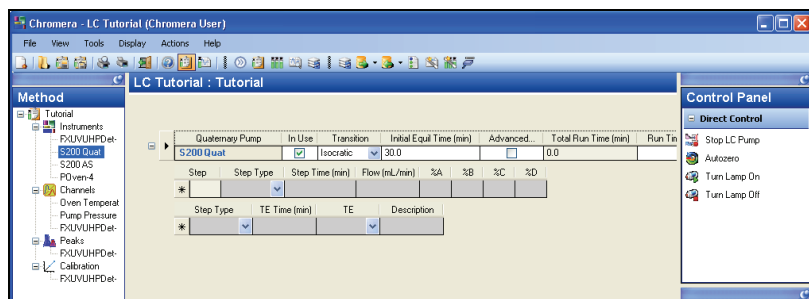
4. Click **S200 Quat** and observe that just the parameters for the **Quaternary Pump** display. Click the plus sign **+** to display all of the pump parameters.




5. Click in the **Step** field and enter the parameters for this step, Step **0**.

76 . Chromera Tutorial for New Users

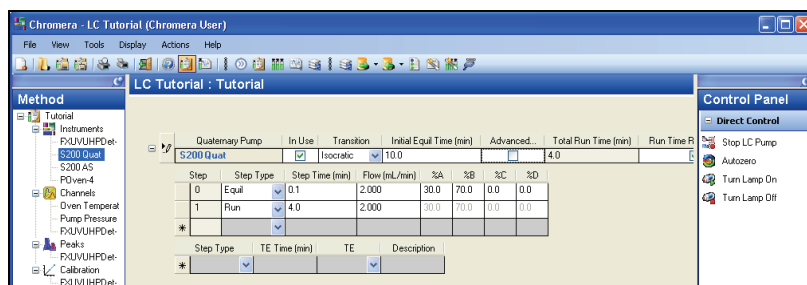When done, click in the field under Step 0 and enter the parameters for Step **1**. Leave the rest of the parameter as the default parameters.




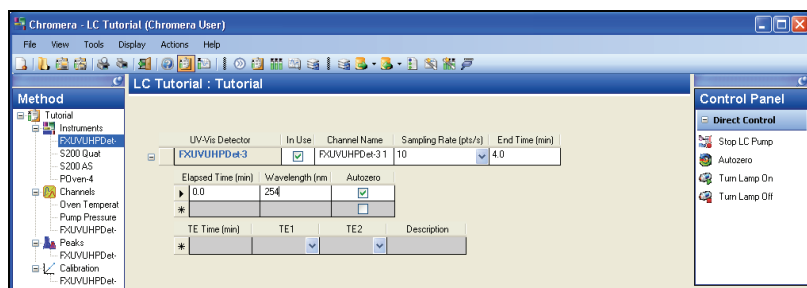- Click **FXUVUHPDet-3** and observe that just the parameters for the **FXUVUHPDet-3** display. Click the plus sign **+** to display all of the pump parameters.




- After entering all of the parameters for all of the devices, select **Save Method** or **Save Method As** from the **File** menu.

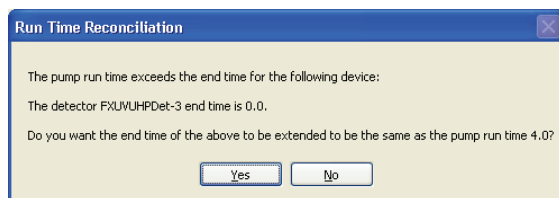You may see a message similar to the following:



Click **Yes** and the method will be updated.

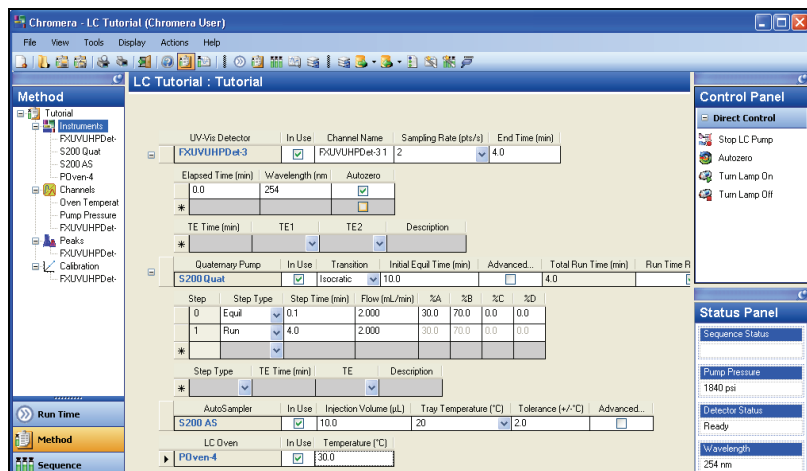
- Click **Instruments** and view your entire method. Click each plus sign **+** to display all of the parameters.



Running the Method in the Single Run Mode


To run the method, go to Run Time and make the method you just saved your Active method.

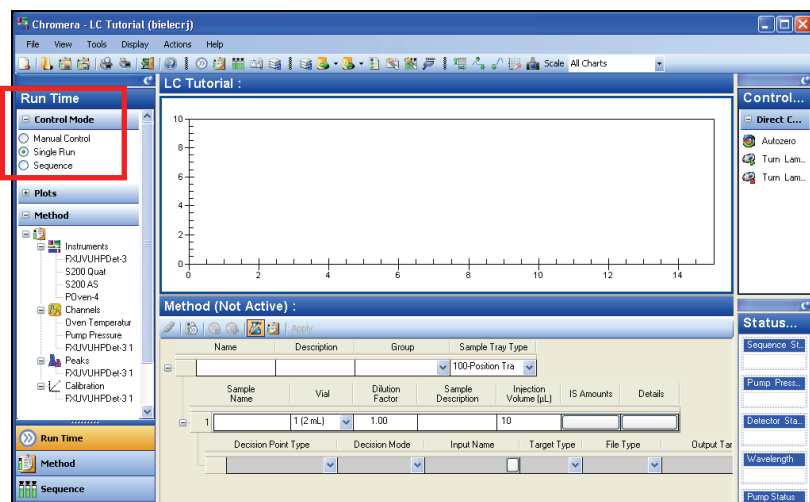
1. Click the **Run Time** button, then select the **Single Run** Control Mode.




Now you need to open the method you saved.

2. Select **Open Method** from the **File** menu.

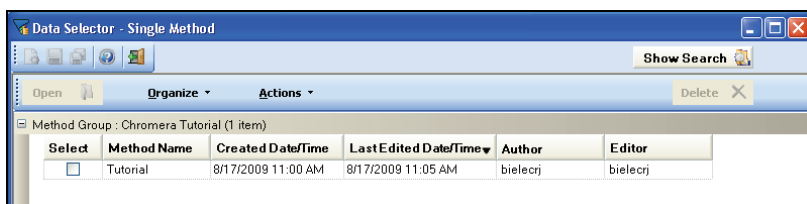
The **Data Selector – Single Method** displays. Click the plus sign  to display the available methods. This example shows on one method – **Tutorial**.



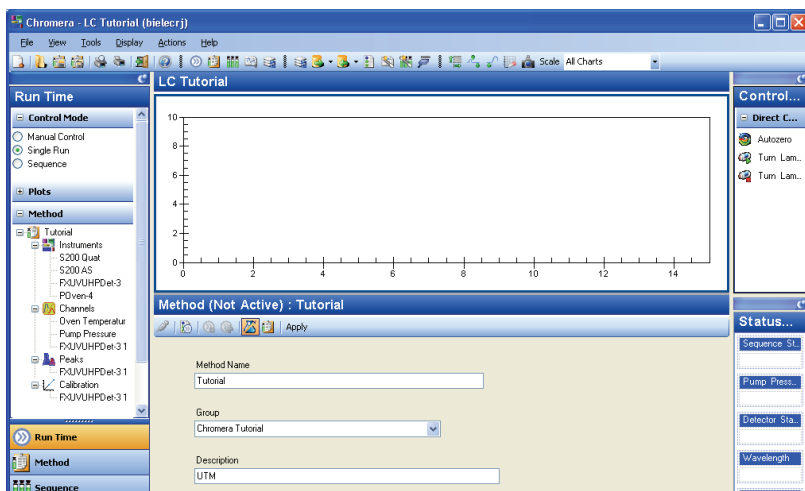3. Click in the **Select** box to the left of the **Tutorial** method, then click **Open**.
The Method appears as a **Not Active** method:



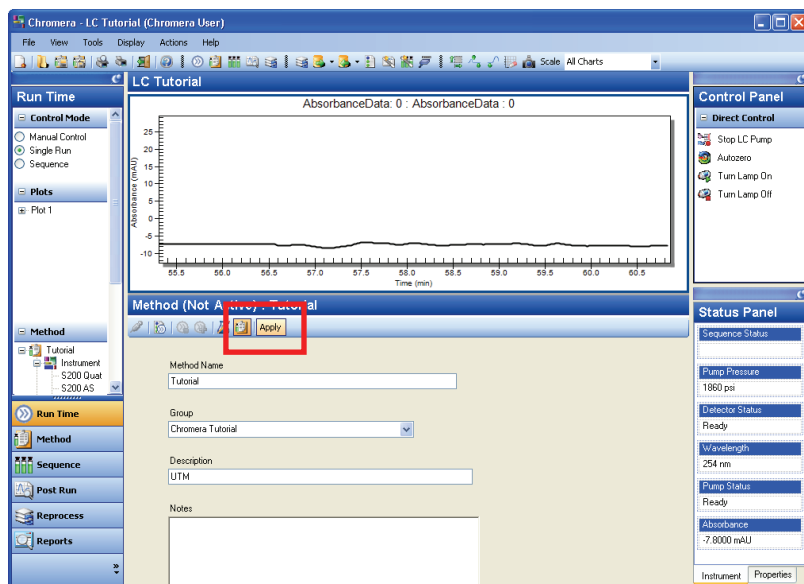- Click **Apply** to activate your method.




An **Active Method** is the method currently in control of the instrument and can be edited or modified as required while it is running.

- Any step of the method not yet begun can be changed.
- Any processing or calibration change can be made.
- Apply the changes to the method without saving them.
- Save the method and have the changes apply to all subsequent injections.

- Click **Autzero** in the **Control Panel**.

- Click the Syringe button  to start the method run.

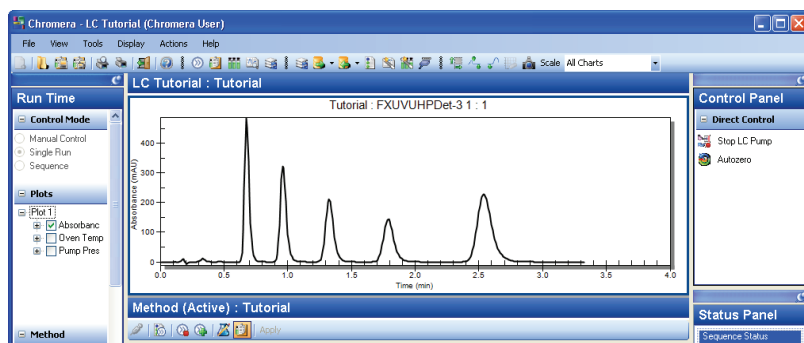The Autosampler starts to inject and the display a Real-Time Plot. The Real-Time Plot lets you view a plot of the data as it is being acquired by an interface. It autoscales to the largest peak.



When running in Single Run every run creates a new batch as each injection becomes a separate entry in a Chromera Batch. Each injection is saved with its complete method, the raw data, and any results generated.

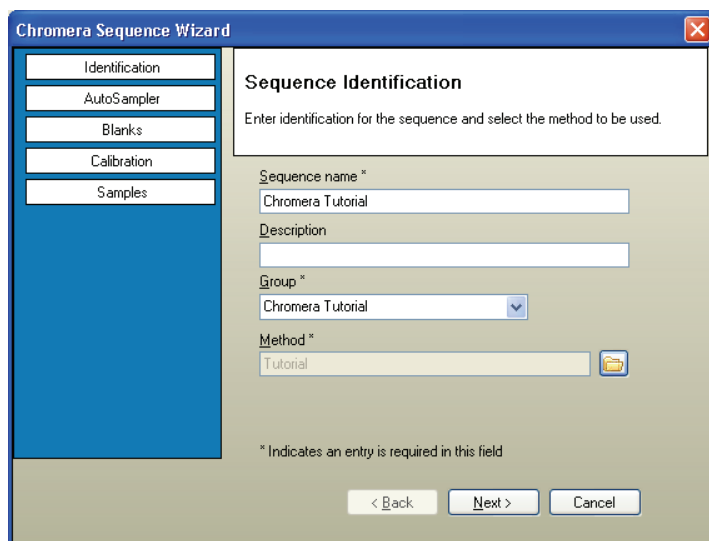
For most labs the proper way to run Chromera is through a Sequence. Now lets build a simple sequence.

Building a Sequence

Sequences are used to control acquisition and analysis of chromatographic data produced by a series of injections. Use the **Sequence Wizard** to create a sequence for running new samples and acquiring data. It is also possible to create sequences through the Sequence Editor.


Using the Sequence Wizard

1. From the **Tools** menu select the **Sequence Wizard** to create a new sequence to acquire new data with the new method.
2. Enter a **Sequence Name, Group** name and select your desired **Method**.
In this example, the method is Tutorial.



The screenshot shows the 'Chromera Sequence Wizard' dialog box, specifically the 'Sequence Identification' step. On the left is a vertical sidebar with buttons for 'Identification', 'AutoSampler', 'Blanks', 'Calibration', and 'Samples'. The 'Identification' button is highlighted. The main area has a title 'Sequence Identification' and a subtitle 'Enter identification for the sequence and select the method to be used.' Below this are four input fields: 'Sequence name *' with the text 'Chromera Tutorial', 'Description' (empty), 'Group *' with a dropdown menu showing 'Chromera Tutorial', and 'Method *' with a dropdown menu showing 'Tutorial' and a folder icon to its right. A note at the bottom states '* Indicates an entry is required in this field'. At the bottom are three buttons: '< Back', 'Next >', and 'Cancel'.

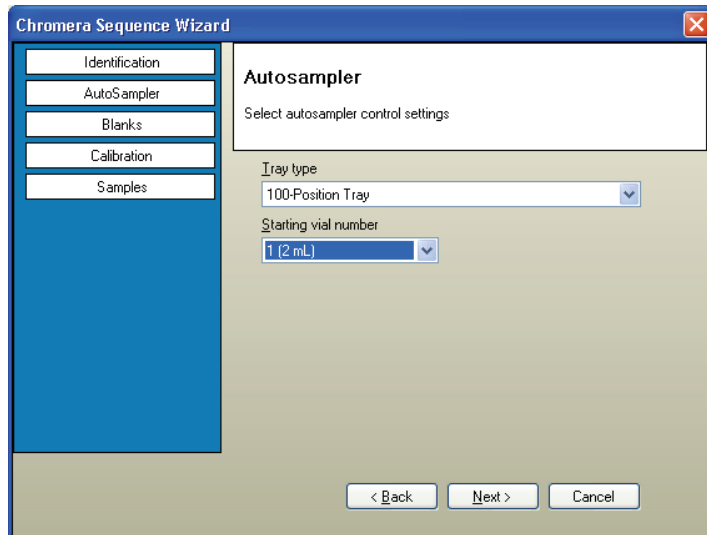3. Click **Next**.
The **Autosampler** screen displays.




The screenshot shows the 'Chromera Sequence Wizard' dialog box, specifically the 'Autosampler' step. The left sidebar is the same, but the 'AutoSampler' button is now highlighted. The main area has a title 'Autosampler' and a subtitle 'Select autosampler control settings'. Below this are two input fields: 'Tray type' with a dropdown menu showing '100-Position Tray', and 'Starting vial number' with a dropdown menu showing '1 (2 mL)'. At the bottom are three buttons: '< Back', 'Next >', and 'Cancel'.

The **Tray Type** selection is restricted to those available for the currently configured autosampler.

4. Click **Next**.

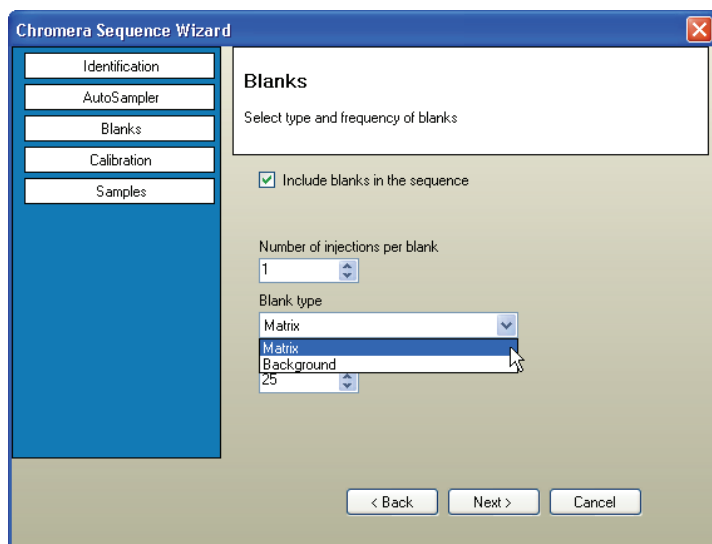
The Blanks screen displays.



The screenshot shows the 'Blanks' screen of the 'Chromera Sequence Wizard'. On the left is a vertical navigation pane with buttons for 'Identification', 'AutoSampler', 'Blanks', 'Calibration', and 'Samples'. The 'Blanks' button is highlighted. The main area has a title bar 'Blanks' and a subtitle 'Select type and frequency of blanks'. It contains a checked checkbox 'Include blanks in the sequence', a spinner for 'Number of injections per blank' set to 1, a dropdown for 'Blank type' with 'Matrix' selected, and a spinner for 'Z5' set to 25. At the bottom are '< Back', 'Next >', and 'Cancel' buttons.

5. Since you will not be using a Blank, click **Next**.

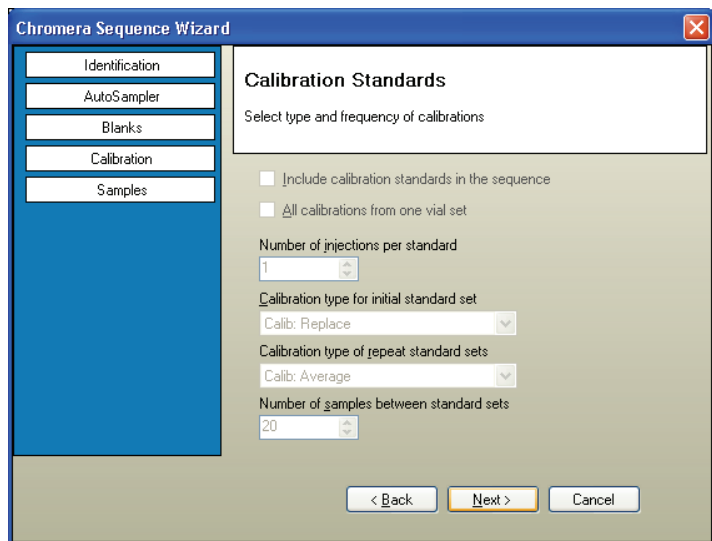The **Calibration Standards** screen displays.



The screenshot shows the 'Calibration Standards' screen of the 'Chromera Sequence Wizard'. The left navigation pane is the same, but the 'Calibration' button is highlighted. The main area has a title bar 'Calibration Standards' and a subtitle 'Select type and frequency of calibrations'. It contains two unchecked checkboxes: 'Include calibration standards in the sequence' and 'All calibrations from one vial set'. Below these are a spinner for 'Number of injections per standard' set to 1, a dropdown for 'Calibration type for initial standard set' with 'Calib: Replace' selected, a dropdown for 'Calibration type of repeat standard sets' with 'Calib: Average' selected, and a spinner for 'Number of samples between standard sets' set to 20. At the bottom are '< Back', 'Next >', and 'Cancel' buttons.

6. Since you will not be using a calibration standard, click **Next**.

The **Samples** screen displays.

You can apply a naming convention to sample naming if desired. If you entered a **Prefix**: Tutorial; **Leading Zeroes**: 2; **Number**: 1 and a **Suffix**: PKI you would create sample names in the series of:

- Tutorial001PKI
- Tutorial002PKI
- Tutorial003PKI
- etc.

7. Be mindful to select to **View sequence** in Editor or in Run Time. Otherwise, it will be created, and saved in the database, but you won't see it until you choose to explicitly open it.

This example suggests to View the sequence **in Run Time**.

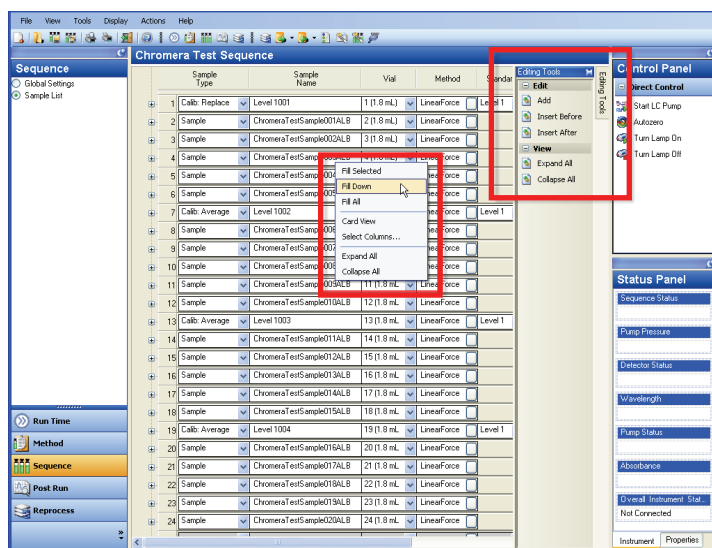
8. When done click **Finish**.

Editing Your Sequence

When sequences are created through the **Sequence Wizard**, they can be automatically opened in the **Sequence Editor** when the wizard completes. Also, if you prefer, new sequences can be created directly from within the Sequence Editor, and existing sequences can be modified.


The editing functions here are the same as those in the **Sequence** mode of **Run Time**, except that the sequence shown in this example is an Active sequence.

The Sequence Editor has a fly-out **Editing Tools** menu containing editing commands.



Decision Points can be set between each row of sequence. They allow the sequence to be paused for a user to perform manual data review:

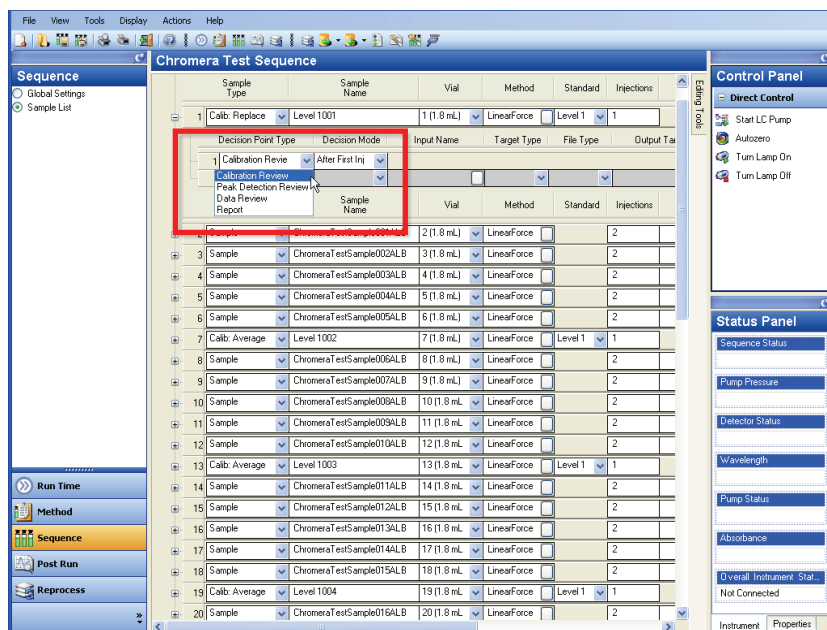
- Review calibration injections and data.
- Inspect and adjust integration parameters for the active method.






Running a Sequence

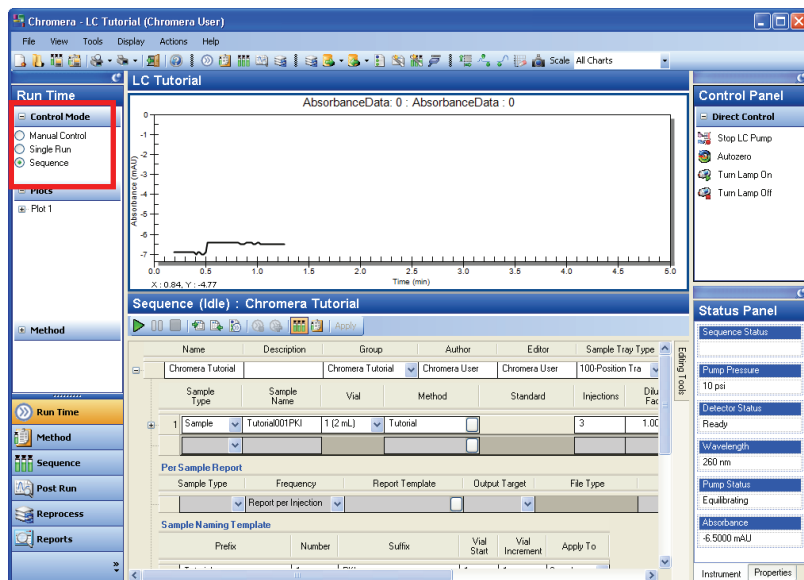
The **Active Sequence** currently controls the instrument, and can be edited or modified as required while it is running. Any row not yet completed can be edited and the changes will be applied, simply by clicking the **Apply** button. Starting and Ending rows for acquisition can be selected from the sequence. Rows that are de-selected will be greyed-out.

1. Click the **Run Time** button, then select the **Sequence Control Mode**.




2. You can use the controls on the Sequence toolbar (mid screen) to **Start** , **Pause** , or **Stop**  the sequence as desired.

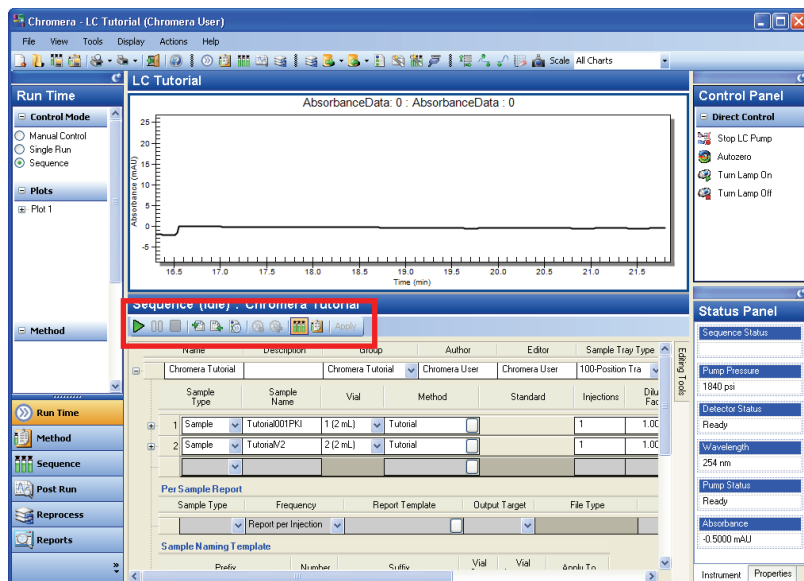You can also select the **Display the Method**  icon to view and/or modify the **Active Method** currently associated with the sequence.




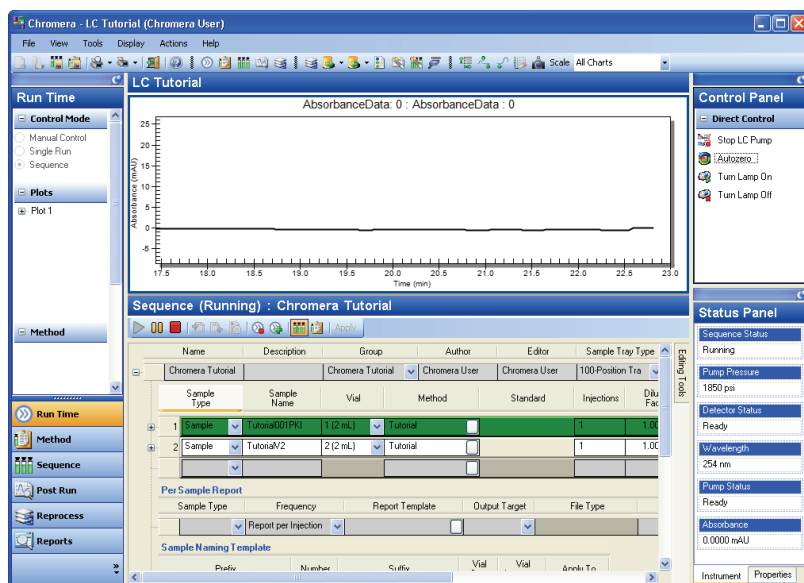3. Click **Start**  to start the sequence.



The sequence starts **Running** and the current row is highlighted green.

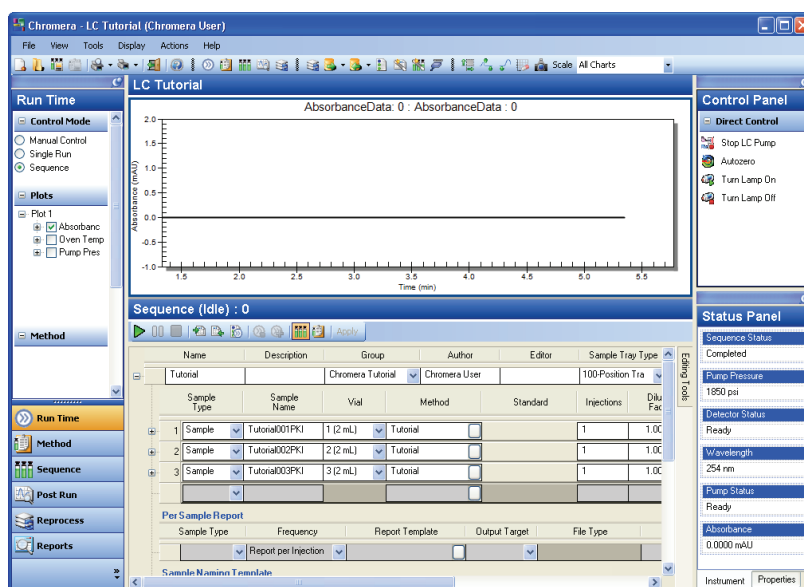- When complete the sequence stops and remains **Idle**.




Congratulations!
You have successfully completed the Chromera Tutorial.

This tutorial was designed to give new users a brief overview of how to setup and explore the main feature of the Chromera software. There are many different features and options available in Chromera that were not discussed in this tutorial. We strongly recommend that all users read the Chromera User's Guide and the on-line Help carefully in order to fully understand this application and how it works.

New users should consider attending a Chromera training class to supplement this tutorial. To find out more information about a training class in your area, please contact 1-800-762-4000 Option 4 or visit PerkinElmer on the web at <http://www.perkinelmer.com/Training>.

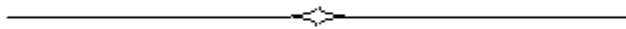What If I Have Problems?



Try This First

Troubleshooting system and software problems can be a real challenge to users who have *extensive* experience in computer systems, networks, and communications. For those of us who have much less experience, and just want to get our samples out the door, the complexity and proliferation of computerized data handling systems and software may seem a bit overwhelming.

The basics of what to do when a problem occurs are pretty much the same. Approaching a system problem methodically BEFORE contacting software support or service personnel can save you time and money. For example, start by answering the following questions:

➤ **Does your system at least meet the MINIMUM requirements as stated in the Chromera Release Notes?**

If it does not, you WILL have problems. These truly are the basic requirements for the version of software.

➤ **Did you check the documentation?**

Check the section of your documentation that explains the commands and procedures you are using to make sure you are operating the software correctly. It may be that what you are trying to do simply cannot be done, or perhaps the feature you are thinking of is actually located in another module of the software.

Do not forget the HELP files (by pressing the F1 key) right at your fingertips, as well as the Release Notes and Known Bugs files available from your Chromera program group. These files contain last minute information we received from our beta testers after manuals were printed. Your question and its resolution may be right here.

➤ **Have you performed good housekeeping procedures on your system?**

Users with MIS staffs overseeing their system needs probably have established system housekeeping procedures. However, if you also pull double duty as the keeper of your computer, you will need to perform weekly maintenance. This includes disk defragmentation, file backup, and removal of any files you no longer need (for example; error log files). Data acquisition systems such as Chromera are incredibly disk intensive, so the cleaner your hard drive, the faster your system.

Also check your system's TEMP directory to make sure it's not cluttered with .TMP files you no longer need. Use a general diagnostics utility once each month to check all other components of the system for signs of failure. Do not forget to check those cables as well to make sure they are seated correctly and are not worn from moving equipment around. A little preventative action now will help ensure that your system is in top shape.

➤ **Are you encountering specific error messages?**

"I think it said something like . . ." may not be much help to people trying to identify what is wrong. Determine whether those errors are coming from your computer, DOS, Windows, or the chromatography software, and write down or make a screen print of exactly what the message says. Each of these error messages should be described in the appropriate manuals supplied with your system or software; this description will help direct you to the problem.

Chromera has a **log** directory that stores error messages encountered by Chromera. Your Chromera Support Specialists may request that you send these log files for further analysis. We recommend that you keep at least a week's worth of these files on hand. This will help us to determine if the system has a history of repeated error messages.

➤ **Is the problem repeatable, or can it be duplicated on another system?**

Fundamental to solving software and hardware mysteries is knowing exactly when they occur and what the conditions are that lead to their occurrence.

- a) If you have multiple Chromera systems available, check whether the problem occurs on other systems as well. If it does not, this may point to a system problem and will help focus the troubleshooting effort.
- b) Can the problem be duplicated with the test data files that came with the system? If not, this may point to a problem with the new files you created.
- c) Are you on a network? See if the problem still exists when you are disconnected from the network to rule out a network problem.

Remember that system operation is dependent on everything communicating in harmony. In the case of data acquisition software, the timing and number of these communications are incredibly complex and demanding on the system, involving both internal (computer, disk drive, printer) and external (instrument control).

Occasionally, something interferes with this communication and simply cannot be traced. Perhaps it is the result of electrical disturbances in your lines from voltage brownouts or surges; a strong magnetic field created by equipment in the next lab; a loose connection; a failing piece of hardware; or just the night cleaning crew plugging their vacuum cleaners into the same circuit as the system analyzing your overnight run. If the problem seems to have disappeared, make a note of it just in case, but you can generally proceed without concern.

➤ **What has changed since the problem started?**

"It was working fine yesterday. No one touched anything...really!"

Software doesn't just break. Ask everyone who uses or maintains your system whether anything has been done to it that you do not know about.

- Has someone installed a new software package or driver that might be interfering?
- Has someone on the night shift altered the Chromera or Windows configuration file since you last operated the system?
- Has someone changed my directory permissions?
- Have my user and/or Chromera process account passwords expired?

It is not impossible that a software file gets corrupted somehow, but before you reinstall the software, check out the obvious culprits first. These can often be diagnosed quickly by rebooting your system with a blank system disk containing basic autoexec and config files. If the problem disappears with the use of this system disk, look to the startup files for a potential conflict.

Hardware on the other hand does break. Even worse, it can break with an insidious style that drives users and support specialists crazy. It can work fine one day and create havoc the next. Often we cannot catch this until it finally breaks for good. Every system should have a good diagnostic software tool that puts its hardware through the paces so hardware faults can be determined or ruled out. Do not forget the obvious here, or maybe not so obvious, too.

- Did someone move equipment? If so are the cables seated properly?
- Is your hard disk full?
- Is your printer connected and turned on?
- Did someone wrap that extra length of cable around the base of an old centrifuge to clean up the lab counter top?

More than half of all the calls to Software Support have nothing to do with the actual Chromera software, but are the result of problems with Windows, network configuration, hardware system configuration, or device drivers for printers and other devices.

Internet Support Help

Additional support information is provided on the PerkinElmer web site:

<http://las.perkinelmer.com/content/onesource/index.html>

This service provides an economical vehicle for users to get the latest information regarding product support issues.

When All Else Fails and You Need Help

If you tried all of the above and still need help or advice, it is time to contact us. PerkinElmer has several options in place for assistance with PerkinElmer software and hardware. Your local sales or service representative will provide you with information regarding how to access providers of assistance in your particular country.

If your shipment contains a companion pamphlet covering software support options in your country, review that pamphlet to see what options are available to you BEFORE that first question or problem comes up. Keep this information along with user manuals, software media, and the license agreement in a safe, accessible location.

Register for a User's Class

Most labs will have a few users that will quickly understand Chromera. As the first month of Chromera implementation nears completion, arrange to have one or more users participate in a product-training course. Do not sign up for the class before you have at least 30 days on the system. Our technical trainers have seen new users take the class before the installation of Chromera, and they just do not have the experience with the software to get everything out of the class they would have if they had waited.

In the U.S., you can request a Technical Training brochure, current schedule, or class registration by calling the For more training contact PerkinElmer at 1-800-762-4000 Option 4 to hear the available training options or visit PerkinElmer on the web at <http://www.perkinelmer.com/Training>. Outside the U.S., contact your PerkinElmer representative for local information.

PerkinElmer Product Specialists

The PerkinElmer field organization has a staff of highly skilled local specialists in GC, LC, and data handling products who are available for on-site assistance with any application need your lab might have. These specialists draw on many years of experience in laboratory environments just like yours. If you would like to use these resources, feel free to contact your PerkinElmer sales staff. In the U.S., call the Sales Administration staff at 1-800-762-4000.

Placing the Call . . .

When you need to talk to a technical representative, the following list can help you get the most from that conversation. Like you, we want to resolve your call as quickly as possible.

- **Have the phone within reach of the system.**
We may have you try several things on your system to determine the source of your problem. Or, we may have you go over the steps you are having problems with, keystroke by keystroke, as we duplicate them on the other side.
- **Have your Serial Number, Software Licenses and/or your original software media available.**
You will need these to access most centralized phone support programs and services.
- **Have your computer hardware, Windows, and software manuals available.**
We may need to refer you to these for additional information.
- **Know the version and service packs of all relevant software packages.**
- **Record any error messages exactly as they appear on the system.**
Make sure you can state exactly where errors occurred and what steps are required to reproduce them.
- **Have available the make and model numbers of your computer and associated hardware.**
- **Have the end user make the call.**

NOTE: ***DO NOT** make calls on behalf of the end user having the problem unless you fully understand both the problem and the system. These calls are seldom resolved quickly because of the unfamiliarity of the caller with the situation.*

Contacting PerkinElmer

	US & Canada Customers	Outside the US
Customer Support	Phone: (800) 762-4000	Local PerkinElmer Sales Office

What Happens to My Problem or Suggestion?

As problems and suggestions regarding software operation, functionality, features, or documentation are reported by registered users, Software Support begins a software investigation. This investigation starts with a search for whether the issue has been previously reported. If it has, you will be notified of the resolution. If a problem has not been reported, Software Support will attempt to duplicate your experience on standard supported computer systems sold by PerkinElmer.

If the problem can be duplicated by Software Support, and it resides in our software, a Software Problem Report will be prepared and forwarded to the Product Department, as well as, Software Development for further action.

What If My Problem Can't Be Resolved Immediately?

Occasionally, everything we do in a remote troubleshooting process comes up empty, and we need to do some research with the product, your system and data, our engineering staff, or a third party. If we cannot reproduce the problem, we will assist you with as much user education, system configuration troubleshooting, and hardware troubleshooting as is possible from a remote location. We will continue until the problem is resolved or until we can advise you what we think is necessary for it to be resolved.


In cases where we cannot make the symptoms appear as described, we may need to get additional information from you regarding the system you are working on and its setup and configuration. We may also need verification that the software is installed according to product documentation. This process may require eliminating potentially interfering third party software or hardware in an effort to return the product to a functioning state. Sometimes this can be a tedious process over the phone, and as Murphy's Law dictates, this always seems to occur when you have to get those 200 samples out the door. Please be patient. You are our eyes, ears, and hands in these situations, and we need you.

If the problem requires on-site assistance, or is the result of operation of the software on non-standard or non-supported hardware or system configurations, Software Support may make several suggestions. One may be to request you involve the local PerkinElmer Service office and get a service engineer on site who can help us evaluate the situation first hand. Software Support will work with the local office to assist in the effort.



Appendix 1: Importing TotalChrom Data and Methods

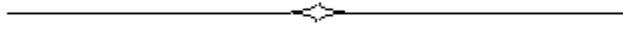(Advanced Task)




How to Import TotalChrom Data and Methods

This exercise will demonstrate how to import TotalChrom data (TotalChrom Index Files Data, TotalChrom raw data, or TotalChrom component names from TotalChrom methods) provided that the RAW data files are in TotalChrom 6.3.x format.


The component names are read into the Chromera **Component Dictionary**. The **Component Dictionary** is a typing aid used in the Chromera **Method Editor**, to allow component names to be assigned to peaks in a new method by selecting them from the dictionary, or by using the “auto complete” mode when typing a component name. This functions much like Microsoft Outlook when typing the name of a recipient into the **To:** or **Cc:** list for a new message.

To Import TotalChrom Data and Methods:

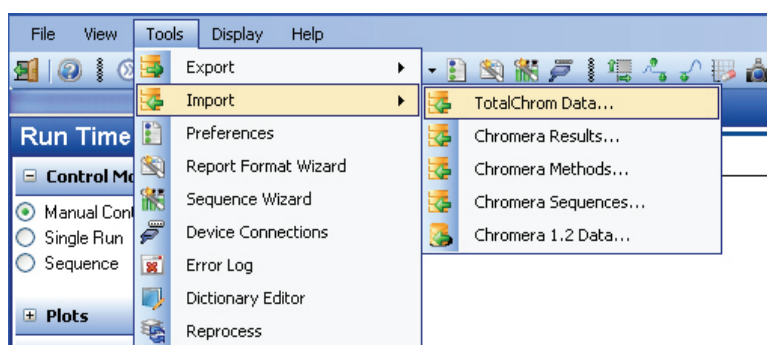- From the **Tools** menu, select **TotalChrom Data** from the **Import** Menu.



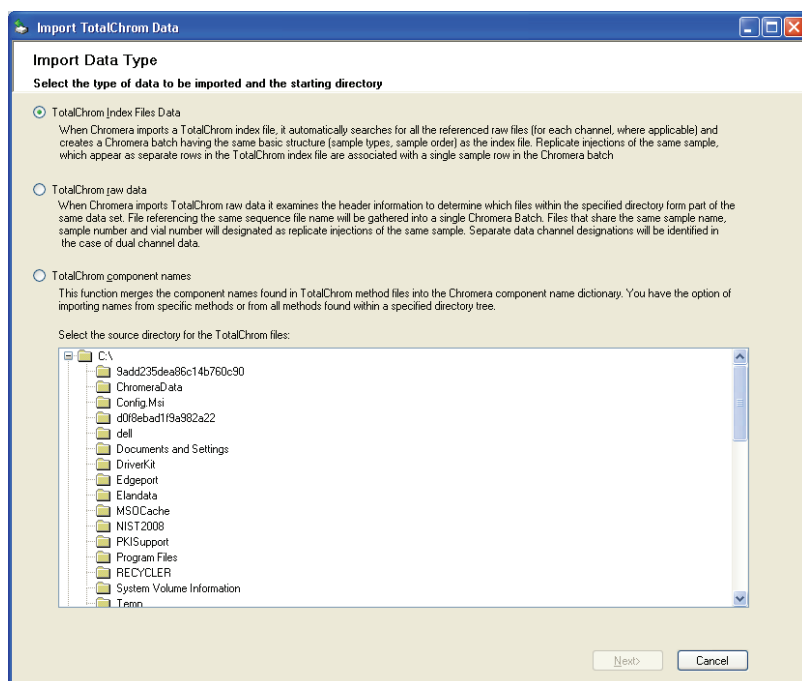The **Import Data Type** screen appears:




This initial page displays for you to select the type of TotalChrom data to be imported (TotalChrom Index Files Data, TotalChrom raw data, or TotalChrom component names) and select the source directory containing the TotalChrom files.

Import TotalChrom Index Data Files

This option enables you to import a group of TotalChrom RAW data files by selecting the Index (IDX) file that points to them. This creates a corresponding Chromera Batch as a consequence.

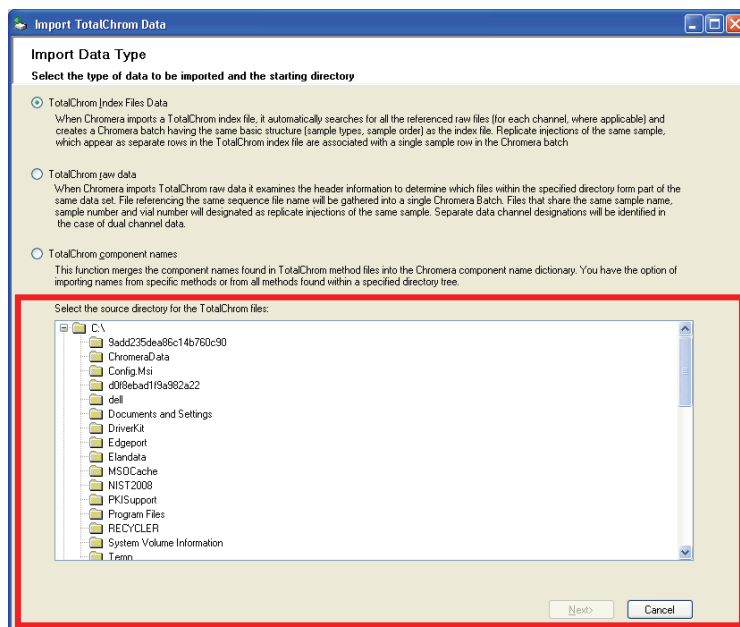1. Select **TotalChrom Index Files Data**.




2. Browse to a folder containing the TotalChrom Index files then click the **Next** button.

Only the TotalChrom files in that specific folder will be displayed for importing. Subfolders of the current one will not be searched automatically.


NOTE: The **Next** button is not enabled until you select a directory containing files corresponding to the selected data type. For example, if you select **TotalChrom Index Files Data**, the directory must contain **IDX** files.

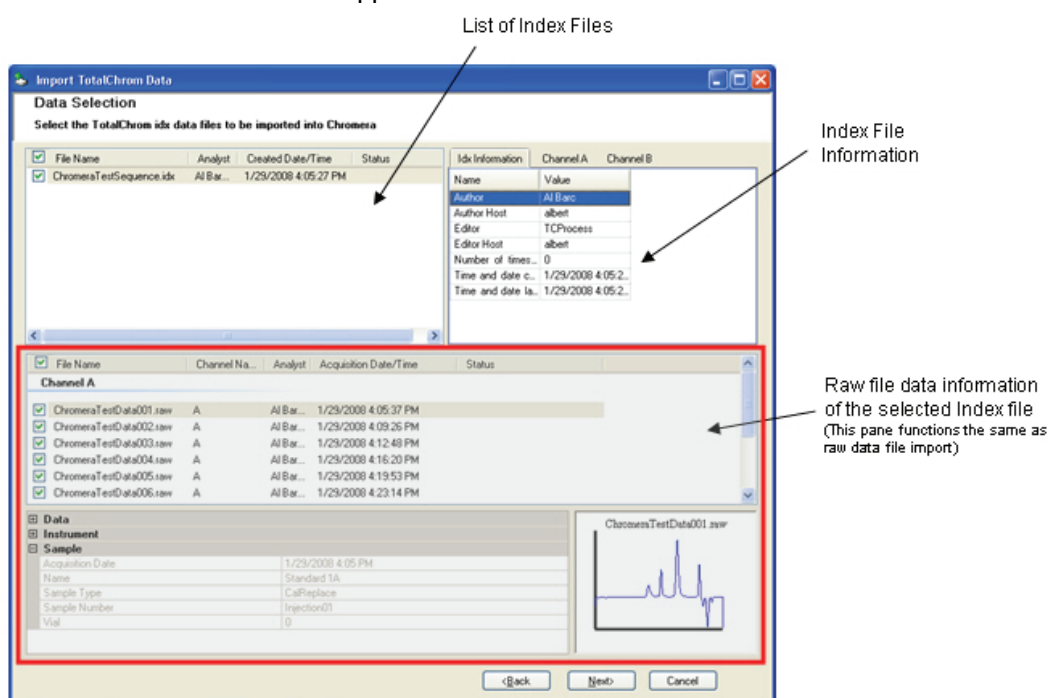The **Data Selection** screen appears:



The software reads each file in the selected directory with an .idx extension and determines if it is a valid TotalChrom IDX file. For the purposes of import into Chromera, "valid" means that its format is the latest TotalChrom IDX file revision. If the file is valid it is included in the upper pane display list. Also, a list of raw files for each idx is displayed in the lower pane (The view is same as seen in raw file import). The first row is selected and highlighted in the IDX file pane, and the first raw file is highlighted with its information and chromatogram displayed in the bottom two panes.

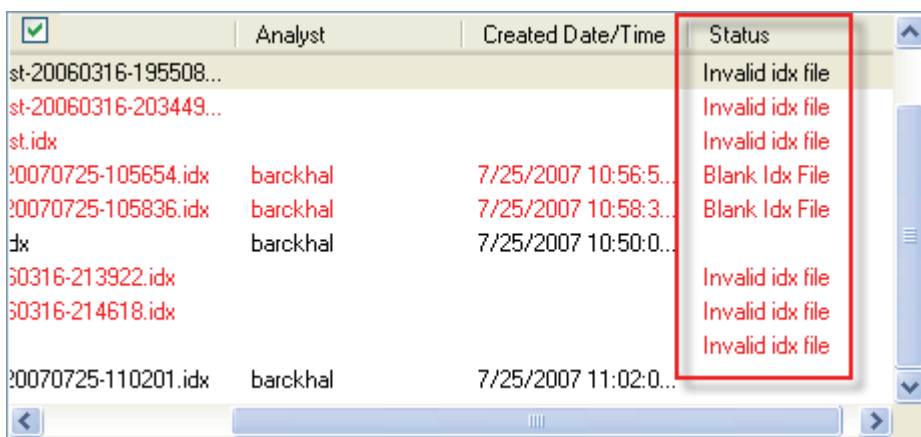If IDX files are highlighted in red it indicates a problem shown in the **Status** field ... which may be scrolled out of sight.




<input checked="" type="checkbox"/>	Analyst	Created Date/Time	Status
<input checked="" type="checkbox"/>			Invalid idx file
<input checked="" type="checkbox"/>			Invalid idx file
<input checked="" type="checkbox"/>			Invalid idx file
<input checked="" type="checkbox"/>	barckhal	7/25/2007 10:56:5...	Blank Idx File
<input checked="" type="checkbox"/>	barckhal	7/25/2007 10:58:3...	Blank Idx File
<input checked="" type="checkbox"/>	barckhal	7/25/2007 10:50:0...	
<input checked="" type="checkbox"/>			Invalid idx file
<input checked="" type="checkbox"/>			Invalid idx file
<input checked="" type="checkbox"/>			Invalid idx file
<input checked="" type="checkbox"/>	barckhal	7/25/2007 11:02:0...	

3. You can select one or more IDX files to read. All valid ones will be selected by default.
4. You can select individual raw files from total available raw files for a select index file, by clicking in each check box for the desired raw files.

NOTE: You can also select all / deselect all files of either type by clicking in the check box next to **File Name** in either the IDX or RAW file list panes. This selects / deselects all valid index files and all corresponding valid raw files. Also, selecting any raw file, and then clicking the **Next** button, will automatically select its corresponding IDX file.

5. After selecting your desired Index and Raw files, click **Next**.

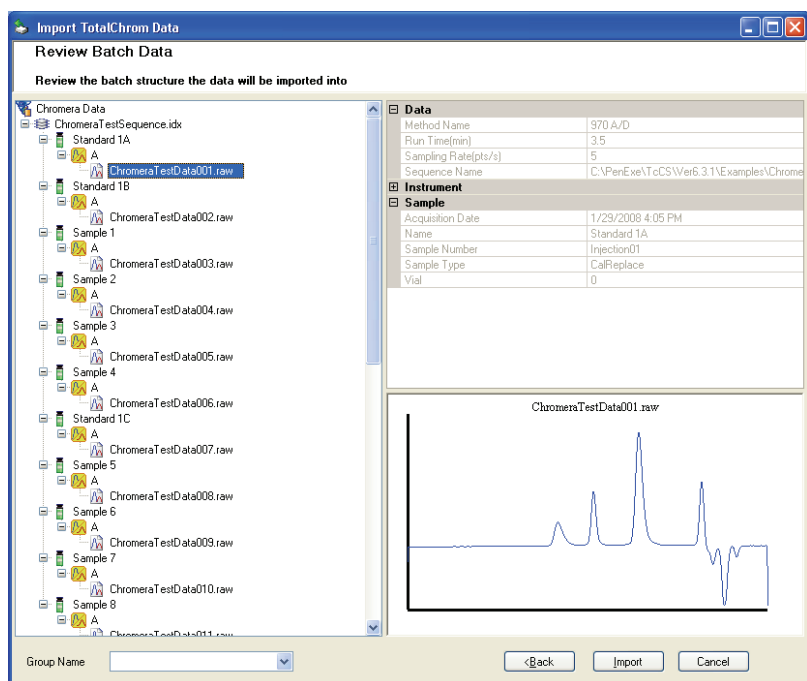The **Review Batch Data** screen appears, displaying the contents of the IDX file in the logical layout of the Batch that will be created.:



NOTE: When reviewing batch data from an index file, you can rename the batch name (for example, ChromeraTestSequence.idx) by right-clicking on it, selecting **Rename Node**, and typing a new name.

6. Select from the tree display the individual files to import, or select the batch name (the .idx file name) to import all the data identified in that IDX file.

Clicking on an entry in the list displays the chromatogram and sample information for that file.

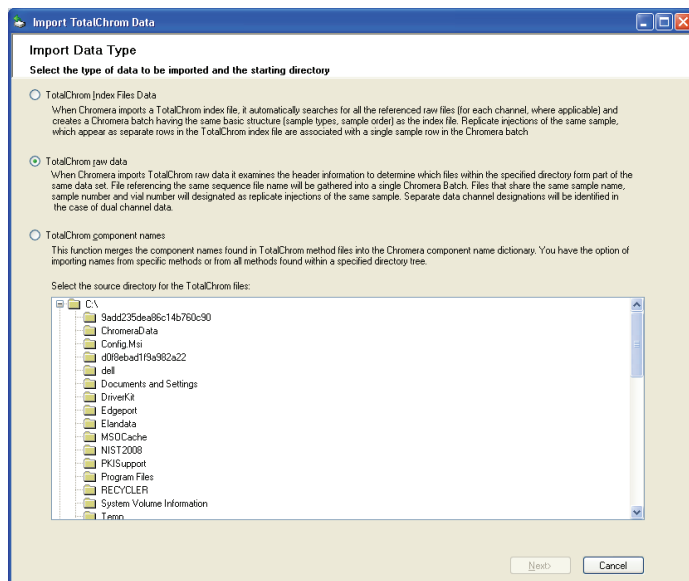
NOTE: The data files pointed to by the IDX file must be accessible to this computer if they are stored elsewhere, such as on a network drive or file server.

7. Type a **Group Name** then click the **Import** button to load the selected data into the Chromera database as a new **Batch**, associated with the **Group Name** entered.

Import TotalChrom raw data Files

This option enables you to import a group of TotalChrom RAW data files by selecting them directly from a folder. A corresponding Chromera Batch will be created as a consequence.

1. Select **TotalChrom raw data** Files.




2. Browse to a folder containing the TotalChrom raw data files then click the **Next** button.

Only the TotalChrom files in that specific folder will be displayed for importing. Subfolders of the current one will not be searched automatically.

NOTE: The **Next** button will not be enabled until you select a directory containing files corresponding to the selected data type. For example, if you select **TotalChrom raw data**, the directory must contain raw files.

The software reads each file in the selected directory with a .raw extension and determines if it is a valid TotalChrom raw file. For the purposes of import into Chromera, "valid" means that its format is the latest TotalChrom raw file revision. If the file is valid it is included in the display list.

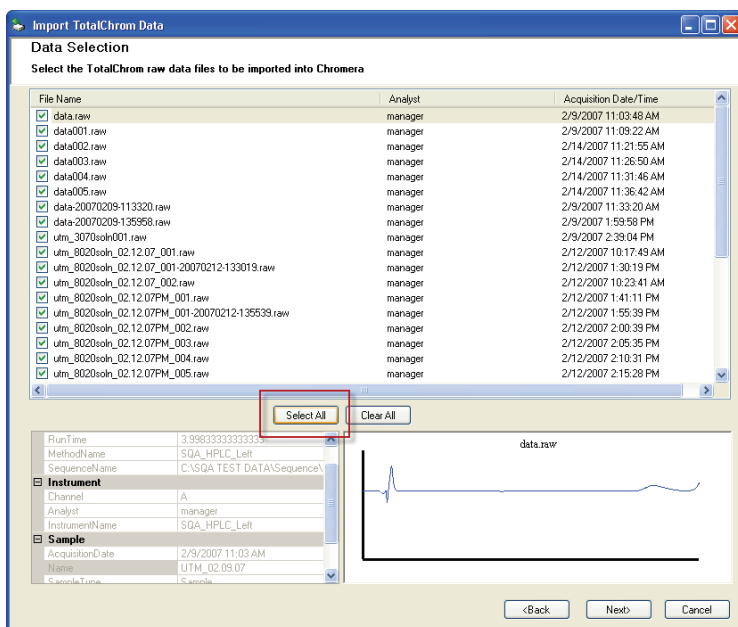
A progress bar on the bottom of the screen displays while files are examined. The progress is updated based on the percentage of raw files contained in the directory that have been examined. The **Select All**, **Clear All** and **Back** buttons are disabled during this process. When complete the **Data Selection** screen appears:



Select the individual files to import, or click **Select All** to import all the data shown in the list. Only the TotalChrom RAW files from the specific folder will be displayed for importing. Subfolders of the current one will not be searched automatically. Clicking on an entry in the list displays the chromatogram and sample information for that file.

3. Click the **Next** button.

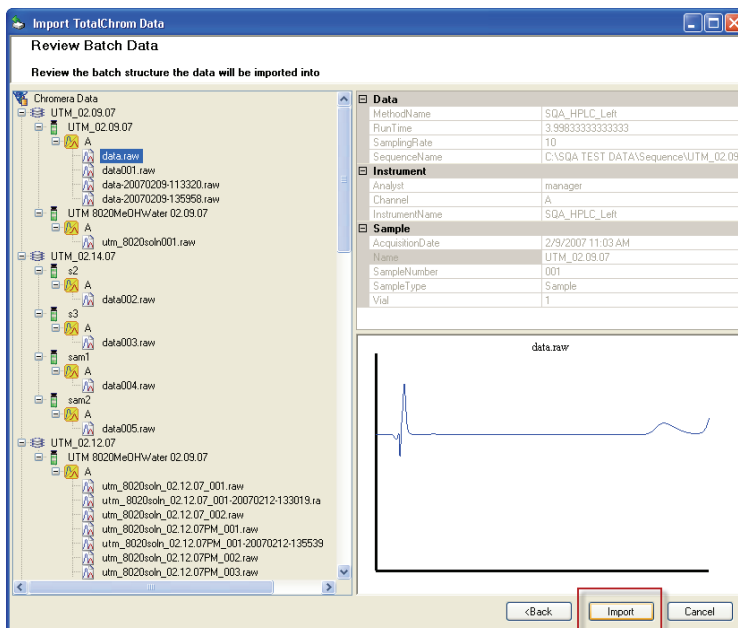The next step inspects the data to see which are valid for importing. Only data from TotalChrom 6.3.x or newer can be imported directly. Any older files must be converted to a newer version through TotalChrom Batch Reprocess using **"Update existing RAW file header"** option.




Any TotalChrom data not of a supported version, or otherwise not suitable for import, will be highlighted in red in the **Batch Data** list.

4. Type a **Group Name** the click the **Import** button to load the selected data into the Chromera database as a new **Batch**, with the name of the TotalChrom sequence from which the RAW data were originally acquired.

Import TotalChrom component names

This option enables you to populate the Chromera **Component Dictionary** with the complete list of component names read from all the TotalChrom methods found in a selected folder and subfolders.

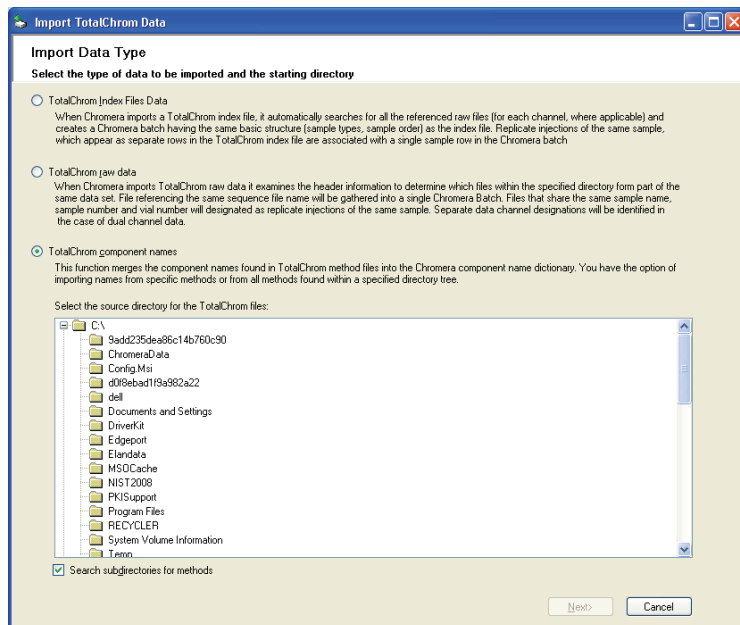1. Select TotalChrom component names.




To import component names from TotalChrom methods, browse to a folder that contains the MTH files to be read. Individual folders can be searched, or all the subfolders below the selected one can be searched as well by selecting **Search subdirectories for methods**.

2. Browse to a folder containing the TotalChrom methods then click the **Next** button.

NOTE: The **Next** button will not be enabled until you select a directory containing files corresponding to the selected data type. For example, if you select TotalChrom component names, the directory must contain method files.

The software reads each file in the selected directory with a .mth extension and determines if it is a valid TotalChrom method file. For the purposes of import into Chromera, "valid" means that its format is the latest TotalChrom method file revision. If the file is valid it is included in the display list.

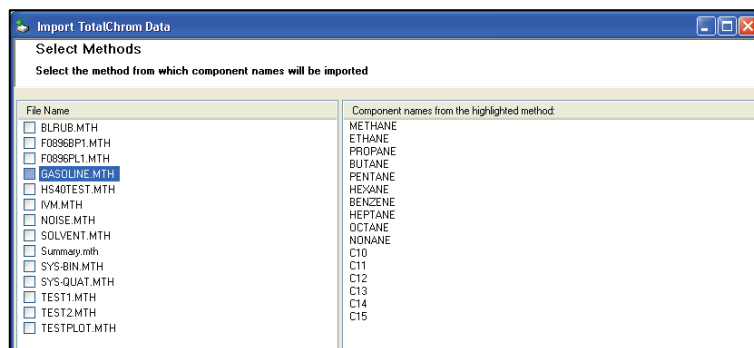3. Click on the method name to highlight it and display the method's components in the right pane.




Highlight each method to see if it contains the components you want.

4. If the method contains your desired components, click in a check box to select a method.

You can continue to highlight methods to see components and select methods that contain the components you want.

5. After selecting the methods that contain your desired components, click **Next**.

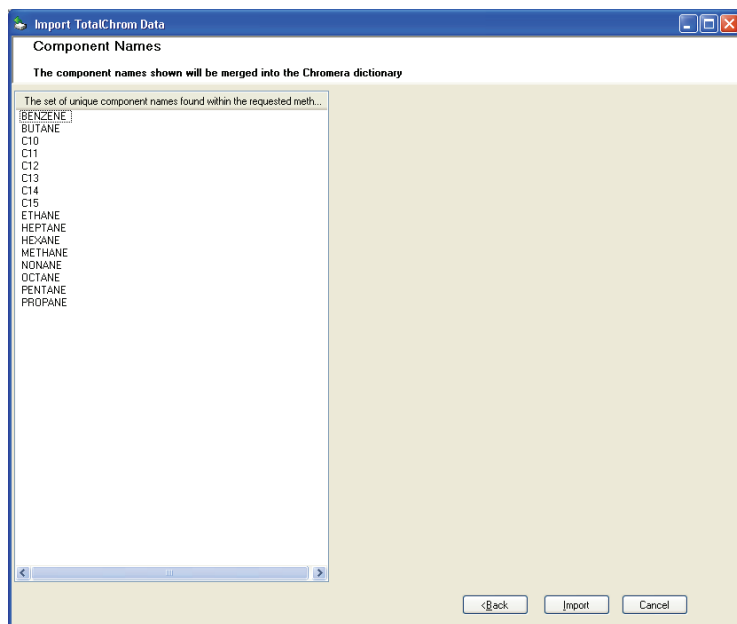A list of components from all of the selected methods are displayed to be merged into the Chromera directory.




6. Click the **Import** button to load the component names into the **Component Dictionary**.


Copying Component Information in Chromera Method Editor

It is also possible to copy detailed component information from an individual TotalChrom method into the Chromera method editor, to populate the **Peaks** and **Calibration** sections of a Chromera method.

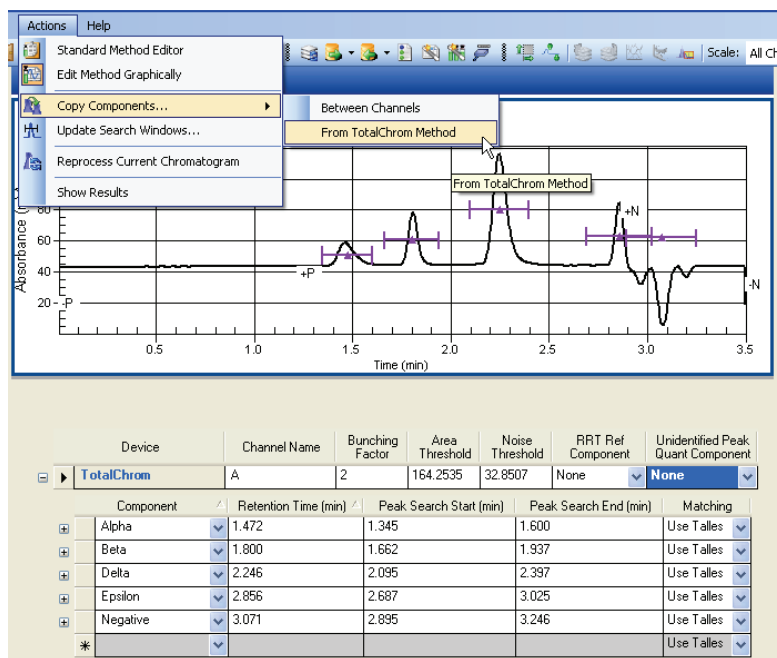1. Select the **Actions** menu. Then select **Copy Components > From TotalChrom Method**.



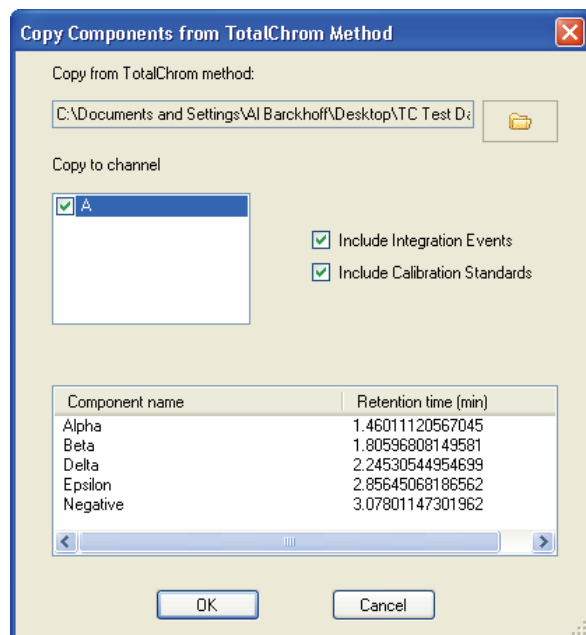This function will copy not just component names, but also retention times, calibration information, and integration timed events to a selected Chromera data channel.



2. In **Copy from TotalChrom method** browse for a method.

Displays the source of the components to be copied. The field is populated by selecting a file from a standard File Open dialog displayed when the browse button is clicked.

3. In **Copy to channel** select the checkbox to set the Chromera method data channel(s) the components will be copied to. The list is populated with all data channels from the method. You can copy the components to more than one channel by selecting the required channels within this control.
 - Select to **Include Integration Events** to copy the integration events from the TotalChrom method to the destination channel.
 - Select **Include Calibration Standards** to copy the calibration standard amounts with the other component information.
4. The display at the bottom lists the components found in the TotalChrom method, together with their retention times.
5. Click **OK** to copy the components.





PerkinElmer, Inc.

710 Bridgeport Avenue
Shelton, CT 06484-4794, U.S.A.

Internet: <http://www.perkinelmer.com>
email: info@perkinelmer.com

PerkinElmer is a registered trademark of PerkinElmer, Inc.