




ANTIARRHYTHMICS HPLC

Pharmaceutical: Antiarrhythmic Verapamil HCl Raw Material (EP 2002 modified) on Brownlee Validated IBD

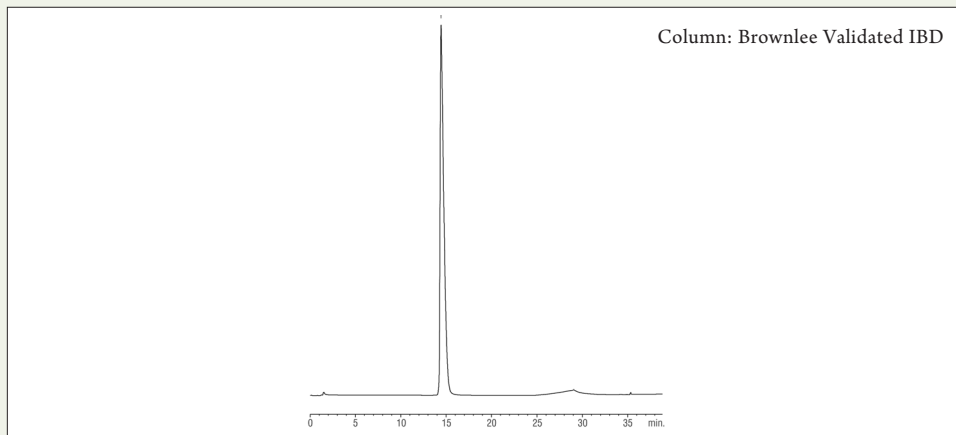

Peak List

Peak List	Retention Time (min)
1. verapamil	15.0

Column

Column name:	Brownlee Validated IBD
Part number:	N9303584
Dimensions:	250 x 4.6 mm
Particle size:	5 μ m
Pore size:	100Å

Conditions

Mobile phase:	A: 40 mM potassium phosphate monobasic, pH 6.25
	B: acetonitrile

Time (min)	%B
0.00	37
22.0	37
27.0	65
35.0	65
36.0	37
46.0	37

Flow:	1.5 mL/min
Temperature:	25 °C
Detector:	UV @ 278 nm

Sample

Injection volume:	20 μ L
Concentration:	25 μ g/mL
Solvent:	initial mobile phase

Notes

Recommended Guard Cartridge	N9303702
Recommended Guard Holder	N9303713

PerkinElmer, Inc.
940 Winter Street
Waltham, MA 02451 USA
P: (800) 762-4000 or
(+1) 203-925-4602
www.perkinelmer.com



For a complete listing of our global offices, visit www.perkinelmer.com/ContactUs

Copyright ©2013 PerkinElmer, Inc. All rights reserved. PerkinElmer® is a registered trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.

007573A_22