




AROMATICS

HPLC

Environmental: Aromatics on Brownlee Analytical C18

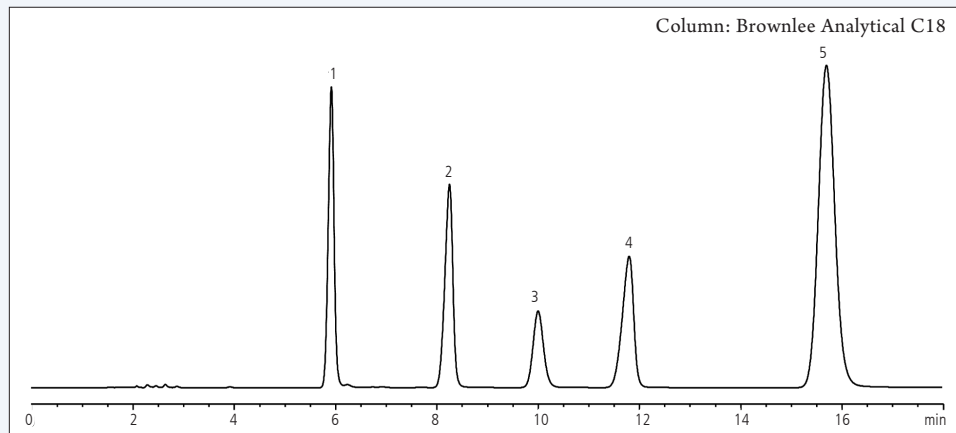

Peak List

Peak List	Concentration (µg/mL)
-----------	-----------------------

1. propylbenzene	15000
2. butylbenzene	15000
3. o-terphenyl	200
4. amylbenzene	15000
5. triphenylene	100

Column

Column name:	Brownlee Analytical C18
Part number:	N9303513
Dimensions:	150 x 4.6 mm
Particle size:	5 µm
Pore size:	110Å

Conditions

Mobile phase:	water:methanol (20:80 v/v)
Flow:	1.0 mL/min
Temperature:	ambient
Detector:	UV @ 254 nm

Sample

Injection volume:	5 µL
Solvent:	methanol

Notes

Recommended Guard Cartridge	N9303677
Recommended Guard Holder	N9303713

PerkinElmer, Inc.
940 Winter Street
Waltham, MA 02451 USA
P: (800) 762-4000 or
(+1) 203-925-4602
www.perkinelmer.com



For a complete listing of our global offices, visit www.perkinelmer.com/ContactUs

Copyright © 2013 PerkinElmer, Inc. All rights reserved. PerkinElmer® is a registered trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.

007571A_05