




DIETARY SUPPLEMENTS/ PLANT EXTRACTS HPLC

Foods and Flavors: Tocopherols in Palm Oil on Brownlee Analytical Silica

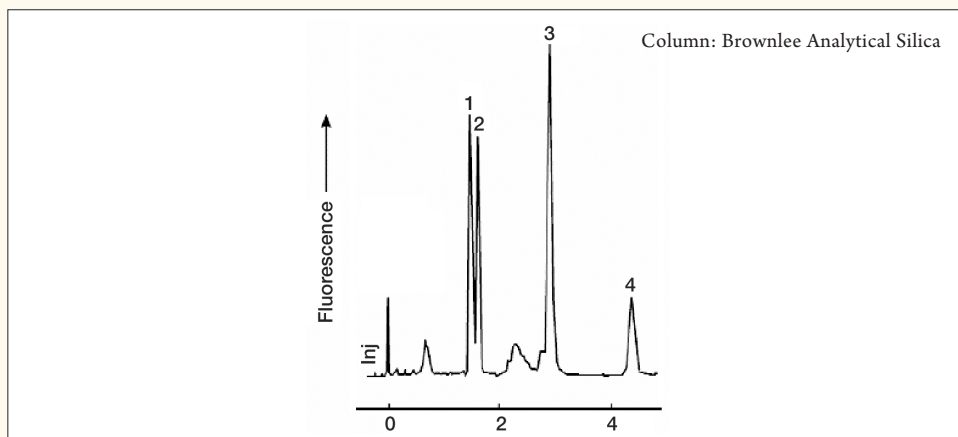

Peak List

1. α -tocopherol
2. α -tocotrienol
3. β -tocotrienol
4. γ -tocotrienol

Column

Column name: Brownlee Analytical Silica
Part number: N9303525
Dimensions: 150 x 4.6 mm
Particle Size: 5 μ m

Conditions

Mobile phase: 0.74% ethanol in hexane
Flow: 2.0 mL/min
Pressure: 900 psi (6.4 MPa)
Detector: FLR @ Ex 298 nm; Em 325 nm

Notes

Recommended Guard Cartridge N9303687
Recommended Guard Holder N9303713

PerkinElmer, Inc.
940 Winter Street
Waltham, MA 02451 USA
P: (800) 762-4000 or
(+1) 203-925-4602
www.perkinelmer.com


PerkinElmer[®]
For the Better

For a complete listing of our global offices, visit www.perkinelmer.com/ContactUs

Copyright ©2013 PerkinElmer, Inc. All rights reserved. PerkinElmer[®] is a registered trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.

007572A_33