




CARBOXYLIC ACID HPLC

Foods and Flavors: Organic Acids on Brownlee Polypore H

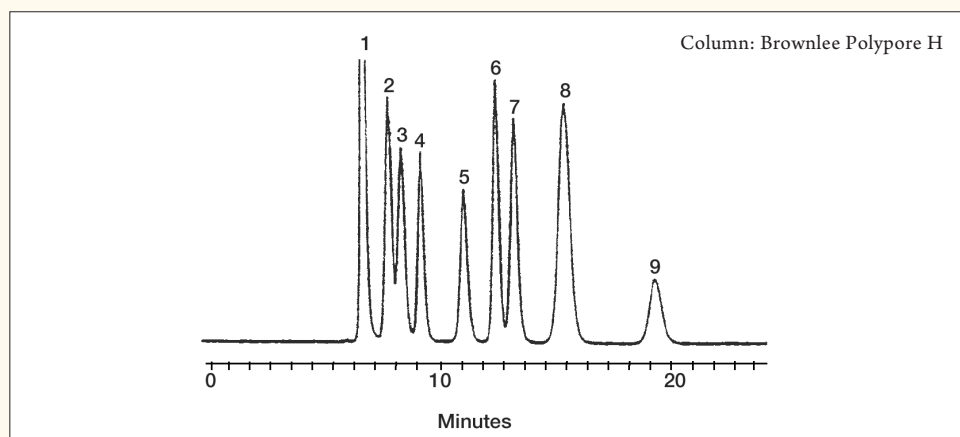

Peak List

1. oxalic acid
2. citric acid
3. α -ketoglutaric acid
4. malic acid
5. succinic acid
6. formic acid
7. acetic acid
8. fumaric acid
9. butyric acid

Conditions

Mobile phase: 0.01 NH_2SO_4
Flow: 0.15 mL/min
Temperature: 25 °C
Detector: UV @ 210 nm

Notes

Recommended Column Holder	07150015
Recommended Guard Cartridge	07110085
Recommended Guard Holder	07150013

Column

Column name: Brownlee Polypore H
Part number: 07110089
Dimensions: 220 x 4.6 mm
Particle Size: 10 μm

PerkinElmer, Inc.
940 Winter Street
Waltham, MA 02451 USA
P: (800) 762-4000 or
(+1) 203-925-4602
www.perkinelmer.com

PerkinElmer
For the Better

For a complete listing of our global offices, visit www.perkinelmer.com/ContactUs

Copyright © 2013 PerkinElmer, Inc. All rights reserved. PerkinElmer® is a registered trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.

007572A_34