




DIETARY SUPPLEMENTS/ PLANT EXTRACTS HPLC

Foods and Flavors: Vitamins A and E on Brownlee Analytical Silica

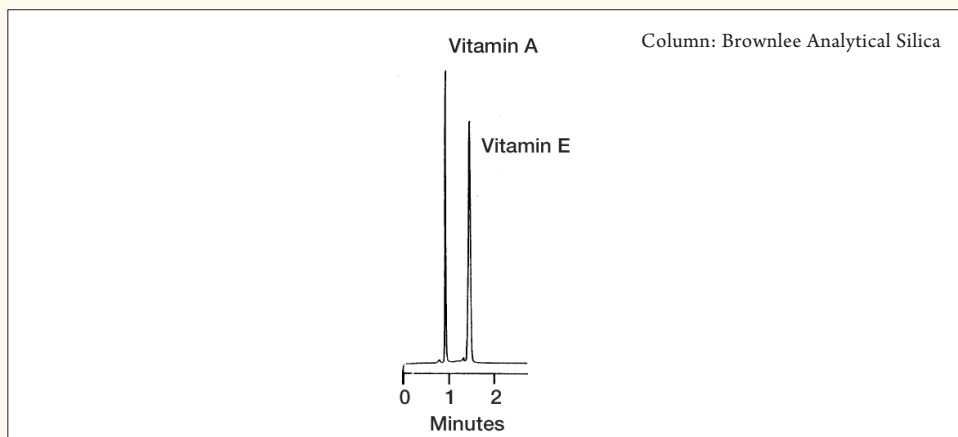

Column

Column name: Brownlee Analytical Silica
Part number: N9303525
Dimensions: 150 x 4.6 mm
Particle Size: 5 μ m

Conditions

Mobile phase: 95% hexane, 5% chloroform
Flow: 2 mL/min
Detector: FLR @
325/480 nm for vitamin A
295/335 nm for vitamin E

Notes

Recommended Guard Cartridge N9303687
Recommended Guard Holder N9303713

PerkinElmer, Inc.
940 Winter Street
Waltham, MA 02451 USA
P: (800) 762-4000 or
(+1) 203-925-4602
www.perkinelmer.com



For a complete listing of our global offices, visit www.perkinelmer.com/ContactUs

Copyright ©2013 PerkinElmer, Inc. All rights reserved. PerkinElmer® is a registered trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.

007572A_36