




PRESERVATIVES HPLC

Foods and Flavors: Preservatives on Brownlee Pecosphere 5 C18

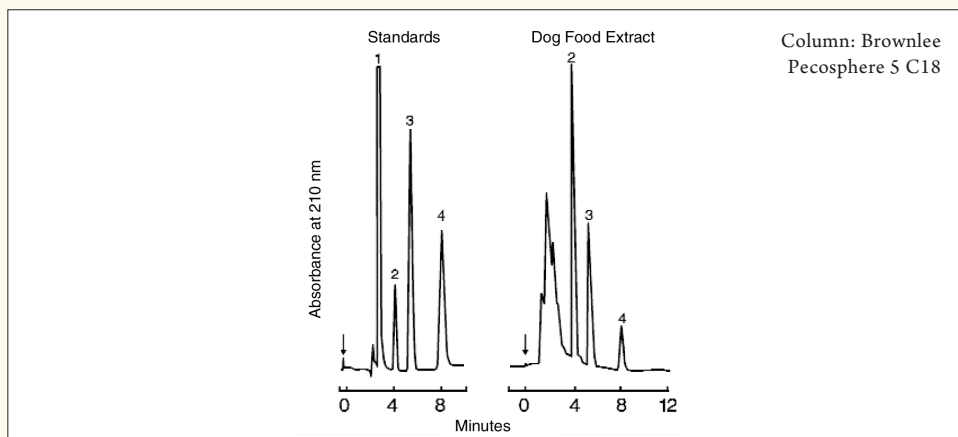

Peak List

1. p-hydroxybenzoic acid
2. sorbic acid
3. p-hydroxybenzoic acid ethyl ester
4. p-hydroxybenzoic acid propyl ester

Column

Column name: Brownlee Pecosphere 5 C18
Part number: 02580169
Dimensions: 150 x 4.6 mm
Particle size: 5 μ m

Conditions

Mobile phase: methanol:water (40:60),
0.2% v/v perchloric acid
Flow: 1.0 mL/min
Detector: UV @ 210 nm

Notes

Recommended Column Holder	07150030
Recommended Guard Cartridge	00402237
Recommended Guard Holder	07150001
Recommended Column Guard Holder	07150030
	07150018
	07150002

PerkinElmer, Inc.
940 Winter Street
Waltham, MA 02451 USA
P: (800) 762-4000 or
(+1) 203-925-4602
www.perkinelmer.com



For a complete listing of our global offices, visit www.perkinelmer.com/ContactUs

Copyright © 2013 PerkinElmer, Inc. All rights reserved. PerkinElmer® is a registered trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.

007572A_44