




TRIGLYCERIDES HPLC

Foods and Flavors: Olive Oil Triglycerides on Two Brownlee Pecosphere 3 C18

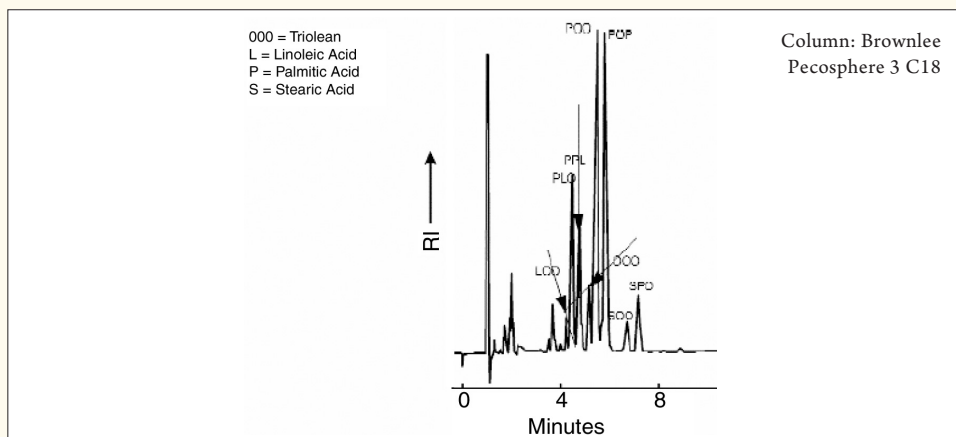

Column

Column name: Brownlee Pecosphere 3 C18
Part number: 02580166
Dimensions: 83 x 4.6 mm
Particle Size: 3 μ m

Conditions

Mobile phase: acetone:acetonitrile (70:30)
Flow: 2.5 mL/min
Detector: RI

Notes

Recommended Column Holder	07150029
Recommended Guard Cartridge	00402237
Recommended Guard Holder	07150001
Recommended Column and Guard Holder	07150030
	07150018
	07150002

PerkinElmer, Inc.
940 Winter Street
Waltham, MA 02451 USA
P: (800) 762-4000 or
(+1) 203-925-4602
www.perkinelmer.com


PerkinElmer[®]
For the Better

For a complete listing of our global offices, visit www.perkinelmer.com/ContactUs

Copyright ©2013 PerkinElmer, Inc. All rights reserved. PerkinElmer[®] is a registered trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.

007572A_49