


PROTEINS, PEPTIDES AND CATECHOLAMINES HPLC

Biochemical: Reversed Phase Test Mix on Brownlee Bio C18

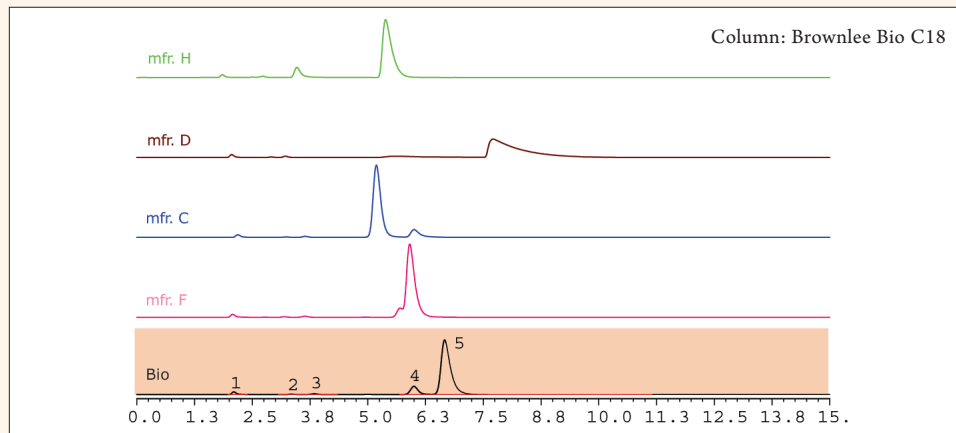

Peak List

1. uracil
2. toluene
3. ethylbenzene
4. quinizarin
5. amitriptyline

Column

Column name: Brownlee Bio C18
Part number: N9303648
Dimensions: 150 x 2.1 mm wide pore
Particle size: 5 μm
Pore size: 300 \AA

Conditions

Mobile phase: water:methanol (25:75 v/v)
Flow: 0.2 mL/min
Detector: UV @ 254 nm

Sample

Injection volume: 10 μL HPLC Performance
(NIST 870) Test Mix

Notes

Recommended Guard Cartridge N9303708
Recommended Guard Holder N9303713

PerkinElmer, Inc.
940 Winter Street
Waltham, MA 02451 USA
P: (800) 762-4000 or
(+1) 203-925-4602
www.perkinelmer.com


PerkinElmer[®]
For the Better

For a complete listing of our global offices, visit www.perkinelmer.com/ContactUs

Copyright ©2013 PerkinElmer, Inc. All rights reserved. PerkinElmer[®] is a registered trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.

007568A_17