


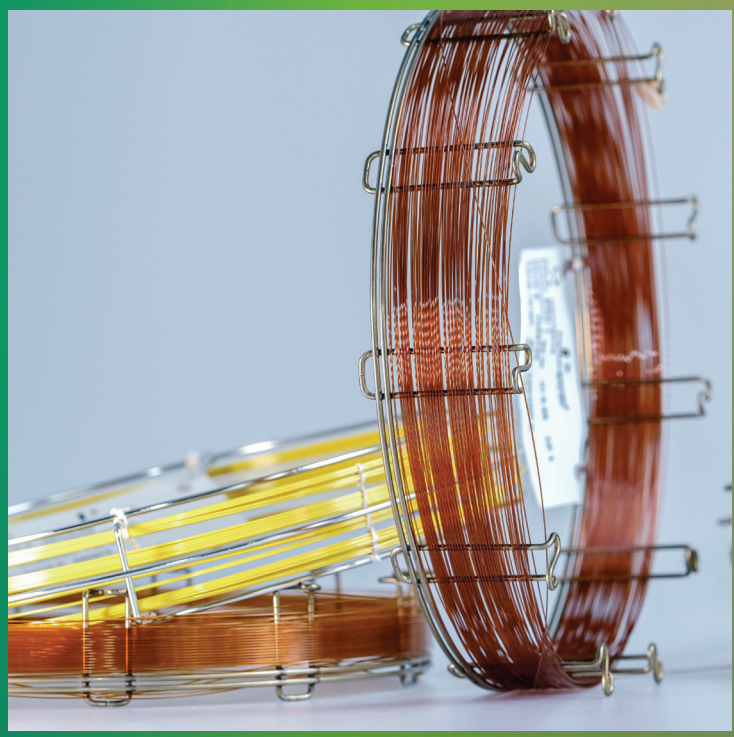

# GC COLUMNS GEARED TOWARDS THE FUTURE

INTERACTIVE PDF

[CLICK HERE TO VIEW THE MAIN INDEX](#)

Gas Chromatography Columns

  
**PerkinElmer**  
*For the Better*

# ORDERING ONLINE: IT'S THAT EASY



Sign Up and Get Started Today at [www.perkinelmer.com/register](http://www.perkinelmer.com/register)

## Quick and Efficient Three Steps - that's how easy it is to buy online

1. Visit [www.perkinelmer.com](http://www.perkinelmer.com)
2. Click on "Login/Register" at the top of the web page to register.  
All you need is your PerkinElmer account number.
3. Search for the product you need and add to your cart

### Benefits of ordering online:

- Real time product availability
- Immediate order processing
- Order confirmations
- Shipment notifications
- Quick re-order
- Save items to favorites list
- Tracking information and updates
- All orders available in "my account" order history

Step 1: Quick registration.

Step 2: Shop by product number, alphabetical listing, or by technology.

Step 3: Checkout via your shopping cart.



The power and convenience of real-time purchasing information are now available online.  
When will the product ship? Shipment date of each item is visible in the shopping cart.

Shop Online today!

[www.perkinelmer.com/category/consumables](http://www.perkinelmer.com/category/consumables)

# Chromatography

Today's advanced and increasingly diverse Advanced Materials laboratories are facing new challenges on a daily basis – starting from raw materials right up to the finished product. PerkinElmer's comprehensive portfolio of analytical solutions is designed to give you the higher accuracy, sensitivity, and ease of use your laboratory demands for examining with confidence, the purity, composition, and performance of your compounds. What's more, a range of complementary services is available to keep your laboratory up and running, meeting the stringent requirements of a variety of environments and working practices.

Better insights for better products.

## CHOOSE PERKINELMER.

### Ordering Information

Place your order online at: [www.perkinelmer.com/shop](http://www.perkinelmer.com/shop)

Online ordering is currently available in Australia, Austria, Belgium, Canada, China, Denmark, Finland, France, Germany, Ireland, Italy, Luxembourg, the Netherlands, Norway, Puerto Rico, Singapore, Spain, Sweden, Switzerland, the United Kingdom and the United States.

Please contact your local PerkinElmer Customer Care Representative.

Or visit [www.perkinelmer.com/contactus/#/offices](http://www.perkinelmer.com/contactus/#/offices).

Visit us online at [www.perkinelmer.com/contactus/#/offices](http://www.perkinelmer.com/contactus/#/offices)

to obtain the correct address for mailing your company Purchase Order.

### Return Goods Policy

Before returning any item to PerkinElmer, please obtain a Return Authorization Number from a customer service representative. Call your local PerkinElmer office and please have the following information ready:

- Our Order Number or your Purchase Order Number.
- The PerkinElmer part numbers of the items that you wish to return.
- Reason for return.

Our customer service representative will provide you with a Return Authorization Number and the return address. We cannot accept any returns without a Return Authorization Number. Policy may vary. Please contact your local PerkinElmer representative for details.

### Technical Support

For technical assistance, please contact your local PerkinElmer Customer Care Representative. See local office listings.

Check with your local PerkinElmer Customer Care center for current pricing.

PerkinElmer and the PerkinElmer logo are registered trademarks of PerkinElmer, Inc. All other product and company names mentioned herein may be the trademarks or registered trademarks of their respective owners.

PerkinElmer is not responsible for typographical or photographic errors that may appear in this catalog. We reserve the right to make adjustments or corrections due to changing market conditions, product discontinuations, or typographical or photographic errors in advertisements or product descriptions. Terms and conditions and item availability are subject to change without notice.

### Terms and Conditions

For details on our terms and conditions, including warranties, please visit: <https://www.perkinelmer.com/corporate/policies/terms-conditions-of-sale.html>

## Table of Contents

|                                                                                                 |             |
|-------------------------------------------------------------------------------------------------|-------------|
| <b>GC COLUMNS</b>                                                                               | <b>4-30</b> |
| ➤ GC Capillary Column Cutter                                                                    | 6           |
| ➤ Cross Reference Charts                                                                        | 7-8         |
| ➤ General Purpose GC Columns                                                                    | 9-19        |
| ➤ Elite-1                                                                                       | 9           |
| ➤ Elite-5                                                                                       | 10          |
| ➤ Elite-17, Elite-35                                                                            | 11-12       |
| ➤ Elite-200, Elite-225                                                                          | 13-14       |
| ➤ Elite-624, Elite-1301                                                                         | 15-16       |
| ➤ Elite-1701                                                                                    | 17          |
| ➤ Elite-WAX, Elite-MWAX, Elite-WAX ETR                                                          | 18-19       |
| ➤ GC Columns for GC/MS                                                                          | 20-23       |
| ➤ Elite-1ms                                                                                     | 20          |
| ➤ Elite-5ms                                                                                     | 21          |
| ➤ Elite-17ms                                                                                    | 22          |
| ➤ Elite-35ms                                                                                    | 22          |
| ➤ Elite-624ms                                                                                   | 23          |
| ➤ High Temperature GC Columns                                                                   | 24          |
| ➤ Speciality Columns                                                                            | 25-30       |
| ➤ Elite-BAC Advantage                                                                           | 25          |
| ➤ Elite-VMS, Elite-XLB                                                                          | 26          |
| ➤ Elite-Volatiles, Elite-CLPesticides                                                           | 27          |
| ➤ Elite-502.2, Elite-RX, Elite-Betacyclodextrin, Elite-SimDist, Elite-608                       | 28          |
| ➤ Elite-TPH, Elite-PONA, Elite-FFAP, Elite-5 Amine, Elite-2330/Elite-23, Elite-MTBE, Elite-2560 | 29          |
| ➤ Elite-Carbon                                                                                  | 29          |
| ➤ PLOT Columns                                                                                  | 30          |
| ➤ Miscellaneous Accessories                                                                     | 30          |
| ➤ ORDERING INFORMATION                                                                          | 31-36       |

# GC Columns


PerkinElmer offers a wide selection of superior quality products designed to work with your PerkinElmer instruments. Our precision designed products deliver the peace of mind that comes from knowing that you'll get the results you need.

## Quick Reference Index

|                                                                                                 | Page  |
|-------------------------------------------------------------------------------------------------|-------|
| ➤ Cross Reference Charts                                                                        | 7-8   |
| ➤ GC Capillary Column Cutter                                                                    | 6     |
| ➤ GC Columns for GC/MS                                                                          | 20-23 |
| ➤ Elite-1ms                                                                                     | 20    |
| ➤ Elite-5ms                                                                                     | 21    |
| ➤ Elite-17ms                                                                                    | 22    |
| ➤ Elite-35ms                                                                                    | 22    |
| ➤ Elite-624ms                                                                                   | 23    |
| ➤ General Purpose GC Columns                                                                    | 9-19  |
| ➤ Elite-1                                                                                       | 9     |
| ➤ Elite-5                                                                                       | 10    |
| ➤ Elite-17, Elite-35                                                                            | 11-12 |
| ➤ Elite-200, Elite-225                                                                          | 13-14 |
| ➤ Elite-624, Elite-1301                                                                         | 15-16 |
| ➤ Elite-1701                                                                                    | 17    |
| ➤ Elite-WAX, Elite-MWAX, Elite-WAX ETR                                                          | 18-19 |
| ➤ High Temperature GC Columns                                                                   | 24    |
| ➤ Miscellaneous Accessories                                                                     | 30    |
| ➤ Speciality Columns                                                                            | 25-30 |
| ➤ Elite-BAC Advantage                                                                           | 25    |
| ➤ Elite-VMS, Elite-XLB                                                                          | 26    |
| ➤ Elite-Volatiles, Elite-CLPesticides                                                           | 27    |
| ➤ Elite-502.2, Elite-RX, Elite-Betacyclodextrin, Elite-SimDist, Elite-608                       | 28    |
| ➤ Elite-TPH, Elite-PONA, Elite-FFAP, Elite-5 Amine, Elite-2330/Elite-23, Elite-MTBE, Elite-2560 | 29    |
| ➤ Elite-Carbon                                                                                  | 29    |
| ➤ PLOT Columns                                                                                  | 30    |


## Blood Alcohol Columns

The Elite-BAC Advantage columns are optimized for selectivities guaranteed to resolve ethanol, internal standards, and frequently encountered interferences.

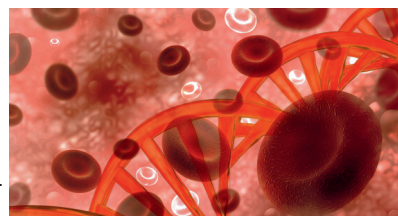
[VIEW PAGE](#)

## Elite MS Columns

The Elite range of MS columns are engineered for extremely low bleed for MS detectors, providing optimum sensitivity. They cover a wide range of polarities and applications.

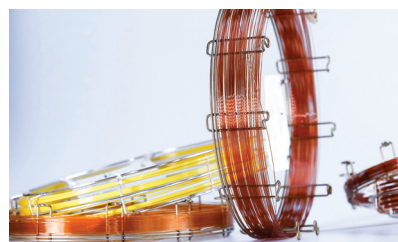
[VIEW PAGE](#)

## GC Column Cutter

PerkinElmer's capillary column cutting tool ensures you make a perfect cut of your GC column, first time, every time. The rotating diamond blade with a built in magnifier to verify a square cut, affords a precise clean cut of fused silica columns.


[VIEW PAGE](#)

## Elite Guard Columns

Using the Elite-Guard or the Elite-Siltek Guard Column lengthens the life of the capillary column and improves the analyte focusing.

[VIEW PAGE](#)



### Clarus® 590/690 GC

Sensitive, high-capacity, high-throughput GC systems delivering the power and functionality needed to meet your analytical goals. A robust autosampler delivers easy access to two injector ports, while the Clarus 690's patented high-performance oven delivers the fastest heat-up and cool-down of any oven in the business.



### TurboMatrix™ Headspace and Headspace Trap

TurboMatrix Headspace and high-sensitivity Headspace Trap samplers provide unparalleled precision and ease of use for numerous GC or GC/MS volatile-analysis applications. Different models offering a range of capacities are available to satisfy the requirements of virtually any laboratory.



### Torion® T-9 Portable GC/MS

The world's smallest portable GC/MS, Torion T-9 allows you to test samples where hazardous events occur: No sample processing, packaging, and prep time. You simply respond, collect, analyze, and identify in the field. And that means remedial action 70 times faster than with conventional labs.



### Clarus® SQ 8 GC/MS

With the flexibility to choose your level of sensitivity and dynamic range, the Clarus SQ 8 GC/MS eliminates background noise, maximizes analyte signals, and enables you to reconfigure between EI and CI, simply and quickly.



### TurboMatrix Thermal Desorbers

The five different TurboMatrix Thermal Desorbers allow you to match throughput and technology to your laboratory and applications needs. Use this clean technique to simplify and speed up a wide range of GC applications.


# Finest Quality High-Strength Fused Silica

## Why Choose Fused Silica?

Many factors influence the quality of a column. Fused silica is considered to be the purest form of glass, with fewer metal oxides (Lewis acid sites) and hydrogen bonding (surface silanol) groups. The stationary phase is cross linked (polymerized) and also bonded to the surface of the column to provide a high degree of stability, resulting in lower bleeding of the stationary phase at elevated temperatures. The superior inertness of the column means that acidic and basic compounds can be analyzed on the same column.

## Selecting the Right Stationary Phase

The inherent efficiency (large number of theoretical plates) of capillary columns allows you to choose from relatively few types of phases, compared to the many varieties of packed columns previously required. Perhaps more importantly, because capillary columns are more efficient, you will see superior resolution resulting in narrower, taller peaks that allow easier integration from your data system. Identification of small peaks are facilitated by a reduced baseline bleed and lower baseline noise. Non-polar Elite-1 columns from PerkinElmer will preferentially retain non-polar compounds, whereas the PerkinElmer Elite-200 column phase provides high selectivity for analytes containing lone pair electrons, such as nitro and carbonyl groups. Elite-WAX polyethylene glycol columns are highly selective toward polar compounds such as alcohols.

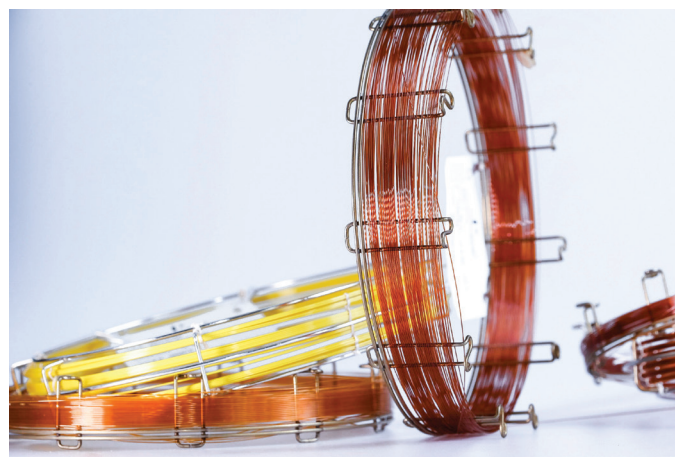

## What Length Do I Need?

Typically capillary columns are available in lengths from 15 to 105 meters. The longer the column the more resolving power, but this also increases the analysis time. Doubling a column length only increases resolution by approximately 40%. Under isothermal conditions, the analysis time will double if using temperature programmed analysis retention times are more dependent on the temperature than on the column length. We provide columns in the most popular lengths of 5, 10, 12, 15, 25, 30, 50, 60, 75, 100 and 105 meters depending upon the column ID.

## Elite Guard Columns

Using the Elite-Guard or the Elite-Siltek Guard Column lengthens the life of the capillary column and improves the analyte focusing. The 5 m length of deactivated uncoated fused silica is connected to the inlet end of the capillary column and traps nonvolatile residues, preventing them from collecting at the head of the analytical column. This length of fused silica contains no stationary phase adding only a minimal amount of time to the analysis.

| Description        | Tubing ID (mm) | Length (m) | Part No. |
|--------------------|----------------|------------|----------|
| Elite Siltek Guard | 0.25           | 5          | N9316607 |
| Elite Siltek Guard | 0.32           | 5          | N9316608 |
| Elite Siltek Guard | 0.53           | 5          | N9316609 |
| Elite Guard        | 0.25           | 5          | N9316603 |
| Elite Guard        | 0.32           | 5          | N9316604 |
| Elite Guard        | 0.53           | 5          | N9316606 |

## GC Capillary Column Cutter

Use PerkinElmer's capillary column cutting tool to ensure a perfect cut of your GC column, first time, every time. The rotating diamond blade with a built in magnifier to verify a square cut, affords a precise clean cut of fused silica columns. Suitable for use with 0.25 mm ID to 0.53 mm ID tubing, (0.78 mm OD maximum).



| Description                                                                                   | Part No. |
|-----------------------------------------------------------------------------------------------|----------|
| GC capillary cutting tool                                                                     | N6107245 |
| Maintenance kit<br>Contains replacement cutting wheel,<br>O-rings and tool to open the cutter | N6107246 |

Cross Reference Chart by Phase

| PerkinElmer   | Phase Composition                                      | USP                     | Agilent®                                                                     | Alltech®         | Machery-Nagel®                  | Ohio Valley® | Phenomenex®       | Quadrex®          | Restek®       | SGE®         | Supelco®        |
|---------------|--------------------------------------------------------|-------------------------|------------------------------------------------------------------------------|------------------|---------------------------------|--------------|-------------------|-------------------|---------------|--------------|-----------------|
| Elite-1       | Dimethyl polysiloxane                                  | G1, G2, G38             | HP-1, DB-1, CP-Sil 5 CB                                                      | 007-1AT-1, EC-1  | OPTIMA 1                        | OV-1         | ZB-1              | 007-1             | Rtx-1, Mtx-1  | BP1          | SPB-1           |
| Elite-1ht     | Dimethyl polysiloxane                                  | G1, G2, G38             | DB-1ht                                                                       | AT-1ht           |                                 |              | ZB-1HTInferno     |                   | Rxi-1HT       |              |                 |
| Elite-1ms     | Dimethyl polysiloxane (low bleed)                      | G1, G2, G38             | HP-1, HP-1ms, HP-1msUI, DB-1, DB-1MS, DB-1msUI, Ultra-1, VF-1ms, CP-Sil 5 CB | AT-1ms           | OPTIMA 1 MS, OPTIMA 1 MS Accent |              | ZB-1, ZB-1ms      | 007-1             | Rxi-1ms       | BP-1         | SPB-1, Equity-1 |
| Elite-5       | Diphenyl dimethyl polysiloxane                         | G27, G36                | HP-5, DB-5, CP-Sil 8 CB                                                      | EC-5, AT-5       | OPTIMA 5                        | OV-5         | ZB-5              | 007-5             | Rtx-5         | BP5          | SPB-5           |
| Elite-5ht     | Diphenyl dimethyl polysiloxane                         | G27, G36                | DB-5ht, VF-5ht                                                               |                  | OPTIMA 5HT                      |              | ZB-5HTInferno     |                   | Rxi-5HT       | HT5          |                 |
| Elite-5ms     | 1,4-bis(dimethylsiloxy)phenylene dimethyl polysiloxane | G27, G36                | DB-5ms, DB-5msUI, VF-5ms, CP Sil 8 CB MS                                     |                  | OPTIMA 5 MS Accent              | OV-5MS       | ZB-5msi           | 007-5MS           | Rxi-5Sil MS   | BPX5         | SLB-5ms         |
| Elite-17      | Phenyl methyl polysiloxane                             | G3                      | DB-17, CP Sil 24 CB                                                          | AT-50            |                                 | OV-17        |                   | 007-17            | Rtx-50        |              | SPB-50          |
| Elite-17ht    | Phenyl methyl polysiloxane                             | G3                      | DB-17ht                                                                      |                  |                                 |              |                   |                   |               |              |                 |
| Elite-17ms    | Phenyl methyl polysiloxane                             | G3                      | HP-50+, DB-17, DB-17ht, DB-608, CP Sil 24 CB                                 |                  | OPTIMA 17                       |              | ZB-50             |                   | Rxi-17        |              | SPB-17          |
| Elite-35      | Diphenyl dimethyl polysiloxane                         | G42                     | HP-35, DB-35, VF-35                                                          | AT-35, AT-35ms   |                                 | OV-35        | ZB-35             | 007-35            | Rtx-35        | BPX35,BPX608 | SPB-35,SPB-608  |
| Elite-35ms    | Diphenyl dimethyl polysiloxane                         | G42                     | DB-35ms, DB35msUI                                                            |                  | OPTIMA 35 MS                    |              | MR2               |                   | Rxi-35Sil MS  | BPX35        |                 |
| Elite-200     | Trifluoropropylmethyl polysiloxane                     | G6                      | DB-210, DB-200, VF-200ms                                                     | AT-210           | OPTIMA 210                      |              |                   |                   | Rtx-200       |              |                 |
| Elite-225     | Cyanopropylmethyl phenylmethyl polysiloxane            | G7, G19                 | DB-225ms, CP Sil 43 CB                                                       | AT-225           | OPTIMA 225                      | OV-225       |                   | 007-225           | Rtx-225       | BP225        | SPB-225         |
| Elite-624     | Cyanopropylphenyl dimethyl polysiloxane                | G43                     | DB-1301, DB-624, VF-624ms, CP-1301                                           | AT-624, AT-1301  | OPTIMA 1301, OPTIMA 624         | OV-624       | ZB-624            | 007-1301, 007-624 | Rtx-624       | BP624        | SPB-624         |
| Elite-624ms   | Unique phase                                           | G43                     | DB-624, VF-624ms, CP-Select 624 CB                                           |                  | OPTIMA 624 LB                   |              | ZB-624            |                   | Rxi-624Sil MS | BP624        |                 |
| Elite-1301    | Cyanopropylphenyl dimethyl polysiloxane                | G43                     | DB-1301, DB-624, VF-1301ms, VF-624ms, CP-1301                                | AT-624, AT-1301  | OPTIMA 1301, OPTIMA 624         | OV-1301      | ZB-624            | 007-1301, 007-624 | Rtx-624       | BP624        | SPB-624         |
| Elite-1701    | Cyanopropylphenyl dimethyl polysiloxane                | G46                     | DB-1701R, DB-1701, CP Sil 19 CB, VF-1701ms, VF-1701 Pesticides               | AT-1701          | OPTIMA 1701                     | OV-1701      | ZB-1701, ZB-1701P | 007-1701          | Rtx-1701      | BP10         | Equity-1701     |
| Elite-WAX     | Polyethylene glycol                                    | G14, G15, G16, G20, G39 | DB-Wax, CP Wax 52 CB, VF WAX                                                 | AT-WAXms, EC-WAX | OPTIMA WAX                      | Carbowax 20M | ZB-Wax            | 007-CW            | Rtx-Wax       | BP20         |                 |
| Elite-WAX ETR | Polyethylene glycol                                    | G14, G15, G16, G20, G39 | HP-INNOWax, CP Wax 52 CB, VF-WAX MS                                          | AT-WAX           | OPTIMA WAX plus                 |              | ZB-WaxPLUS        |                   | Stabilwax     |              | Supelcowax-10   |

Cross Reference Chart by Application

| PerkinElmer           | Applications                                                               | Agilent®                                  | Alltech®               | Machery-Nagel®                                | Ohio Valley® | Phenomenex® | Quadrex® | Restek®            | SGE®    | Supelco®                       |
|-----------------------|----------------------------------------------------------------------------|-------------------------------------------|------------------------|-----------------------------------------------|--------------|-------------|----------|--------------------|---------|--------------------------------|
| Elite-23              | cis/trans FAMES and Dioxins                                                | VF-23ms                                   | AT-Silar90             |                                               |              |             | 007-23   |                    | BPX70   | SP-2330, SP-2331, SP-2380      |
| Elite-502             | Volatile analytes by EPA Method 502.2                                      | DB-502.2                                  |                        |                                               |              |             |          | Rtx-502.2          |         | VOCOL                          |
| Elite-608             | Semivolatile pesticides by EPA Method 608                                  | DB-608, HP-608                            |                        |                                               |              |             | 007-608  |                    |         | SPB-608                        |
| Elite-2560            | cis/transFAMES                                                             | HP-88, CP Sil 88                          |                        |                                               |              |             |          | Rt-2560            |         | SPB-2560                       |
| Elite-Alumina PLOT    | Light hydrocarbons                                                         | Alumina PLOT                              |                        |                                               |              |             |          |                    |         |                                |
| Elite-BAC 1 Advantage | Blood alcohol testing                                                      | DB-ALC1                                   |                        |                                               |              | ZB-BAC1     |          | Rtx-BAC Plus 1     |         |                                |
| Elite-BAC 2 Advantage | Blood alcohol testing                                                      | DB-ALC2                                   |                        |                                               |              | ZB-BAC2     |          | Rtx-BAC Plus 2     |         |                                |
| Elite-CLPesticides    | Organochlorine pesticides by EPA Methods 504, 608, 8081, 8082, and CLP     | DB-CLP1                                   |                        |                                               |              |             |          | Rtx-CLPesticides   |         |                                |
| Elite-CLPesticides2   | Organochlorine pesticides by EPA Methods 504, 608, 8081, 8082, and CLP     | DB-CLP2                                   |                        |                                               |              |             |          | Rtx-CLPesticides 2 | Rtx-200 |                                |
| Elite-Cyclosil B      | Chiral separations                                                         |                                           |                        |                                               |              |             |          |                    |         |                                |
| Elite-FFAP            | Free fatty acids                                                           | HP-FFAP, DB-FFAP, CP WAX58 CB, CP-FFAP CB | AT-AquaWax DA, AT-1000 | PERMABOND FFAP, OPTIMA FFAP, OPTIMA FFAP Plus | OV-351       | ZB-FFAP     |          | Stabilwax-DA       | BP-21   | Nukol                          |
| Elite-Molesieve PLOT  | Permanent gases                                                            |                                           |                        |                                               |              |             |          |                    |         |                                |
| Elite-PONA            | Detailed analysis of petroleum naphtha                                     | HP-PONA, DB-Petro, CP Sil PONA CB         |                        |                                               |              |             |          | Rtx-DHA            | BP1PONA | Petrocol DH                    |
| Elite-Carbon          | Permanent gases and light hydrocarbons                                     |                                           |                        |                                               |              |             |          |                    |         |                                |
| Elite-SimDist         | Simulated Distillation and Hydrocarbons – ASTM 2887                        | DB-2887, CP SimDist                       | AT-2887                |                                               |              |             |          | Rtx-2887           |         | Petrocol 2887, Petrocol EX2887 |
| Elite-THP             | Total petroleum hydrocarbons                                               |                                           |                        |                                               |              |             |          |                    |         |                                |
| Elite-VMS             | Volatiles Organic Pollutants by GC-MS for EPA Methods 8260,624,524         | Unique Phase                              |                        |                                               |              |             |          |                    |         |                                |
| Elite-VRX             | Volatile analytes by EPA Methods 502.2, 601, 602, 8010, 8020               | DB-VRX                                    |                        |                                               |              |             |          |                    |         |                                |
| Elite-XLB             | Polychlorinated biphenyl analytes by EPA Methods 8082, 6008, PCB congeners | DB-XLB, VF-XMS                            |                        |                                               |              | MR1, ZB-XLB |          | Rxi-XLB            |         |                                |

# Elite-1

The Elite-1 100% dimethyl polysiloxane columns is a highly versatile, non-polar, cross-linked general purpose phase that is rugged, exhibiting long column lifetime, low bleed, and high maximum operating temperatures.


## Features

- Temperature Range: -60 °C to 350 °C
- Equivalent to USP G1, G2, and G38 phases

## Applications

- Ideal for analysis of non-polar petrochemical samples
- Also excellent for solvents, chemicals, flavors & fragrances, air toxins and pesticides

Elite-1 Structure

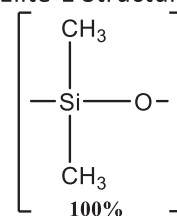

| ID (mm) | df (μm) | Temp Limits (°C) | 5 m Part No. | 15 m Part No. | 30 m Part No.         | 60 m Part No. | 105 m Part No. |
|---------|---------|------------------|--------------|---------------|-----------------------|---------------|----------------|
| 0.25    | 0.10    | -60 to 330/350   |              | N9316006      | N9316009              | N9316012      |                |
|         | 0.25    | -60 to 330/350   |              | N9316007      | N9316010              | N9316013      |                |
|         | 0.50    | -60 to 330/350   |              | N9316686      | N9316685              |               |                |
|         | 1.00    | -60 to 320/340   |              | N9316008      | N9316011              | N9316014      |                |
| 0.32    | 0.10    | -60 to 330/350   |              | N9316016      | N9316022              | N9316027      |                |
|         | 0.25    | -60 to 330/350   | N9316596     | N9316017      | N9316023              | N9316028      |                |
|         | 0.50    | -60 to 330/350   |              |               | N9316021 <sup>1</sup> | N9316691      |                |
|         | 1.00    | -60 to 320/340   |              | N9316018      | N9316024              | N9316029      |                |
|         | 1.50    | -60 to 310/330   |              |               | N9316050              | N9316580      |                |
|         | 3.00    | -60 to 280/300   |              | N9316019      | N9316025              | N9316030      |                |
|         | 5.00    | -60 to 260/280   |              | N9316020      | N9316026              | N9316031      |                |
| 0.45    | 0.13    | -60 to 330/350   |              | N9316032      |                       |               |                |
|         | 0.42    | -60 to 310/330   |              | N9316037      | N9316041              |               |                |
|         | 1.27    | -60 to 310/330   |              | N9316034      | N9316038              | N9316042      |                |
|         | 2.55    | -60 to 270/290   |              | N9316035      | N9316039              |               | N9316043       |
|         | 4.25    | -60 to 260/280   | N9316032     | N9316036      | N9316040              |               |                |
| 0.53    | 0.15    | -60 to 320/340   |              | N9316045      |                       |               |                |
|         | 0.50    | -60 to 310/330   |              | N9316049      | N9316053              |               |                |
|         | 1.50    | -60 to 310/330   |              | N9316046      | N9316050              | N9316054      |                |
|         | 3.00    | -60 to 270/290   |              | N9316047      | N9316051              | N9315499      | N9316692       |
|         | 5.00    | -60 to 270/290   | N9316044     | N9316048      | N9316052              |               |                |

| ID (mm) | df (μm) | Temp Limits (°C) | 10 m Part No. | 12 m Part No. | 20 m Part No. | 25 m Part No. | 50 m Part No.         |
|---------|---------|------------------|---------------|---------------|---------------|---------------|-----------------------|
| 0.05    | 0.05    | -60 to 330/350   | N9316056      |               |               |               |                       |
|         | 0.20    | -60 to 330/350   | N9316057      |               |               |               |                       |
| 0.10    | 0.10    | -60 to 330/350   | N9316058      |               |               |               |                       |
|         | 0.40    | -60 to 320/340   |               |               | N9316061      |               |                       |
| 0.18    | 0.18    | -60 to 330/350   | N9316001      |               | N9316003      |               |                       |
|         | 0.40    | -60 to 320/340   | N9316002      |               | N9316004      |               | N9316005 <sup>2</sup> |
| 0.20    | 0.33    | -60 to 330/350   |               | N9316062      |               | N9316063      | N9316064              |

<sup>1</sup> N9316021: Elite-1, 25M x 0.32 mm x 0.52 μm. <sup>2</sup> The length of N9316005 is 40 m.

## Elite-5

The Elite-5 is a 5% diphenyl/95% dimethyl polysiloxane stationary phase. It is regarded as a general purpose, low polarity phase that is the most popular GC stationary phase used for a wide variety of applications. A crosslinked phase in which all residual catalysts and low molecular weight fragments have been removed providing a tight mono-modal distribution and extremely low bleed.


### Features

- Temperature Range: -60 °C to 350 °C
- Equivalent to USP G27 and G36 phases

### Applications

- Drugs, pesticides and solvent impurities
- Hydrocarbons and PCBs
- Essential oils and semivolatiles

Elite-5 Structure

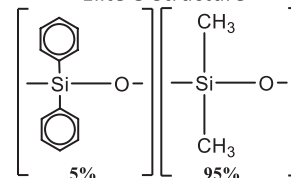


| ID (mm) | df (μm)         | Temp Limits (°C) | 10 m Part No.               | 15 m Part No.   | 20 m Part No.               | 30 m Part No.   | 50 m Part No.               | 60 m Part No.   |
|---------|-----------------|------------------|-----------------------------|-----------------|-----------------------------|-----------------|-----------------------------|-----------------|
| 0.05    | 0.05            | -60 to 325/350   | <b>N9316104</b>             |                 |                             |                 |                             |                 |
|         | 0.20            | -60 to 325/350   | <b>N9316105</b>             |                 |                             |                 |                             |                 |
| 0.10    | 0.10            | -60 to 330/350   |                             |                 | <b>N9316108</b>             |                 |                             |                 |
|         | 0.40            | -60 to 320/340   |                             |                 | <b>N9316109</b>             |                 |                             |                 |
| 0.18    | 0.18            | -60 to 330/350   | <b>N9316066</b>             |                 | <b>N9316068</b>             |                 |                             |                 |
|         | 0.4             | -60 to 320/340   | <b>N9316067</b>             |                 | <b>N9316069</b>             |                 | <b>N9316071<sup>2</sup></b> |                 |
| 0.20    | 0.33            | -60 to 330/350   | <b>N9316110<sup>1</sup></b> |                 | <b>N9316111<sup>1</sup></b> |                 | <b>N9316112</b>             |                 |
| 0.25    | 0.10            | -60 to 330/350   |                             | <b>N9316072</b> |                             | <b>N9316075</b> |                             | <b>N9316078</b> |
|         | 0.25            | -60 to 330/350   |                             | <b>N9316073</b> |                             | <b>N9316076</b> |                             | <b>N9316079</b> |
|         | 0.50            | -60 to 330/350   |                             |                 |                             | <b>N9316107</b> | <b>N6107815</b>             |                 |
|         | 1.00            | -60 to 320/340   |                             | <b>N9316074</b> |                             | <b>N9316077</b> |                             | <b>N9316080</b> |
| 0.32    | <b>NEW</b> 0.10 | -60 to 330/350   |                             | <b>N9316081</b> |                             | <b>N9316085</b> |                             | <b>N9316089</b> |
|         | 0.25            | -60 to 330/350   |                             | <b>N9316082</b> |                             | <b>N9316086</b> |                             | <b>N9316090</b> |
|         | 1.00            | -60 to 320/340   |                             | <b>N9316083</b> |                             | <b>N9316087</b> |                             | <b>N9316091</b> |
| 0.45    | 0.13            | -60 to 340/350   |                             |                 |                             |                 |                             | <b>N9316097</b> |
|         | 0.42            | -60 to 310/330   |                             | <b>N9316093</b> |                             | <b>N9316096</b> |                             |                 |
|         | 1.27            | -60 to 310/330   |                             | <b>N9316092</b> |                             | <b>N9316094</b> |                             |                 |
|         | 4.25            | -60 to 260/280   |                             |                 |                             | <b>N9316095</b> |                             |                 |
| 0.53    | 0.50            | -60 to 310/330   |                             | <b>N9316099</b> |                             | <b>N9316102</b> |                             |                 |
|         | 1.50            | -60 to 310/330   |                             | <b>N9316098</b> |                             | <b>N9316100</b> |                             | <b>N9316103</b> |
|         | 5.00            | -60 to 270/290   |                             |                 |                             | <b>N9316101</b> |                             |                 |

<sup>1</sup> The lengths of **N9316110** and **N9316111** are 12 m and 25 m respectively. <sup>2</sup> **N9316071** use 40 m length column.

### Industrial

#### Determination of sulfur compounds in air by online TD-GC/FPD.

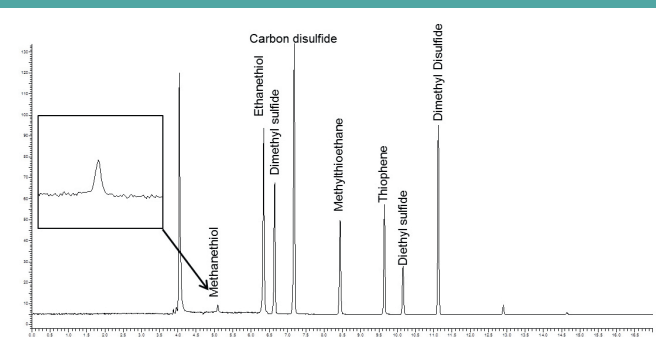

Recommended Column: Elite-5, 60 m x 0.32 mm x 0.25 mm, Part No. **N9316090**

# Elite-17

The Elite-17 columns are general purpose, mid-polarity, (50%-phenyl)-methylpolysiloxane phases and incorporates a crosslinking technology for very low bleed and long column lifetimes.


## Features

- Temperature Range: 40 °C to 330 °C
- Equivalent to USP G3 phase

## Applications

- Herbicides and pesticides
- Phthalate esters, sterols and rosin acids

Elite-17 Structure

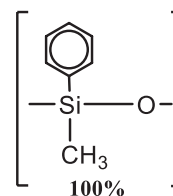

| ID (mm) | df (μm) | Temp Limits (°C) | 5 m Part No. | 10 m Part No. | 15 m Part No. | 20 m Part No. | 30 m Part No. | 60 m Part No. |
|---------|---------|------------------|--------------|---------------|---------------|---------------|---------------|---------------|
| 0.05    | 0.05    | 40 to 280/300    |              | N9316138      |               |               |               |               |
|         | 0.10    | 40 to 280/300    |              | N9316139      |               |               |               |               |
| 0.10    | 0.02    | 40 to 280/300    |              | N9316141      |               |               |               |               |
|         | 0.10    | 40 to 280/300    |              | N9316140      |               | N9316142      |               |               |
|         | 0.20    | 40 to 280/300    |              |               |               | N9316143      |               |               |
| 0.18    | 0.18    | 40 to 310/330    |              | N9316113      |               | N9316115      |               |               |
|         | 0.3     | 40 to 300/320    |              | N9316114      |               | N9316116      |               |               |
| 0.25    | 0.15    | 40 to 300/320    |              |               | N9316117      |               | N9316120      |               |
|         | 0.25    | 40 to 300/320    |              |               | N9316118      |               | N9316121      | N9316123      |
|         | 0.50    | 40 to 290/310    |              |               | N9316119      |               | N9316122      |               |
| 0.32    | 0.15    | 40 to 300/320    |              |               | N9316124      |               | N9316127      |               |
|         | 0.25    | 40 to 300/320    |              |               | N9316125      |               | N9316128      |               |
|         | 0.50    | 40 to 290/310    |              |               | N9316126      |               | N9316129      |               |
| 0.45    | 0.85    | 40 to 270/290    |              |               | N9316131      |               | N9316132      | N9316133      |
| 0.53    | 1.00    | 40 to 260/280    |              |               | N9316135      |               | N9316136      | N9316137      |
|         | 2.00    | 40 to 250/270    | N9316134     |               |               |               |               |               |



## Elite-35

The Elite-35 columns are general purpose, mid-polarity columns that are coated with a crosslinked, (35%-diphenyl)-dimethylpolysiloxane commonly used for organochlorine pesticides, PDB congeners. It is a popular confirmation column for pesticides and herbicides, in conjunction with an Elite-5 or Elite-1701. The higher phenyl content results in useful elution order and retention time changes.


### Features

- Temperature Range: 40 °C to 300/320 °C
- Equivalent to USP G42 phase

### Applications

- Pesticides and herbicides
- Pharmaceuticals, sterols, rosin acids and phthalate esters


Elite-35 Structure

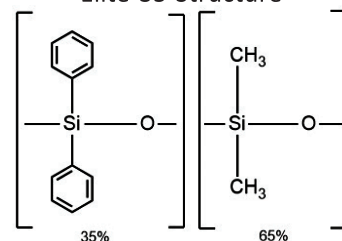

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.   | 30 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|
| 0.25    | 0.15    | 40 to 300/320    |                 | <b>N9316144</b> |
|         | 0.25    | 40 to 300/320    |                 | <b>N9316145</b> |
| 0.32    | 0.25    | 40 to 300/320    |                 | <b>N9316146</b> |
|         | 0.50    | 40 to 290/310    |                 | <b>N9316147</b> |
| 0.45    | 0.42    | 40 to 290/310    |                 | <b>N9316150</b> |
|         | 0.85    | 40 to 280/300    | <b>N9316148</b> | <b>N9316149</b> |
| 0.53    | 0.50    | 40 to 260/280    | <b>N9303929</b> | <b>N9316153</b> |
|         | 1.00    | 40 to 260/280    | <b>N9316151</b> | <b>N9316152</b> |

### Food and Flavor

#### Analysis of ethyl carbamate using GC/MS.

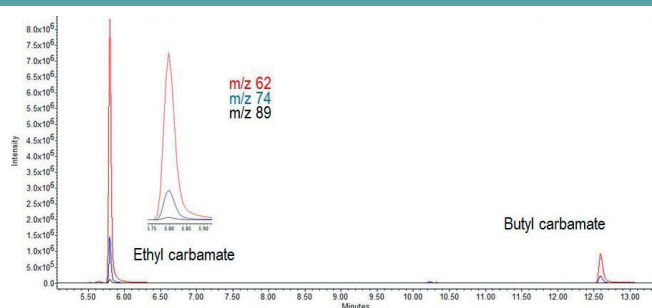

Recommended Column: Elite-35, 30 m x 0.25 mm x 0.25 μm  
(Part No. **N9316145**)

Liner: 2 mm ID glass (no glass wool)

## Elite-200

Elite-200 columns are comprised of a (trifluoropropyl)-methylpolysiloxane stationary phase that has a unique selectivity which changes elution orders and resolves compounds that phenyl, cyano, or Carbowax® phase cannot. These columns have accomplished many difficult separations not possible on any other bonded stationary phase. It offers exceptional thermal stability, low bleed, and superior inertness – even for active compounds such as phenols, and with sensitive detectors such as ECDs, NPDs, and MSDs. It is often used for the confirmation of phenols, nitrosamines, organochlorine pesticides, chlorinated hydrocarbons, and chlorophenoxy herbicides when paired with an Elite-5 column.


### Features

- Temperature Range: 40 °C to 320/340 °C
- Equivalent to USP G6 phase

### Applications

- Solvents, fluorocarbons, ketones and phenols
- Alcohols and drugs of abuse
- Chlorinated herbicides and pesticides

Elite-200 Structure

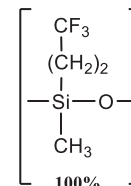

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No. | 30 m Part No. | 60 m Part No. |
|---------|---------|------------------|---------------|---------------|---------------|
| 0.25    | 0.10    | -20 to 320/340   | N9316616      | N9316617      |               |
|         | 0.25    | -20 to 320/340   | N9316618      | N9316619      |               |
|         | 0.50    | -20 to 310/330   | N9316620      | N9316621      |               |
|         | 1.00    | -20 to 290/310   | N9316622      | N9316623      | N9316624      |
| 0.32    | 0.25    | -20 to 320/340   | N9316625      | N9316626      |               |
|         | 0.50    | -20 to 310/330   | N9316627      | N9316628      |               |
|         | 1.00    | -20 to 290/310   | N9316629      | N9316630      | N9316631      |
|         | 1.50    | -20 to 280/300   | N9316632      | N9316633      | N9316634      |
| 0.53    | 0.25    | -20 to 310/330   | N9316635      | N9316636      | N9316637      |
|         | 0.50    | -20 to 300/320   | N9316638      | N9316639      | N9316640      |
|         | 1.00    | -20 to 290/310   | N9316641      | N9316642      | N9316643      |
|         | 1.50    | -20 to 280/300   | N9316644      | N9316645      | N9316646      |
|         | 3.00    | -20 to 260/280   | N9316647      | N9316648      | N9316649      |

## Elite-225


The Elite-225 is a general purpose column for the analysis of FAMES, carbohydrates, sterols and flavor compounds. The cyanopropyl-containing Elite-225 phase is slightly less polar than bonded polyethylene glycol (PEG) phases, but it can be used for many of the same applications. Improvements to the Elite-225 polymer have increased thermal stability, reduced bleed, and improved inertness. In other similar columns, the Carbowax® deactivation layer is not fully compatible with the cyanopropyl siloxane polymer, which can cause tailing of active compounds, and lower efficiency.

### Features

- Temperature Range: 40 °C to 220/240 °C
- Equivalent to USP G7, G19 phases

### Applications

- FAMES and carbohydrates
- Sterols and flavor compounds

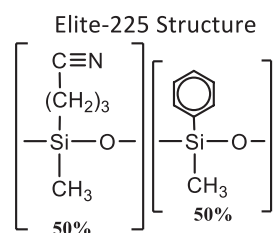

| ID<br>(mm) | df<br>(μm) | Temp Limits<br>(°C) | 10 m<br>Part No. | 15 m<br>Part No. | 20 m<br>Part No. | 30 m<br>Part No. | 60 m<br>Part No. |
|------------|------------|---------------------|------------------|------------------|------------------|------------------|------------------|
| 0.05       | 0.05       | 45 to 220/240       | N9316186         |                  |                  |                  |                  |
| 0.10       | 0.10       | 45 to 220/240       |                  |                  | N9316187         |                  |                  |
| 0.18       | 0.2        | 45 to 220/240       | N9316172         |                  | N9316173         |                  |                  |
| 0.25       | 0.15       | 45 to 220/240       |                  | N9316174         |                  | N9316176         | N9305631         |
|            | 0.25       | 45 to 220/240       |                  | N9316175         |                  | N9316177         |                  |
| 0.32       | 0.15       | 45 to 220/240       |                  | N9316178         |                  | N9316180         |                  |
|            | 0.25       | 45 to 220/240       |                  | N9316179         |                  | N9316181         |                  |
| 0.45       | 0.85       | 40 to 200/220       |                  | N9316182         |                  | N9316183         |                  |
| 0.53       | 1.00       | 40 to 200/220       |                  | N9316184         |                  | N9316185         |                  |

## Elite-624

The Elite-624 column is a specially engineered, low to mid-polarity (6%-cyanopropylphenyl)-dimethylpolysiloxane phase. The unique polarity of this phase makes it ideal for analyzing volatile organic pollutants and it is recommended in U.S. EPA methods. The Elite-624 phase produces greater than 90% resolution of the first six gases in EPA Methods 8260 and 524.2. This stationary phase is especially well-suited for EPA Method 524.2 since it resolves 2-nitropropane from 1,1-dichloropropanone, which share quantification ion m/z 43 and must be separated chromatographically.


### Features

- Temperature Range: -20 °C to 240 °C
- Equivalent to USP G43 phase

### Applications

- Volatile organic pollutants
- EPA methods 524.2 and 8260


### Elite-624 Structure

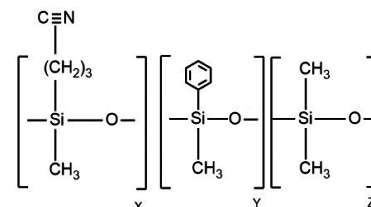

| ID (mm) | df (μm) | Temp Limits (°C) | 20 m Part No.   | 25 m Part No.   | 30 m Part No.   | 60 m Part No.   | 75 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|-----------------|-----------------|-----------------|
| 0.18    | 1.00    | -20 to 240       | <b>N9316200</b> |                 |                 |                 |                 |
| 0.20    | 1.12    | -20 to 240       |                 | <b>N9316209</b> |                 |                 |                 |
| 0.25    | 1.40    | -20 to 240       |                 |                 | <b>N9316201</b> | <b>N9316202</b> |                 |
| 0.32    | 1.80    | -20 to 240       |                 |                 | <b>N9316203</b> | <b>N9316204</b> |                 |
| 0.45    | 2.55    | -20 to 240       |                 |                 | <b>N9316205</b> |                 | <b>N9316206</b> |
| 0.53    | 3.00    | -20 to 240       |                 |                 | <b>N9316207</b> | <b>N9305699</b> | <b>N9316208</b> |

### Pharma

#### Residual solvents in pharmaceuticals by USP 467.

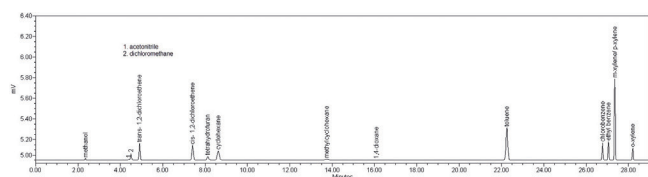


Elite-624, 30 m x 0.32 mm x 1.8 μm, Part No. **N9316203**


### Environmental

#### Analysis of Volatile Organic Compounds (VOCs) in air using US EPA Method TO-17.


Vinyl Chloride 0.05 μg/m<sup>3</sup>  
Single to noise: 255:1

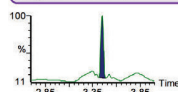

Chloroform 0.05 μg/m<sup>3</sup>  
Single to noise: 250:1

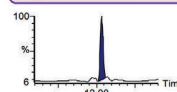


Benzene 0.05 μg/m<sup>3</sup>  
Single to noise: 450:1

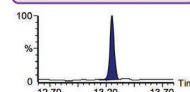


Trichloroethylene 0.05 μg/m<sup>3</sup>  
Single to noise: 670:1



Tetrachloroethylene 0.05 μg/m<sup>3</sup>  
Single to noise: 600:1

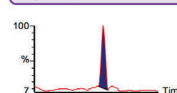

Naphthalene 0.05 μg/m<sup>3</sup>  
Single to noise: 710:1

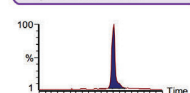

Recommended Column: Elite-624, 60 m x 0.25 mm x 1.4 μm, Part No. **N9316006**

# Elite-1301

The Elite-1301 column is a general purpose low to mid-polarity phase commonly used for the analysis of residual solvents, alcohols, oxygenates and volatile organic compounds. Our polymer is fully characterized to ensure long-term reproducibility, column-to-column consistency, and low bleed – even with sensitive detectors such as ECD and MS.


## Features

- Temperature Range: -20 °C to 280 °C
- Equivalent to USP G43 phase

## Applications

- Residual solvents, alcohols
- Oxygenates, VOCs

Elite-1301 Structure

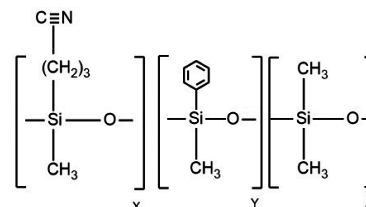

| ID (mm) | df (μm) | Temp Limits (°C) | 10 m Part No. | 15 m Part No. | 20 m Part No. | 30 m Part No. | 60 m Part No. |
|---------|---------|------------------|---------------|---------------|---------------|---------------|---------------|
| 0.18    | 0.40    | -20 to 280       | N9316210      |               | N9316211      |               |               |
| 0.25    | 0.25    | -20 to 280       |               | N9316212      |               | N9316214      | N9316216      |
|         | 1.00    | -20 to 260/280   |               |               |               | N9316215      | N9316217      |
| 0.32    | 0.25    | -20 to 280       |               | N9316218      |               | N9316220      | N9316222      |
|         | 1.00    | -20 to 260/280   |               | N9316219      |               | N9316221      | N9316223      |
| 0.45    | 0.85    | -20 to 260/280   |               | N9316224      |               | N9316225      |               |
| 0.53    | 1.00    | -20 to 260/280   |               | N9316226      |               | N9316227      |               |



## Elite-1701

The Elite-1701 has a stationary phase of (14%-cyanopropylphenyl)-methylpolysiloxane. It is regarded as a good general purpose column for the analysis of alcohols, oxygenates, and pesticides. The mix of cyano and phenyl functional groups increases the polarity and offers a different elution order relative to less polar Elite-1 or Elite-5 columns. An Elite-1701 column is ideal for confirmation analysis in combination with an Elite-35 or Elite-5 column. The polymer is fully characterized to ensure long-term reproducibility, column-to-column consistency, and low bleed – even with sensitive detectors such as ECD and MS.


### Features

- Temperature Range: -20 °C to 280 °C
- Equivalent to USP G46 phase

### Applications

- Alcohols, oxygenates
- PCB congeners, pesticides

Elite-1701 Structure

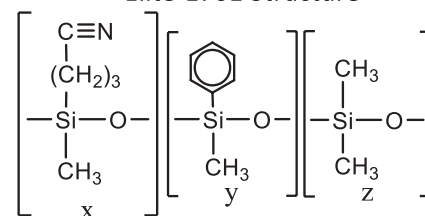

| ID (mm) | df (μm) | Temp Limits (°C) | 10 m Part No. | 15 m Part No. | 20 m Part No. | 30 m Part No. | 50 m Part No. | 60 m Part No. |
|---------|---------|------------------|---------------|---------------|---------------|---------------|---------------|---------------|
| 0.05    | 0.05    | -20 to 280       | N9316257      |               |               |               |               |               |
|         | 0.20    | -20 to 280       | N9316258      |               |               |               |               |               |
| 0.10    | 0.10    | -20 to 280       | N9316259      |               |               |               |               |               |
| 0.18    | 0.4     | -20 to 270/280   | N9316228      |               | N9316229      |               |               |               |
| 0.25    | 0.15    | -20 to 280       | N9316230      |               |               |               |               | N9316236      |
|         | 0.25    | -20 to 280       | N9316231      |               | N9316234      |               | N9316237      |               |
|         | 1.00    | -20 to 260/280   | N9316232      |               | N9316235      |               | N9316238      |               |
| 0.32    | 0.15    | -20 to 280       | N9316239      |               | N9316242      |               | N9316246      |               |
|         | 0.25    | -20 to 280       | N9316240      |               | N9316243      |               | N9316247      |               |
|         | 1.00    | -20 to 260/280   | N9326141      |               | N9316244      |               | N9316245      | N9316248      |
| 0.45    | 0.42    | -20 to 260/270   | N9316250      |               | N9316252      |               |               |               |
|         | 0.85    | -20 to 250/270   | N9316249      |               | N9316251      |               |               |               |
| 0.53    | 0.50    | -20 to 260/270   | N9316254      |               | N9316256      |               |               |               |
|         | 1.00    | -20 to 250/270   | N9316253      |               | N9316255      |               |               |               |

## Elite-WAX

The Elite-WAX column, a Polar Polyethylene Glycol (PEG) stationary phase column, is a general purpose polar PEG phase commonly used for the analysis of polar compounds like alkenols, glycols and aldehydes. The operating temperature range up to 250 °C facilitates the analysis of compounds that have a wide volatility range. Selectivity of the Elite-WAX is comparable to other Carbowax® columns for compounds of intermediate to high polarity.


### Features

- Temperature Range: 20 °C to 250 °C
- Equivalent to USP G14, G15, G16, G20 and G39 phases

### Applications

- FAMES, Glycols
- Alkenols, aldehydes, solvents

### Elite-WAX Structure

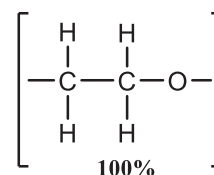


| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No. | 30 m Part No. | 60 m Part No. |
|---------|---------|------------------|---------------|---------------|---------------|
| 0.25    | 0.15    | 20 to 250        | N9316399      |               | N9316405      |
|         | 0.25    | 20 to 250        | N9316400      | N9316403      | N9316406      |
|         | 0.50    | 20 to 250        | N9316401      | N9316404      | N9316407      |
| 0.32    | 0.15    | 20 to 250        | N9316408      | N9316411      |               |
|         | 0.25    | 20 to 250        | N9316409      | N9316412      | N9316416      |
|         | 0.50    | 20 to 250        | N9316410      | N9316413      | N9316417      |
| 0.45    | 0.42    | 20 to 250        | N9316420      | N9316422      |               |
|         | 0.85    | 20 to 240/250    | N9316419      | N9316421      | N9316423      |
|         | 1.70    | 50 to 230        | N9316418      |               |               |
| 0.53    | 0.50    | 20 to 250        | N9316426      | N9316428      |               |
|         | 1.00    | 20 to 240/250    | N9316425      | N9316427      | N9316429      |

### Elite-MWAX: Metal Column

| ID (mm) | df (μm) | Temp Limits (°C) | 30 m Part No. |
|---------|---------|------------------|---------------|
| 0.53    | 1.00    | 20 to 240/250    | N9316478      |

### Food and Flavor

The determination of low levels of benzene, toluene, ethylbenzene, xylenes and styrene in olive oil.

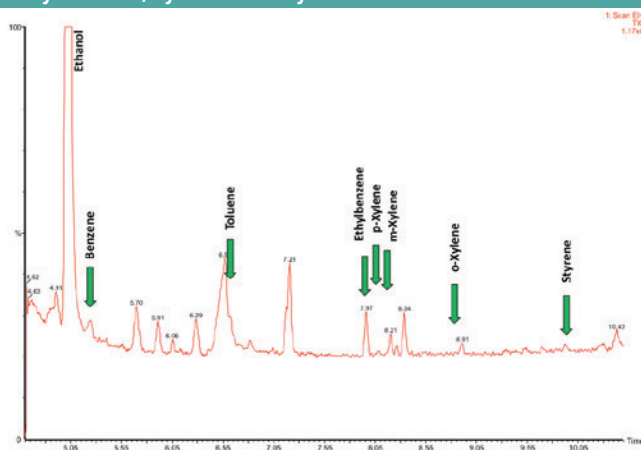

Recommended Column: Elite-WAX, 30 m x 0.25 mm x 1.0 μm (Part No. N9316485)

## Elite-WAX ETR

The Elite-WAX ETR (Extended Temperature Range) columns are manufactured with a special bonding process that binds the Carbowax® polymer to the polar deactivated silica. This results in a low bleed WAX column that exhibits extended lifetimes even when repeatedly heated to 250 °C. The bonding mechanism makes this column rugged enough to stand up to repeated water injections and allows solvent washing to rejuvenate the column. The Elite-WAX ETR has a wide applicability including FAMES, flavor compounds, acrolein/ acrylonitrile (EPA 603), oxygenated compounds, and impurities in water matrices.

### Features

- Temperature Range: 40 °C to 250 °C
- Equivalent to USP G14, G15, G16, G20 and G39 phases


### Applications

- FAMES, flavor compounds, essential oils
- Solvents, aromatics, alcoholic beverages
- EPA method 603

| ID (mm) | df (μm) | Temp Limits (°C) | 5 m Part No. | 15 m Part No. | 30 m Part No. | 50 m Part No. | 60 m Part No. |
|---------|---------|------------------|--------------|---------------|---------------|---------------|---------------|
| 0.25    | 0.25    | 40 to 250        |              | N9316547      | N9316549      |               | N9316551      |
|         | 0.50    | 40 to 250        |              | N9316548      | N9316550      |               |               |
| 0.32    | 0.25    | 40 to 250        |              | N9316552      | N9316555      |               | N9316559      |
|         | 0.50    | 40 to 250        |              | N9316553      | N9316556      |               | N9316560      |
|         | 1.00    | 40 to 240/250    |              | N9316554      | N9316557      | N9316558      | N9316561      |
| 0.45    | 0.85    | 40 to 250        |              | N9316563      | N9316564      |               | N9316565      |
|         | 1.70    | 40 to 230/250    | N9316562     |               |               |               |               |
| 0.53    | 1.00    | 40 to 240/250    |              | N9316567      | N9316569      |               | N9316571      |
|         | 2.00    | 40 to 220/250    | N9316566     | N9316568      | N9316570      |               |               |

### Food and Flavor

#### Determination of omega-3 (n-3) and omega-6 (n-6) fatty acid composition in evening primrose oil and flax seed oil.

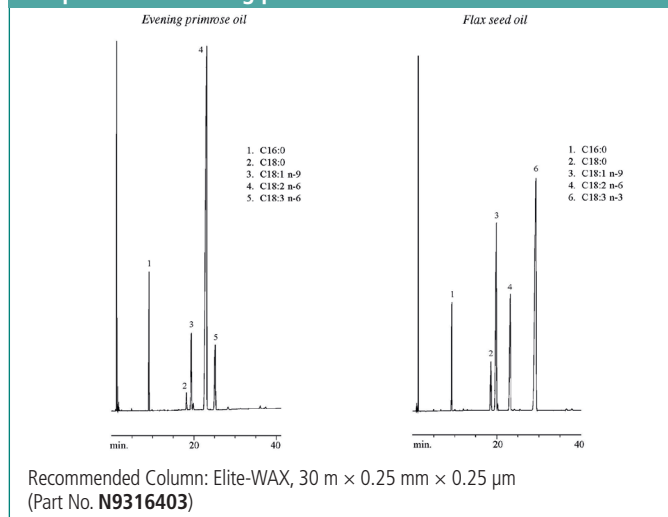

## GC Columns for GC/MS

The Elite range of MS columns are engineered for extremely low bleed for MS detectors, providing optimum sensitivity. They cover a wide range of polarities and applications.

### Elite-1ms

The Elite-1ms phase is a non-polar phase, (crosslinked dimethyl polysiloxane) designed to be robust for MS applications. With improved thermal stability and ultra low bleed provides increased sensitivity. It is regarded as a good general-purpose columns for arson accelerants, essential oils, hydrocarbons, pesticides, PCB congeners (e.g., Aroclor mixes), sulfur compounds, amines and solvent impurities.


#### Features

- Temperature range: -60 °C to 330/350 °C.
- Equivalent to USP G1, G2, and G38 phases

#### Applications

- Ideal for analysis of non polar petrochemical samples
- Also excellent for solvents, chemicals, flavors & fragrances, air toxins and pesticides

#### Elite-1ms Structure

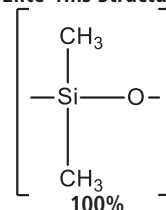

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.   | 20 m Part No.   | 30 m Part No.   | 60 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|-----------------|-----------------|
| 0.18    | 0.18    | -60 to 330/350   |                 | <b>N9305635</b> |                 |                 |
|         | 2.00    | -60 to 330/350   |                 | <b>N9305636</b> |                 |                 |
| 0.25    | 0.25    | -60 to 330/350   | <b>N9305637</b> |                 | <b>N9305638</b> | <b>N9305639</b> |
|         | 0.50    | -60 to 330/350   | <b>N9305640</b> |                 | <b>N9305641</b> | <b>N9305642</b> |
|         | 1.00    | -60 to 330/350   | <b>N9305643</b> |                 | <b>N9305644</b> | <b>N9305645</b> |
| 0.32    | 0.25    | -60 to 330/350   | <b>N9305646</b> |                 | <b>N9305647</b> | <b>N9305648</b> |
|         | 0.50    | -60 to 330/350   | <b>N9305649</b> |                 | <b>N9305650</b> | <b>N9305651</b> |
|         | 1.00    | -60 to 330/350   |                 |                 | <b>N9305652</b> | <b>N9305653</b> |
|         | 4.00    | -60 to 330/350   |                 |                 | <b>N9305654</b> |                 |

## Elite-5ms

The Elite-5ms phase (1,4-bis(dimethylsiloxy)phenylene dimethyl polysiloxane) incorporates a phenyl group in the polymer backbone to improve thermal stability, reduce bleed and make the phase less prone to oxidation. This results in a phase that is inert to active compounds with extremely low bleed to meet the requirements of sensitive MS detectors. It is a general purpose column ideal for GC/MS analysis of semivolatiles, PAHs, chlorinated hydrocarbons, phthalates, phenols, amines, organochlorine and organophosphorus pesticides, drugs and solvent impurities.


### Features

- Temperature Range: -60 °C to 350 °C
- Similar to USP G27 and G36 phases

### Applications

- Drugs, pesticides and solvent impurities
- Hydrocarbons and PCBs
- Essential oils and semivolatiles

Elite-5ms Structure

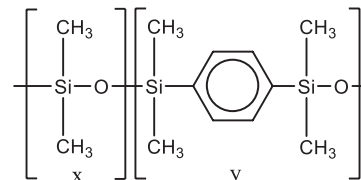

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.                | 30 m Part No.                | 60 m Part No.                |
|---------|---------|------------------|------------------------------|------------------------------|------------------------------|
| 0.18    | 0.18    | -60 to 325/340   |                              | <b>N9316276</b> <sup>1</sup> | <b>N9316277</b> <sup>1</sup> |
| 0.20    | 0.33    | -60 to 330/350   | <b>N9316301</b> <sup>2</sup> | <b>N9316302</b> <sup>2</sup> | <b>N9316303</b> <sup>2</sup> |
| 0.25    | 0.25    | -60 to 330/350   | <b>N9316279</b>              | <b>N9316282</b>              | <b>N9316286</b>              |
|         | 0.50    | -60 to 330/350   |                              | <b>N9316284</b>              |                              |
|         | 1.00    | -60 to 325/350   | <b>N9316280</b>              | <b>N9316283</b>              | <b>N9316287</b>              |
| 0.32    | 0.25    | -60 to 330/350   | <b>N9316289</b>              | <b>N9316293</b>              | <b>N9316297</b>              |
|         | 0.50    | -60 to 330/350   |                              | <b>N9316295</b>              |                              |
|         | 0.52    | -60 to 330/350   |                              | <b>N9316291</b> <sup>3</sup> |                              |
|         | 1.00    | -60 to 325/350   | <b>N9316290</b>              | <b>N9316294</b>              | <b>N9316298</b>              |
| 0.53    | 1.50    | -60 to 310/330   | <b>N9316299</b>              | <b>N9316300</b>              |                              |


<sup>1</sup> The lengths of **N9316276** and **N9316277** are 20 m and 40 m, respectively.

<sup>2</sup> The lengths of **N9316301**, **N9316302** and **N9316303** are 12 m, 25 m and 50 m, respectively.

<sup>3</sup> The length of **N9316291** is 25 m.

### Food and Flavor


The qualitative characterization of fruit juice flavor using a TurboMatrix Hs Trap and a Clarus SQ 8 GC/MS.

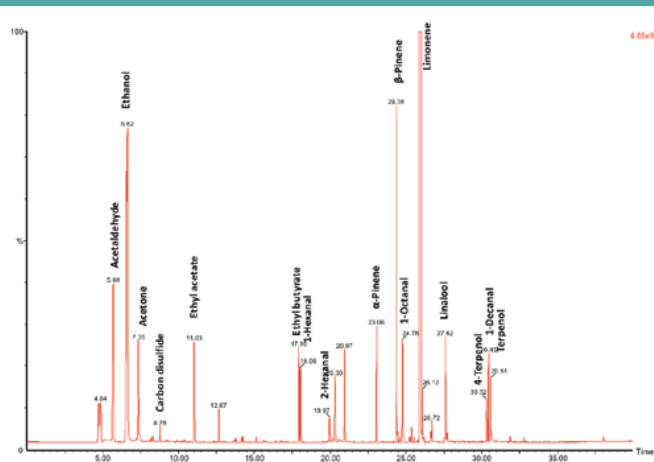

Recommended Column: Elite-5ms, 60 m x 0.25 mm x 1.0 μm, Part No. **N9316287**

### Food

The preparation and analysis of polycyclic aromatic hydrocarbons in meat by GC/MS.

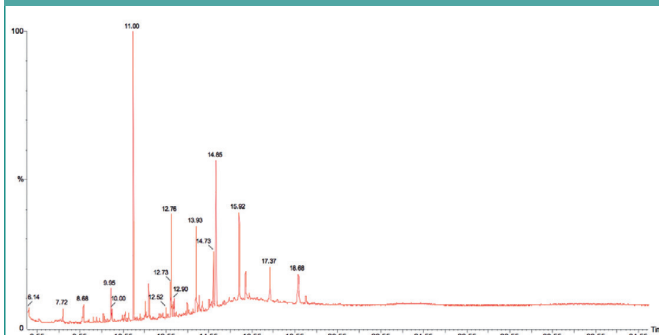

Column: Elite-5ms column, 30 m x 0.25 mm x 0.25 μm, Part No. **N9316282**

Liner Deactivated Liner, Part No. **N6502002**

## Elite-17ms


The Elite-17ms columns are general purpose, mid-polarity columns that are coated with a crosslinked, (50%-diphenyl)-dimethylpolysiloxane engineered for very low bleed to meet the requirements of sensitive MS detectors.

### Features

- Temperature Range: 40 °C to 300/340 °C
- Equivalent to USP G3 phase

### Applications

- Herbicides and pesticides
- Phthalate esters, sterols and rosin acids

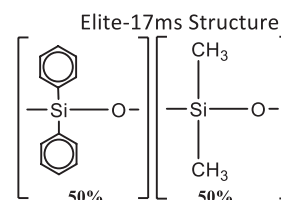

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.   | 30 m Part No.   | 60 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|-----------------|
| 0.18    | 0.18    | 40 to 300/340    | <b>N9316534</b> |                 |                 |
| 0.25    | 0.15    | 40 to 300/320    | <b>N9316535</b> | <b>N9316537</b> |                 |
|         | 0.25    | 40 to 300/320    | <b>N9316536</b> | <b>N9316538</b> | <b>N9316539</b> |
| 0.32    | 0.15    | 40 to 300/320    | <b>N9316540</b> | <b>N9316542</b> |                 |
|         | 0.25    | 40 to 300/320    | <b>N9316541</b> |                 |                 |

## Elite-35ms


The Elite-35ms columns are general purpose, mid-polarity columns offering extremely low bleed at higher temperatures. They are coated with a unique blend of linked dimethyl polysiloxanes and diphenyl polysiloxanes that are inert and selective for substituted polar compounds, such as drugs, pesticides, herbicides, PCBs and phenyls, while maintaining a similar selectivity and polarity as traditional Elite-35 phases.

### Features

- Temperature Range: 50 °C to 340/ 360 °C
- Equivalent to USP G42 phase

### Applications

- Pesticides and herbicides
- PCBs

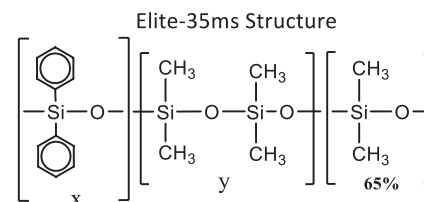

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.   | 30 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|
| 0.25    | 0.25    | 50 to 340/360    | <b>N9305686</b> | <b>N9305687</b> |
|         | 0.50    | 50 to 340/360    | <b>N9305688</b> | <b>N9305689</b> |
|         | 1.00    | 50 to 320/340    | <b>N9305690</b> | <b>N9305691</b> |
| 0.32    | 0.25    | 40 to 340/360    | <b>N9305692</b> | <b>N9305693</b> |
|         | 0.50    | 40 to 340/360    | <b>N9305694</b> | <b>N9305695</b> |
|         | 1.00    | 40 to 320/340    | <b>N9305696</b> | <b>N9305697</b> |

# Elite-624ms

The Elite-624ms incorporates a unique proprietary blend of cyanopropyl and methyl siloxanes that results in a very inert, extremely low bleed and high thermal stability column. This column provides excellent peak shape for a wide range of compounds and is highly selective for residual solvents making it a great choice for USP<467>. These columns are manufactured for column-to-column reproducibility, so they are well suited for validated methods.


## Features

- Temperature Range: -20 °C to 300/320 °C
- Similar to USP G43 phase
- Ideal choice for USP method 467

## Applications

- Residual solvents
- Suitable for USP 467


Elite-624ms Structure

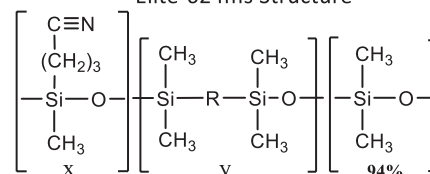

| ID (mm) | df (μm) | Temp Limits (°C) | 20 m Part No.   | 30 m Part No.   | 60 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|-----------------|
| 0.18    | 1.00    | -20 to 300/320   | <b>N9315067</b> |                 |                 |
| 0.25    | 1.40    | -20 to 300/320   |                 | <b>N9315068</b> | <b>N9315066</b> |
| 0.32    | 1.80    | -20 to 300/320   |                 | <b>N9315069</b> | <b>N9315070</b> |

## Industrial

### Determination of nine carbonates in lithium ion battery electrolyte by GC/MS.

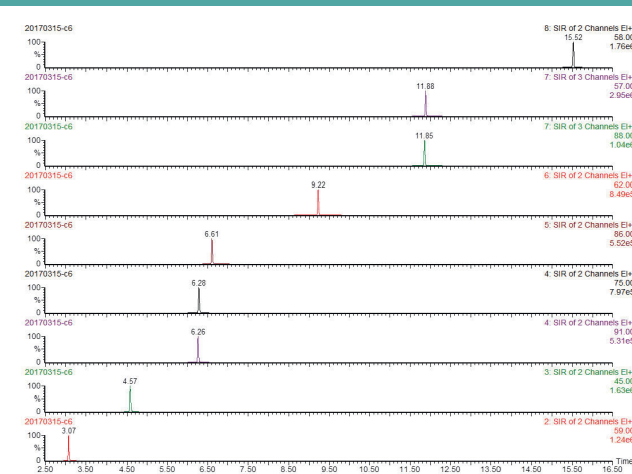

Recommended Column: Elite-35ms, 30 m x 0.25 mm x 0.25 μm (Part No. **N9316438**)

Liner: Capillary splitless deactivated glass liners with deactivated wool (**N9306235**)

## High Temperature Columns

Available in a range of phases with varying polarity, the high temperature (ht) columns are specifically designed for reduced bleed when operating at higher temperatures, up to 400 °C. The optimum higher operating temperatures varies by phase.

### Elite-1ht

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.   | 30 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|
| 0.25    | 0.10    | -60 to 380/400   | <b>N9316268</b> | <b>N9316269</b> |
| 0.32    | 0.10    | -60 to 380/400   | <b>N9316270</b> | <b>N9316271</b> |

### Elite-17ht

| ID (mm) | df (μm) | Temp Limits (°C) | 30 m Part No.   |
|---------|---------|------------------|-----------------|
| 0.25    | 0.15    | 40 to 300/320    | <b>N9316264</b> |
| 0.32    | 0.15    | 40 to 300/320    | <b>N9316266</b> |

### Elite-5ht

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.   | 30 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|
| 0.25    | 0.10    | -60 to 400       | <b>N9316272</b> | <b>N9316273</b> |
| 0.32    | 0.10    | -60 to 400       | <b>N9316274</b> | <b>N9316275</b> |

### Elite-SimDist ht

Application: High-temperature simulated distillation

Phase: Metal Column, 100% dimethylpolysiloxane, non-polar

| ID (mm) | df (μm) | Temp Limits (°C) | 5 m Part No.    | 6 m Part No.    |
|---------|---------|------------------|-----------------|-----------------|
| 0.53    | 0.15    | -60 to 400       |                 | <b>N6107191</b> |
| 0.53    | 0.10    | -60 to 450       | <b>NR213314</b> |                 |

## Metal High Temperature Columns

Using our new metal capillary columns eliminates the risk of column breakage at higher operating temperatures. A range of phases are offered, covering common applications. Rugged up to 450 °C although the exact upper temperature limits varies depending on phase and column configurations.

|                      | Temp Limits (°C) | Inner Diameter (mm) | df (μm) | 15 m Part No.   | 30 m Part No.   | 60 m Part No.   |
|----------------------|------------------|---------------------|---------|-----------------|-----------------|-----------------|
| <b>Elite-1mht</b>    | -60 to 340/430   | 0.25                | 0.1     |                 | <b>N9303453</b> |                 |
| <b>Elite-1mht</b>    | -60 to 340/430   | 0.25                | 0.25    | <b>N9303454</b> | <b>N9303455</b> | <b>N9303456</b> |
| <b>Elite-5mht</b>    | -60 to 330/430   | 0.25                | 0.1     | <b>N9303457</b> | <b>N9303458</b> |                 |
| <b>Elite-5mht</b>    | -60 to 330/430   | 0.25                | 0.25    | <b>N9303459</b> | <b>N9303460</b> | <b>N9303461</b> |
| <b>Elite-1301mht</b> | -20 to 280       | 0.53                | 3.0     | <b>N9303462</b> | <b>N9303463</b> |                 |
| <b>Elite-1701mht</b> | -20 to 260       | 0.53                | 1.0     | <b>N9303464</b> | <b>N9303465</b> |                 |
| <b>Elite-1701mht</b> | -20 to 250       | 0.53                | 1.5     |                 | <b>N9303466</b> |                 |

# Elite-BAC Advantage: Blood Alcohol Content


The Elite-BAC Advantage columns are optimized for selectivities guaranteed to resolve ethanol, internal standards, and frequently encountered interferences. These application-specific columns for blood alcohol analysis baseline separate all critical compounds, including ethanol, methanol, acetone, *tert*-butanol, acetaldehyde, isopropanol, and *n*-propanol, in less than 2 minutes. Every Elite-BAC 1 Advantage and Elite-BAC 2 Advantage column is qualified with a test mix containing these important BAC target compounds to ensure reproducibility. These columns, baseline separate all blood alcohol compounds in blood, breath, or urine, in less than 2 minutes, under isothermal conditions. Isothermal analysis increases productivity by eliminating the need for oven cycling. Confirmation is easily achieved with this tandem set because there are two elution order changes between the columns.

## Features and Benefits

- Robust and reproducible
- Baseline separation of all components in less than 2 minutes
- Stable to 260 °C


## Clinical

### Elite-BAC 1 advantage.

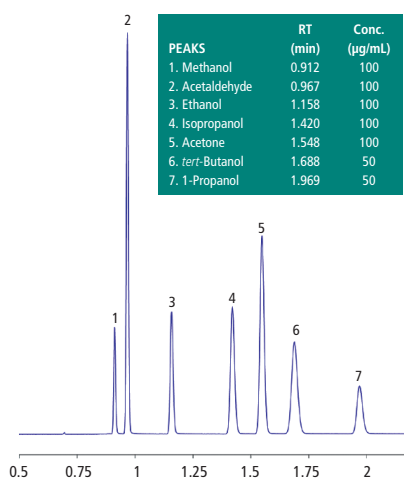

## Clinical

### Elite-BAC 2 advantage.

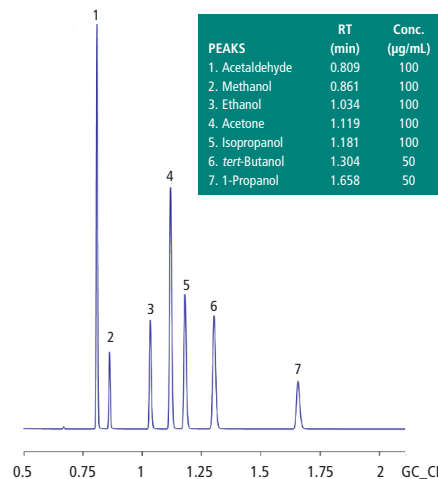

- Baseline resolution
- Symmetrical peaks
- Run time under 2 min

| Column Type           | ID (mm) | df (µm) | Temp Limits (°C) | 10 m Part No. | 30 m Part No. |
|-----------------------|---------|---------|------------------|---------------|---------------|
| Elite-BAC 1 Advantage | 0.18    | 1.00    | -20 to 240/260   | N9315075      |               |
|                       | 0.32    | 1.80    | -20 to 240/260   |               | N9315071      |
|                       | 0.53    | 3.00    | -20 to 240/260   |               | N9315072      |
| Elite-BAC 2 Advantage | 0.18    | 0.34    | -20 to 240/260   | N9315076      |               |
|                       | 0.32    | 0.60    | -20 to 240/260   |               | N9315073      |
|                       | 0.53    | 1.00    | -20 to 240/260   |               | N9315074      |

## Elite-VMS

Elite-VMS columns offer lower bleed, better selectivity, and overall faster analysis for separating volatile organic compounds. The stationary phase is a highly stable polymer that provides outstanding analysis of volatile compounds on MS detectors. The 0.18 and 0.25 mm ID columns allow sample splitting at the injection port, eliminating the added expense and maintenance of a jet separator. A 0.45 mm or 0.53 mm ID column can be directly connected to the purge-and-trap transfer line in a system equipped with a jet separator.

| ID (mm) | df (μm) | Temp Limits (°C) | 30 m Part No.               | 60 m Part No.               |
|---------|---------|------------------|-----------------------------|-----------------------------|
| 0.18    | 1.00    | -40 to 240/260   | <b>N9316650<sup>1</sup></b> | <b>N9316651<sup>1</sup></b> |
| 0.25    | 1.40    | -40 to 240/260   | <b>N9316652</b>             | <b>N9316653</b>             |
| 0.32    | 1.80    | -40 to 240/260   | <b>N9316654</b>             | <b>N9316655</b>             |
| 0.45    | 2.55    | -40 to 240/260   | <b>N9316656</b>             | <b>N9316657</b>             |
| 0.53    | 3.00    | -40 to 240/260   | <b>N9316658</b>             | <b>N9316659</b>             |

<sup>1</sup> The lengths of N9316650 and N9316651 are 20 m and 40 m, respectively

### Features and Benefits

- Temperature Range: -40 °C to 240/260 °C
- No known equivalent phases

### Applications

- Ideal for analysis of volatile organic pollutants by GC/MS
- Suitable for EPA Method 8260B

## Elite-XLB

The Elite-XLB phase is a proprietary low-polarity, very inert and exceptionally low bleed column for GC/MS analysis of pesticides, PCB congeners (e.g., Aroclor mixes) and PAHs. Improvements in polymer synthesis and tubing deactivation enable us to make inert, stable Elite-XLB columns especially well-suited for analyzing active, high molecular weight compounds with sensitive GC-MS systems, including ion trap detectors.

### Features and Benefits

- Temperature Range: 30 °C to 340/360 °C
- No known equivalent phases
- Exceptionally low bleed for GC/MS

### Applications

- Pesticides, PCB congeners
- Semi volatiles in drinking water
- Suitable for EPA Method 525

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.               | 30 m Part No.               | 60 m Part No.   |
|---------|---------|------------------|-----------------------------|-----------------------------|-----------------|
| 0.18    | 0.18    | 40 to 340/360    |                             | <b>N9316480<sup>1</sup></b> |                 |
| 0.20    | 0.33    | 40 to 340/360    | <b>N9316496<sup>2</sup></b> | <b>N9316497<sup>2</sup></b> |                 |
| 0.25    | 0.10    | 40 to 340/360    |                             | <b>N9316483</b>             |                 |
|         | 0.25    | 40 to 340/360    | <b>N9316481</b>             | <b>N9316484</b>             | <b>N9316487</b> |
|         | 1.00    | 40 to 340/360    | <b>N9318482</b>             | <b>N9316485</b>             |                 |
| 0.32    | 0.10    | 40 to 340/360    |                             | <b>N9316489</b>             |                 |
|         | 0.25    | 40 to 340/360    | <b>N9316488</b>             | <b>N9316490</b>             | <b>N9316493</b> |
|         | 0.50    | 40 to 340/360    |                             | <b>N9316492</b>             |                 |
|         | 1.00    | 40 to 340/360    |                             | <b>N9316491</b>             |                 |
| 0.53    | 1.50    | 40 to 320/340    | <b>N9316494</b>             | <b>N9316495</b>             |                 |

<sup>1</sup> The length of N9316480 is 20 m.

<sup>2</sup> The lengths of N9316496 and N9316497 are 12 m and 25 m, respectively.

## Elite-Volatiles

The Elite-Volatiles stationary phase and optimized column dimensions provide low bleed, excellent resolution, and fast analysis times for volatile organic pollutants.

### Features

- Temperature Range: -20 °C to 240 °C
- Proprietary phase
- Ideal for EPA Method 8021

### Applications

- Volatile organic pollutants
- Suitable for EPA Method 8021

| ID (mm) | df (μm) | Temp Limits (°C) | 30 m Part No.   | 60 m Part No.   | 75 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|-----------------|
| 0.25    | 1.40    | -20 to 240       | <b>N9316388</b> | <b>N9316389</b> |                 |
| 0.32    | 1.80    | -20 to 240       | <b>N9316390</b> | <b>N9316391</b> |                 |
| 0.45    | 2.55    | -20 to 240       | <b>N9316392</b> |                 | <b>N9316393</b> |

## Elite-CLPesticides: Chlorinated Pesticides

Elite-CLPesticides is specially designed to overcome the coelutions and analyte breakdown typically encountered in chlorinated pesticide analytes for U.S. EPA methods 8081, 608, and CLP. Column bleed measured by ECD is extremely low at temperatures greater than 300 °C, which is critical for baking out the column to remove high-boiling compounds commonly found in pesticide/PCB extracts.

### Features and Benefits

- Thermally stable to 340 °C
- Low column bleed – ideal for ECD or GC/MS analysis
- Exceeds performance criteria for U.S. EPA Methods 8081, 608 and CLP
- Baseline separation in less than 15 minutes

### Applications

- Chlorinated pesticides and herbicides
- Exceeds performance criteria for U.S. EPA Methods 8081, 608 and CLP

| Column Type          | ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.   | 30 m Part No.   |
|----------------------|---------|---------|------------------|-----------------|-----------------|
| Elite-CLPesticides   | 0.25    | 0.25    | -60 to 320/340   | <b>N9316661</b> | <b>N9316662</b> |
|                      | 0.32    | 0.50    | -60 to 320/340   | <b>N9316663</b> | <b>N9316664</b> |
|                      | 0.53    | 0.50    | -60 to 300/320   | <b>N9316665</b> | <b>N9316666</b> |
| Elite-CLPesticides 2 | 0.25    | 0.20    | -20 to 240/260   | <b>N9316667</b> | <b>N9316668</b> |
|                      | 0.32    | 0.25    | -20 to 240/260   | <b>N9316669</b> | <b>N9316670</b> |
|                      | 0.53    | 0.42    | -20 to 240/260   | <b>N9316671</b> | <b>N9316672</b> |

## Elite-502.2: U.S. EPA Method 502.2

**Application:** Analysis of volatiles by U.S. EPA method 502.2

**Phase:** Proprietary Dimethyl-diphenyl polysiloxane, low-polarity

| ID (mm) | df (μm) | Temp Limits (°C) | 60 m Part No. | 75 m Part No. | 105 m Part No. |
|---------|---------|------------------|---------------|---------------|----------------|
| 0.25    | 1.40    | 0 to 250/270     | N9316498      |               |                |
| 0.45    | 2.55    | 0 to 250/270     |               | N9316188      | N9316189       |
| 0.53    | 3.00    | 0 to 250/270     |               |               | N9316190       |

## Elite-RX: Drugs of Abuse

**Application:** Analysis of drugs of abuse

| Phase        | ID (mm) | df (μm) | Temp Limits (°C) | 12 m Part No. | 25 m Part No. |
|--------------|---------|---------|------------------|---------------|---------------|
| Elite-1 RX   | 0.20    | 0.33    | -60 to 330/350   | N9316345      | N9316346      |
| Elite-5ms RX | 0.20    | 0.33    | -60 to 330/350   | N9316349      | N9316350      |
| Elite-17 RX  | 0.20    | 0.33    | 40 to 300/320    | N9316347      | N9316348      |

## Elite-Betecyclodextrin: Chiral Separations

**Application:** General-purpose chiral, Chiral compounds in essential oils

| Column Type      | ID (mm) | df (μm) | Temp Limits (°C) | 30 m Part No. |
|------------------|---------|---------|------------------|---------------|
| Elite-Betacydex  | 0.25    | 0.25    | 40 to 230        | N9316319      |
| Elite-Cyclosil B | 0.25    | 0.25    | 40 to 230        | N9316545      |

## Elite-SimDist

**Application:** Simulated distillation

**Phase:** Specially processed dimethylpolysiloxane, non-polar

| ID (mm) | df (μm) | Temp Limits (°C) | 10 m Part No. |
|---------|---------|------------------|---------------|
| 0.45    | 2.55    | -60 to 360       | N9316261      |
| 0.53    | 3.00    | -60 to 360       | N9316262      |

## Elite-608

**Application:** Analysis of semivolatile pesticides by U.S. EPA method 608

**Phase:** Phenyl methyl polysiloxane, mid-polarity

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No. | 30 m Part No. |
|---------|---------|------------------|---------------|---------------|
| 0.32    | 0.50    | 40 to 290/310    |               | N9316191      |
| 0.45    | 0.42    | 40 to 270/290    | N9316194      | N9316195      |
|         | 0.70    | 40 to 260/280    | N9316192      | N9316193      |
| 0.53    | 0.50    | 40 to 270/290    | N9316198      | N9316199      |
|         | 0.83    | 40 to 260/280    | N9316196      | N9316197      |

## Elite-TPH

**Application:** Analysis of total petroleum hydrocarbons

**Phase:** (5%-diphenyl)-dimethylpolysiloxane, low polarity

| ID (mm) | df (μm) | Temp Limits (°C) | 30 m Part No.   |
|---------|---------|------------------|-----------------|
| 0.32    | 0.25    | -10 to 320       | <b>N9316386</b> |
| 0.45    | 1.00    | -10 to 290       | <b>N9316387</b> |

## Elite-PONA

**Application:** Detailed analysis of petroleum naphtha

**Phase:** Specially processed dimethylpolysiloxane, non-polar

| ID (mm) | df (μm) | Temp Limits (°C) | 50 m Part No.   | 100 m Part No.  |
|---------|---------|------------------|-----------------|-----------------|
| 0.20    | 0.50    | -60 to 300/320   | <b>N9316065</b> |                 |
| 0.25    | 0.50    | -60 to 300/320   |                 | <b>N9316015</b> |

## Elite-FFAP

**Application:** Free fatty acids

**Phase:** Nitroterephthalic acid modified PEG (bonded), polar

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.   | 30 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|
| 0.25    | 0.25    | 40 to 250        | <b>N9316351</b> | <b>N9316352</b> |
| 0.32    | 0.25    | 40 to 250        | <b>N9316353</b> | <b>N9316354</b> |
| 0.45    | 0.85    | 40 to 240/250    | <b>N9316355</b> | <b>N9316356</b> |
| 0.53    | 1.00    | 40 to 240/250    | <b>N9316357</b> | <b>N9316358</b> |

## Elite-5 Amine

**Application:** Amines and other basic compounds including alkylamines and di/triamines

| ID (mm) | df (μm) | Temp Limits (°C) | 15 m Part No.   | 30 m Part No.   |
|---------|---------|------------------|-----------------|-----------------|
| 0.25    | 0.50    | -60 to 300/315   | <b>N9316684</b> | <b>N9316673</b> |
|         | 1.00    | -60 to 300/315   | <b>N9316674</b> | <b>N9316675</b> |
| 0.32    | 1.00    | -60 to 300/315   | <b>N9316676</b> | <b>N9316677</b> |
|         | 1.50    | -60 to 290/305   | <b>N9316678</b> | <b>N9316679</b> |
| 0.53    | 1.00    | -60 to 290/305   |                 | <b>N9316680</b> |
|         | 3.00    | -60 to 280/295   | <b>N9316681</b> | <b>N9316682</b> |

## Elite-2330/Elite-23

**Application:** Analysis of cis/trans isomers in FAMES and dioxin isomers. Equivalent to USP G8 and G48

**Phase:** Biscyanopropyl cyanopropylphenyl polysiloxane, highly polar

| ID (mm) | df (μm) | Temp Limits (°C) | 60 m Part No.   | Column Type |
|---------|---------|------------------|-----------------|-------------|
| 0.25    | 0.10    | 0 to 275         | <b>N6107813</b> | Elite-2330  |
| 0.25    | 0.20    | 0 to 275         | <b>N6107814</b> | Elite-2330  |
| 0.25    | 0.25    | 40 to 250/260    | <b>N9316508</b> | Elite-23    |

## Elite-MTBE

**Application:** Analysis of methyl t-butylether and other oxygenates

**Phase:** Proprietary low polarity phase

| ID (mm) | df (μm) | Temp Limits (°C) | 30 m Part No.   |
|---------|---------|------------------|-----------------|
| 0.45    | 2.55    | 10 to 250        | <b>N9316520</b> |
| 0.53    | 3.00    | 10 to 250        | <b>N9316521</b> |

## Elite-2560

**Application:** Application-specific column for cis/trans FAMES

**Phase:** Biscyanopropylpolysiloxane, highly polar

| ID (mm) | df (μm) | Temp Limits (°C) | 100 m Part No.  |
|---------|---------|------------------|-----------------|
| 0.25    | 0.20    | 20 to 250        | <b>N9311570</b> |

## Elite-Carbon Columns

### For Volatiles in Hydrocarbon Streams

The Elite-Carbon columns offer rapid separation of permanent gas/light hydrocarbon mixtures; including carbon monoxide and carbon dioxide without cryogenic cooling. They are preconditioned and thus take less than 30 minutes to stabilize. They are used in conjunction with a molecular sieve column (Molecular sieve 5 Å, 50 m, 0.53 mm, 50 μm Part No. **NR201108**).

| ID (mm) | Length (m) | Mesh Size | Temp Limits (°C) | Part No.        |
|---------|------------|-----------|------------------|-----------------|
| 1.0     | 1          | 100/120   | Up to 300        | <b>N9303927</b> |
| 1.0     | 2          | 100/120   | Up to 300        | <b>N9303926</b> |

Fittings for the micropacked Elite-Carbon columns need to be ordered separately.

| Description                                                  | Part No.        |
|--------------------------------------------------------------|-----------------|
| Installation kit for 1 mm ID columns; for valve applications | <b>N9303450</b> |
| Installation kit for 1 mm ID columns; for direct injections  | <b>N9303451</b> |

## Elite-Alumina/KCl\* PLOT

| ID (mm) | Length (m) | Film Thickness (μm) | Temp Limits (°C) | Part No.        |
|---------|------------|---------------------|------------------|-----------------|
| 0.53    | 50         | 10                  | up to 200        | <b>N9316544</b> |

\* Lower Polarity than Elite-Alumina.

Elite-Alumina/Na<sub>2</sub>SO<sub>4</sub> PLOT

| ID (mm) | Length (m) | Film Thickness (μm) | Temp Limits (°C) | Part No.        |
|---------|------------|---------------------|------------------|-----------------|
| 0.32    | 50         | 5                   | up to 200        | <b>N6107777</b> |

## Elite-Alumina PLOT

## Phase for Analysis of Light Hydrocarbons

| ID (mm) | Length (m) | Film Thickness (μm) | Temp Limits (°C) | Part No.        |
|---------|------------|---------------------|------------------|-----------------|
| 0.53    | 30         | 6                   | -60 to 200       | <b>N9316304</b> |
| 0.53    | 50         | 10                  | -60 to 200       | <b>N9316305</b> |

Note: -60 °C is the lowest temperature used on this phase in our lab. Lower temperatures may be used depending on the sample.

## Elite-Cyclosil B PLOT

## For Chiral Separations

| ID (mm) | Length (m) | Film Thickness (μm) | Temp Limits (°C) | Part No.        |
|---------|------------|---------------------|------------------|-----------------|
| 0.25    | 30         | 0.25                | 35 to 230        | <b>N9316545</b> |
| 0.32    | 30         | 0.25                | 35 to 230        | <b>N9316546</b> |

## Elite-Molesieve PLOT

## Phase for Analysis of Permanent Gases

| ID (mm) | Length (m) | Film Thickness (μm) | Temp Limits (°C) | Part No.        |
|---------|------------|---------------------|------------------|-----------------|
| 0.53    | 30         | —                   | -60 to 300       | <b>N9316361</b> |

Note: -60 °C is the lowest temperature used on this phase in our lab. Lower temperatures may be used depending on the sample.

## Elite-Q PLOT

## Phase for Analysis of Light Gases and Hydrocarbons

| ID (mm) | Length (m) | Film Thickness (μm) | Temp Limits (°C) | Part No.        |
|---------|------------|---------------------|------------------|-----------------|
| 0.32    | 30         | 10                  | -60 to 250       | <b>N9316359</b> |
| 0.53    | 30         | 20                  | -60 to 250       | <b>N9316360</b> |

Note: -60 °C is the lowest temperature used on this phase in our lab. Lower temperatures may be used depending on the sample.

## Miscellaneous Accessories

| Description                                                                                                                  | Part No.        |
|------------------------------------------------------------------------------------------------------------------------------|-----------------|
| 2 oz. Replacement Charcoal (30/60 mesh)                                                                                      | <b>03300904</b> |
| Liner Removal Tool                                                                                                           | <b>N6100102</b> |
| Injector/Detector Adapter 1/4 in.<br>Adapter fits on injector and detector outlet (inside oven) for use with 1/4 in. columns | <b>00080100</b> |
| Silanized Glass Wool (2 oz.)                                                                                                 | <b>03300905</b> |

## Wide-Bore Adapter Kit

Contains all the parts necessary to adapt to packed column injectors quickly and easily for use with wide-bore capillary columns. Includes 0–20 mL/min flow controller element, wide-bore adapter with 1/16 in. fitting, wide-bore glass liner and column support hanger.



| Description                       | Part No.        |
|-----------------------------------|-----------------|
| 0.53 Capillary Column Adapter Kit | <b>N6120001</b> |

## Wafer Scribes

The PerkinElmer ceramic wafer scribe is inexpensive and ideal for cutting polyimide fused silica capillary columns and guard columns. The scribe is easy to hold and simple to use. All four sides can be used as a cutting tool.



| Description             | Part No.        |
|-------------------------|-----------------|
| Wafer Scribes (pkg. 10) | <b>N9301376</b> |

## Connectors

## Universal Connectors



| Description                                                  | Part No.        |
|--------------------------------------------------------------|-----------------|
| Deactivated Presstight Column Connectors (pkg. 5)            | <b>N9302149</b> |
| Metal Universal Connectors: 0.25 mm ID (pkg. 10)             | <b>N9301167</b> |
| Deactivated Y Splitter Presstight Column Connectors (pkg. 3) | <b>N9306380</b> |
| Polyimide Sealing Resin (5 g)                                | <b>N9301343</b> |

## New Investments – Local Support – Proven Quality!

For over 80 years, PerkinElmer has helped our global laboratory partners meet their most demanding application measurement challenges. Our experience in Chromatography, Spectroscopy and Materials Characterization gives us a unique advantage in providing holistic and comprehensive solutions for laboratories in universities and key market segments like Environmental, Energy, and Food and Product Safety.

When you order our precision designed, genuine PerkinElmer consumables and accessories, you can enjoy the peace of mind, ease of ordering and best in class service. You'll get the results you need – accurately and on time. PerkinElmer selects only the highest quality products for our consumables portfolio and invests heavily to test and validate the quality of our products.

| THE AMERICAS       |                   |                          |                                         |
|--------------------|-------------------|--------------------------|-----------------------------------------|
| Country            | Telephone         | Facsimile                | Email                                   |
| ARGENTINA          | 54 11 4554 4004   | 54 11 4554 2807          | consumiblesarg@perkinelmer.com          |
| BELIZE             | 506 2 2724700     | 506 2 2729700            | info@inanbio.co.cr                      |
| BOLIVIA            | (591 2) 241 5446  | (591 2) 241 0639         | consumiblesarg@perkinelmer.com          |
| BRAZIL             | 55 11 3868 6200   | 55 11 3868 6210          | consumiveis@perkinelmer.com             |
| CANADA             | 800 561 4646      | 905 851 1814             | Orders.LSCanada@perkinelmer.com         |
| CHILE              | 56 2 784 7500     | 56 2 236 6821            | consumibleschile@perkinelmer.com        |
| COLUMBIA           | 57 320 2670667    |                          | ventas@ccvgrupo.com                     |
| COSTA RICA         | (506) 2 272 4700  | (506) 2 272 9700         | info@inanbio.co.cr                      |
| DOMINICAN REPUBLIC | (809) 508 6500    | (809) 532 0741           |                                         |
| ECUADOR            | (593) 2 2456 556  | (593) 2 2456 556 ext 102 | info@proinstra.com                      |
| EL SALVADOR        | (503) 2535 0700   | (503) 2535 0777          | coresa@coresaelsalvador.com             |
| GUATEMALA          | (502) 25 02 47 00 | (502) 25 02 47 01        |                                         |
| HONDURAS           | (504) 9982-3628   |                          | Lvallejoh@rhandals.hn                   |
| MEXICO             | 52 55 5337 0700   | 52 55 559 36223          | aftermarket.consumables@perkinelmer.com |
| NICARAGUA          | (505) 2 278 0177  | (505) 278 0279           |                                         |
| PANAMA             | 507 3941 695      | 507 3941699              | Pjroyer@comerrsa.com                    |
| PARAGUAY           | +595 21 228 068   | +595 21 228 223          | Insitu@insitu.com.py                    |
| PERU               | 51-989287824      |                          | ciansac@cientifica-andina.com.pe        |
| PUERTO RICO        | 787 783 3067      | 787 793 7198             |                                         |
| TRINIDAD           | (868) 671 1525    | (868) 665 3703           |                                         |
| UNITED STATES      | 800 762 4000      | 203 944 4904             | CustomerCareUS@perkinelmer.com          |
| URUGUAY            | 59 8 94385747     |                          | gpastoriza@biodiagnostico.com.uy        |
| VENEZUELA          | 58 212 7614757    |                          | ventas@ccvgrupo.com                     |

# ORDERING INFORMATION

| EUROPE             |                   |                   |                            |
|--------------------|-------------------|-------------------|----------------------------|
| Country            | Telephone         | Facsimile         | Email                      |
| ALBANIA            | 0800 111 933      | 0800 111 934      | cc.export@perkinelmer.com  |
| AUSTRIA            | 0800 111 933      | 0800 111 934      | cc.austria@perkinelmer.com |
| BELARUS            | +7 495 935 8888   | +7 495 564 8787   | info@scheltec.ru           |
| BELGIUM            | 0800 40858        | 0800 40859        | cc.benelux@perkinelmer.com |
| BOSNIA HERZEGOVINA | +385 1 2001 767   | +385 1 2001 724   | hebe@hebe.hr               |
| BULGARIA           | +359 2 952 0301   | +359 2 952 3784   | chroma@spnet.net           |
| CROATIA            | +385 1 2001 767   | +385 1 2001 724   | hebe@hebe.hr               |
| CYPRUS             | +357 22 33 82 34  | +357 22 33 82 32  | sales@kriticos.com.cy      |
| CZECH REPUBLIC     | +420 241 430 534  | +420 241 430 535  | info@pesystems.cz          |
| DENMARK            | 8088 4236         | 8088 4237         | cc.nordic@perkinelmer.com  |
| ESTONIA            | +372 7 383 330    | +372 7 383 360    | kent@lanlab.ee             |
| FINLAND            | 0800 117 186      | 0800 117 185      | cc.nordic@perkinelmer.com  |
| FRANCE             | 0805 111 333      | 0805 111 334      | cc.france@perkinelmer.com  |
| GERMANY            | 0800 181 0032     | 0800 181 0031     | cc.germany@perkinelmer.com |
| GREECE             | +30 2310 322525   | +30 2310 321912   | antisel@antisel.gr         |
| HUNGARY            | +36 1 251 1116    | +36 1 251 1461    | perform@per-form.hu        |
| ICELAND            | +354 412 7000     | +354 412 7099     | medor@medor.is             |
| IRELAND            | 1800 932 886      | 1800 932 884      | cc.ie@perkinelmer.com      |
| ISRAEL             | 1700 500 512      | 1700 555 915      | cc.israel@perkinelmer.com  |
| ITALY              | +39 800 906 642   | +39 02 3601 2508  | cc.italy@perkinelmer.com   |
| LATVIA             | +371 6755 1894    | +371 6755 1976    | info@adrona.lv             |
| LITHUANIA          | +370 5 2638748    | +370 5 2638749    | info@linealibera.it        |
| LUXEMBOURG         | 0800 26588        | 0800 26589        | cc.benelux@perkinelmer.com |
| MACEDONIA          | +381 11 22 22 244 | +381 11 22 22 244 | perkinelmer@superlab.rs    |
| MALTA              | +39 081 7013094   | +39 081 7021977   | info@areachem.it           |
| MOLDAVIA           | +40 21 316 2877   | +40 21 316 2876   | office@cromatec.ro         |
| MONTENEGRO         | +381 11 22 22 244 | +381 11 22 22 244 | perkinelmer@superlab.rs    |
| NETHERLANDS        | 0800 0234490      | 0800 0234491      | cc.benelux@perkinelmer.com |
| NORWAY             | 800 18 854        | 800 18 855        | cc.nordic@perkinelmer.com  |
| POLAND             | 0800 88 08 88     | +48 22 310 88 01  | cc.pl@perkinelmer.com      |
| PORTUGAL           | +351 21 782 60 30 | +351 21 782 60 39 | ilc@ilc.pt                 |
| ROMANIA            | +40 21 316 2877   | +40 21 316 2876   | office@cromatec.ro         |

# ORDERING INFORMATION

## EUROPE

| Country                   | Telephone            | Facsimile         | Email                             |
|---------------------------|----------------------|-------------------|-----------------------------------|
| <b>RUSSIAN FEDERATION</b> | +7 495 935 8888      | +7 495 564 8787   | info@scheltec.ru                  |
| <b>SERBIA</b>             | +381 11 22 22 244    | +381 11 22 22 244 | perkinelmer@superlab.rs           |
| <b>SLOVAKIA</b>           | +421 918 506782      | +421 48 2301628   | mLa@mediclabor.sk                 |
| <b>SLOVENIA</b>           | +386 1 427 3213/3650 | +386 1 427 3191   | omega@omega.si                    |
| <b>SPAIN</b>              | 800 099 164          | 800 099 165       | atencionalcliente@perkinelmer.com |
| <b>SWEDEN</b>             | 0200 887 520         | 0200 887 521      | cc.nordic@perkinelmer.com         |
| <b>SWITZERLAND</b>        | 0800 000 015         | 0800 000 016      | cc.switzerland@perkinelmer.com    |
| <b>TURKEY</b>             | +90 212 312 1100     | +90 212 312 1111  | cc.turkey@perkinelmer.com         |
| <b>UKRAINE</b>            | +7 495 935 8888      | +7 495 564 8787   | info@scheltec.ru                  |
| <b>UNITED KINGDOM</b>     | 0800 89 60 46        | 0800 89 17 14     | cc.uk@perkinelmer.com             |

## ASIA

| Country               | Telephone                       | Facsimile        | Email                                    |
|-----------------------|---------------------------------|------------------|------------------------------------------|
| <b>AFGANISTAN</b>     | +92 21 5345581                  | +92 21 5345582   | sales@amspvtltd.com                      |
| <b>ARMENIA</b>        | +7 495 935 8888                 | +7 495 564 8787  | info@scheltec.ru                         |
| <b>AUSTRALIA</b>      | 61 3 9212 8500 or 1 800 033 391 | 61 3 9212 8595   | www.perkinelmer.com/category/consumables |
| <b>AZERBAIJAN</b>     | +90 212 212 5566                | +90 212 212 2829 | istanbul@tetratex.com.tr                 |
| <b>CHINA</b>          | 86 8008205046                   | 86 2160645666    | consumable.china@perkinelmer.com         |
| <b>GEORGIA</b>        | +7 495 935 8888                 | +7 495 564 8787  | info@scheltec.ru                         |
| <b>HONG KONG</b>      | 852 2793 3218                   | 852 2793 3225    | info@aldenor.com.uy                      |
| <b>INDIA</b>          | +91 22 33261700                 | +91 22 33261790  | marketing.india@perkinelmer.com          |
| <b>INDONESIA</b>      | 62 21 5694 5040 to 43           | 62 21 5694 5024  |                                          |
| <b>JAPAN</b> Osaka    | +81 6 6394 5360                 | +81 6 6394 5363  |                                          |
| <b>JAPAN</b> Tokyo    | +81 3 3866 2647                 | +81 3 3866 2652  |                                          |
| <b>JAPAN</b> Yokahama | +81 45 339 5866                 | +81 45 339 5877  |                                          |
| <b>KAZAKHSTAN</b>     | +7 495 935 8888                 | +7 495 564 8787  | info@scheltec.ru                         |
| <b>KOREA</b>          | 82 2 868 5376                   | 82 2 868 4343    | contact.KR@perkinelmer.com               |
| <b>KYRGYZSTAN</b>     | +7 495 935 8888                 | +7 495 564 8787  | info@scheltec.ru                         |
| <b>MALAYSIA</b>       | 603 7949 1118                   | 603 7949 1119    | www.perkinelmer.com/category/consumables |
| <b>MONGOLIA</b>       | +7 495 935 8888                 | +7 495 564 8787  | info@scheltec.ru                         |
| <b>NEW ZEALAND</b>    | 64 21 215 1198 or 0800 724 663  | 64 3 384 1222    | www.perkinelmer.com/category/consumables |
| <b>PAKISTAN</b>       | +92 21 5345581                  | +92 21 5345582   | sales@amspvtltd.com                      |

# ORDERING INFORMATION

## ASIA

| Country                         | Telephone          | Facsimile        | Email                                                                                                  |
|---------------------------------|--------------------|------------------|--------------------------------------------------------------------------------------------------------|
| <b>PHILIPPINES</b>              | 632 822 0511 to 16 | 632 822 0517     | <a href="http://www.perkinelmer.com/category/consumables">www.perkinelmer.com/category/consumables</a> |
| <b>RUSSIAN FEDERATION</b>       | +7 495 935 8888    | +7 495 564 8787  | <a href="mailto:info@scheltec.ru">info@scheltec.ru</a>                                                 |
| <b>SINGAPORE</b>                | 65 6868 1650       | 65 6872 6595     | <a href="http://www.perkinelmer.com/category/consumables">www.perkinelmer.com/category/consumables</a> |
| <b>TAIWAN</b> Taipei            | 886 2 87912589     | 886 2 87919085   | <a href="http://www.perkinelmer.com/category/consumables">www.perkinelmer.com/category/consumables</a> |
| <b>TAIWAN</b> Kaohsiung         | 886 7 5521030      | 886 7 5543402    | <a href="http://www.perkinelmer.com/category/consumables">www.perkinelmer.com/category/consumables</a> |
| <b>TAJIKISTAN</b>               | +7 495 935 8888    | +7 495 564 8787  | <a href="mailto:info@scheltec.ru">info@scheltec.ru</a>                                                 |
| <b>THAILAND</b>                 | 662 319 7901       | 662 319 7900     | <a href="http://www.perkinelmer.com/category/consumables">www.perkinelmer.com/category/consumables</a> |
| <b>TURKEY</b>                   | +90 212 312 1100   | +90 212 312 1111 | <a href="mailto:cc.turkey@perkinelmer.com">cc.turkey@perkinelmer.com</a>                               |
| <b>TURKMENISTAN</b>             | +7 495 935 8888    | +7 495 564 8787  | <a href="mailto:info@scheltec.ru">info@scheltec.ru</a>                                                 |
| <b>VIETNAM</b> Hanoi            | 84 4 3747 2258     | 84 4 3747 2260   | <a href="http://www.perkinelmer.com/category/consumables">www.perkinelmer.com/category/consumables</a> |
| <b>VIETNAM</b> Ho Chi Minh City | 84 8 5431 8877     | 84 8 5431 8570   | <a href="http://www.perkinelmer.com/category/consumables">www.perkinelmer.com/category/consumables</a> |
| <b>UZBEKISTAN</b>               | +7 495 935 8888    | +7 495 564 8787  | <a href="mailto:info@scheltec.ru">info@scheltec.ru</a>                                                 |

## MIDDLE EAST AND GULF

| Country                        | Telephone           | Facsimile           | Email                                                                    |
|--------------------------------|---------------------|---------------------|--------------------------------------------------------------------------|
| <b>BAHRAIN</b>                 | +973 17276176       | +973 17275819       | <a href="mailto:info@ymh.com.bh">info@ymh.com.bh</a>                     |
| <b>GAZA STRIP</b>              | +202 22905306       | +202 22905307       | <a href="mailto:info@htds.fr">info@htds.fr</a>                           |
| <b>IRAQ</b>                    | +202 22905306       | +202 22905307       | <a href="mailto:info@htds.fr">info@htds.fr</a>                           |
| <b>ISRAEL</b>                  | 1700 500512         | 1700 555 915        | <a href="mailto:cc.israel@perkinelmer.com">cc.israel@perkinelmer.com</a> |
| <b>JORDAN</b>                  | +962 (0) 6 534 6523 | +962 (0) 6 534 6527 | <a href="mailto:info@htds.fr">info@htds.fr</a>                           |
| <b>KUWAIT</b>                  | +965 22419797       | +965 22419798       | <a href="mailto:las@emphor.biz">las@emphor.biz</a>                       |
| <b>LEBANON</b>                 | +202 22905306       | +202 22905307       | <a href="mailto:info@htds.fr">info@htds.fr</a>                           |
| <b>OMAN</b>                    | +968 248 14752      | +968 248 13924      | <a href="mailto:big@bigllcoman.com">big@bigllcoman.com</a>               |
| <b>PALESTINIAN TERRITORIES</b> | +202 22905306       | +202 22905307       | <a href="mailto:info@htds.fr">info@htds.fr</a>                           |
| <b>QATAR</b>                   | +974 4654563        | +974 4652256        | <a href="mailto:las@emphor.biz">las@emphor.biz</a>                       |
| <b>SAUDI ARABIA</b>            | +966 11477 7932     | +966 11477 1314     | <a href="mailto:info@aralgosaibico.com">info@aralgosaibico.com</a>       |
| <b>UNITED ARAB EMIRATES</b>    | +971 4 883 0233     | +971 4 883 0133     | <a href="mailto:las@emphor.biz">las@emphor.biz</a>                       |
| <b>WEST BANK</b>               | +202 22905306       | +202 22905307       | <a href="mailto:info@htds.fr">info@htds.fr</a>                           |
| <b>YEMEN</b>                   | +967 234 5073       | +967 234 0504       | <a href="mailto:direct@alfalahye.net">direct@alfalahye.net</a>           |

| AFRICA               |                   |                   |                         |
|----------------------|-------------------|-------------------|-------------------------|
| Country              | Telephone         | Facsimile         | Email                   |
| ALGERIA              | +213 23 23 84 01  | +213 23 23 84 00  | info@htds.fr            |
| ANGOLA               | +27 11 564 2400   | +27 11 564 2440   | info.za@perkinelmer.com |
| BENIN                | +225 78 78 69 32  | +33 1 69 07 69 54 | info@htds.fr            |
| BOTSWANA             | +27 11 564 2400   | +27 11 564 2440   | info.za@perkinelmer.com |
| BRITISH INDIAN OCEAN | +202 22905306     | +202 22905307     | info@htds.fr            |
| BURKINA FASO         | +225 78 78 69 32  | +33 1 69 07 69 54 | info@htds.fr            |
| BURUNDI              | +243 990 086 063  | +33 1 69 07 69 54 | info@htds.fr            |
| CAMEROON             | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |
| CAPE VERDE           | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |
| CENTRAL AFRICA REP.  | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |
| CHAD                 | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |
| COMOROS              | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |
| CONGO                | +243 990 086 063  | +33 1 69 07 69 54 | info@htds.fr            |
| DEMO. REP. CONGO     | +243 990 086 063  | +33 1 69 07 69 54 | info@htds.fr            |
| DSCHIBUTI REP.       | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |
| EGYPT                | +202 22905306     | +202 22905307     | info@htds.fr            |
| EQUATORIAL GUINEA    | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |
| ERITREA              | +202 22905306     | +202 22905307     | info@htds.fr            |
| ETHIOPIA             | +202 22905306     | +202 22905307     | info@htds.fr            |
| GABON                | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |
| GAMBIA               | +225 78 78 69 32  | +33 1 69 07 69 54 | info@htds.fr            |
| GHANA                | +202 22905306     | +202 22905307     | info@htds.fr            |
| GUINEA               | +225 78 78 69 32  | +33 1 69 07 69 54 | info@htds.fr            |
| GUINEA – Bissau      | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |
| IVORY COAST          | +225 78 78 69 32  | +33 1 69 07 69 54 | info@htds.fr            |
| KENYA                | +202 22905306     | +202 22905307     | info@htds.fr            |
| LESOTHO              | +27 11 564 2400   | +27 11 564 2440   | info.za@perkinelmer.com |
| LIBERIA              | +202 22905306     | +202 22905307     | info@htds.fr            |
| LIBYA                | +218 91 695 0708  | +218/21 713 21 12 | info@htds.fr            |
| MADAGASCAR           | +33 1 64 86 28 28 | +33 1 69 07 69 54 | info@htds.fr            |

| AFRICA                          |                      |                      |                         |
|---------------------------------|----------------------|----------------------|-------------------------|
| Country                         | Telephone            | Facsimile            | Email                   |
| <b>MALAWI</b>                   | +27 11 564 2400      | +27 11 564 2440      | info.za@perkinelmer.com |
| <b>MALI</b>                     | +225 78 78 69 32     | +33 1 69 07 69 54    | info@htds.fr            |
| <b>MAURITANIA</b>               | +33 1 64 86 28 28    | +33 1 69 07 69 54    | info@htds.fr            |
| <b>MAURITIUS</b>                | +27 11 564 2400      | +27 11 564 2440      | info.za@perkinelmer.com |
| <b>MOROCCO</b>                  | +212/(0)522 27 49 59 | +212/(0)522 20 83 74 | info@htds.fr            |
| <b>MOZAMBIQUE</b>               | +27 11 564 2400      | +27 11 564 2440      | info.za@perkinelmer.com |
| <b>NAMIBIA</b>                  | +27 11 564 2400      | +27 11 564 2440      | info.za@perkinelmer.com |
| <b>NIGER</b>                    | +225 78 78 69 32     | +33 1 69 07 69 54    | info@htds.fr            |
| <b>NIGERIA</b>                  | +202 22905306        | +202 22905307        | info@htds.fr            |
| <b>REUNION ISLANDS</b>          | +33 1 64 86 28 28    | +33 1 69 07 69 54    | info@htds.fr            |
| <b>RWANDA</b>                   | +33 1 64 86 28 28    | +33 1 69 07 69 54    | info@htds.fr            |
| <b>SENEGAL</b>                  | +225 78 78 69 32     | +33 1 69 07 69 54    | info@htds.fr            |
| <b>SEYCHELLES</b>               | +33 1 64 86 28 28    | +33 1 69 07 69 54    | info@htds.fr            |
| <b>SIERRA LEONE</b>             | +202 22905306        | +202 22905307        | info@htds.fr            |
| <b>SOMALIA</b>                  | +202 22905306        | +202 22905307        | info@htds.fr            |
| <b>REPUBLIC OF SOUTH AFRICA</b> | +27 11 564 0600      | +27 11 564 0899      | info.za@perkinelmer.com |
| <b>SWAZILAND</b>                | +27 11 564 2400      | +27 11 564 2440      | info.za@perkinelmer.com |
| <b>TANZANIA</b>                 | +202 22905306        | +202 22905307        | info@htds.fr            |
| <b>TOGO</b>                     | +225 78 78 69 32     | +33 1 69 07 69 54    | info@htds.fr            |
| <b>TUNISIA</b>                  | +216 70 836 961      | +216 70 836 561      | info@htds.fr            |
| <b>UGANDA</b>                   | +202 22905306        | +202 22905307        | info@htds.fr            |
| <b>WESTERN SAHARA</b>           | +33 1 64 86 28 28    | +33 1 69 07 69 54    | info@htds.fr            |
| <b>ZAMBIA</b>                   | +27 11 564 2400      | +27 11 564 2440      | info.za@perkinelmer.com |
| <b>ZIMBABWE</b>                 | +27 11 564 2400      | +27 11 564 2440      | info.za@perkinelmer.com |

**PerkinElmer, Inc.**  
940 Winter Street  
Waltham, MA 02451 USA  
P: (800) 762-4000 or  
(+1) 203-925-4602  
[www.perkinelmer.com](http://www.perkinelmer.com)



---

For a complete listing of our global offices, visit [www.perkinelmer.com/ContactUs](http://www.perkinelmer.com/ContactUs)

Copyright ©2022, PerkinElmer, Inc. All rights reserved. PerkinElmer is a trademark of PerkinElmer, Inc. All other trademarks are the property of their respective owners.

547522